

ELISA for Human IFN- α (pan specific)

Product Code: 3425-1A-6

CONTENTS, development kit for 6 plates:

Vial 1 (yellow top)

Monoclonal antibodies MT1/3/5 a combination of three different antibodies (300 μ l)

Concentration: 1 mg/ml

Vial 2 (red top)

Biotinylated monoclonal antibodies MT2/4/6 (80 μ l)

Concentration: 1 mg/ml

Vial 3 (white top)

Streptavidin-Alkaline Phosphatase (80 μ l)

Vial 4

Lyophilised native human IFN- α standard (0.5 μ g)

To ensure total recovery of stated quantity, vials have been overfilled.

STORAGE:

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored at -20°C.

General

Intended use: For quantitative determination of native and recombinant human IFN- α in solution, e.g. cell culture supernatant. The system will detect native and recombinant human IFN- α subtypes 1/13, 2, 4, 5, 6, 7, 8, 10, 14, 16 and 17.


Serum/plasma samples: Please note that determination of analyte in human serum/plasma samples by this kit requires the use of ELISA diluent (product code: 3652-D2) for dilution of samples, standard and detection antibody. The diluent prevents false positive read-outs which may be caused by interference of heterophilic antibodies commonly found in human plasma and serum. The ELISA diluent has been validated using serum/plasma from normal healthy human blood donors. Please note that heterophilic antibody interference in samples from human subjects with various diseases or other conditions has not been assessed. Please contact Mabtech for further information.

Reagents: Antibodies are supplied in sterile-filtered (0.2 μ m) PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.15% Kathon CG.

Standard range: 7-700 pg/ml

Standard calibration: 1 ng of supplied standard equals 400 U of 94/784 NIBSC*-standard according to repeated calibrations. Calibration is batch-specific.

*National Institute of Biological Standards and Control, UK.

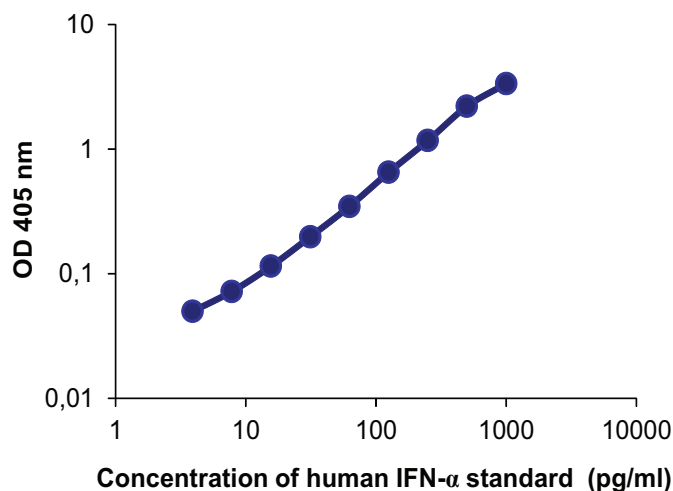

Guidelines for Human IFN- α

(pan specific) ELISA

- Day 1**
1. Coat a high protein binding ELISA plate with mAbs MT1/3/5, diluted to 4 $\mu\text{g}/\text{ml}$ in PBS, pH 7.4, by adding 100 μl /well. Incubate overnight at 4-8°C.
- Day 2**
2. Wash twice with PBS (200 μl /well).
 3. Block plate by adding 200 μl /well of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
 4. Wash five times with PBS containing 0.05% Tween20
 5. Prepare hIFN- α standard by reconstituting contents of vial 4 in 1 ml PBS with 0.1% BSA to make up a stock solution of 500 ng/ml. Use immediately or store in aliquots at -20°C for future use. We recommend the aliquots not to be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
 6. Add 100 μl /well of samples or standards diluted in incubation buffer or ELISA diluent for serum/plasma samples and incubate for 2 hours at room temperature.
 7. Wash as in step 4.
 8. Add 100 μl /well of mAbs MT2/4/6-biotin diluted to 1 $\mu\text{g}/\text{ml}$ in incubation buffer or ELISA diluent for serum/plasma samples. Incubate for 1 hour at room temperature.
 9. Wash as in step 4.
 10. Add 100 μl /well of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
 11. Wash as in step 4.
 12. Add 100 μl /well of appropriate substrate solution e.g. p-nitro-phenyl-phosphate (pNPP).
 13. Measure the optical density (405 nm for pNPP) in an ELISA reader after suitable developing time.

NOTE; for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



MABTECH AB
Box 1233
SE-131 28 Nacka Strand
Sweden
Tel: +46 8 716 27 00
Fax: +46 8 716 27 01
E-mail: mabtech@mabtech.com
www.mabtech.com

MABTECH Inc
M.E.B. 220
3814 West Street
Cincinnati, OH 45227
USA
Tel: +1 513 871 4500
Fax: +1 513 871 7353
E-mail: mabtech.usa@mabtech.com

MABTECH AB Büro Deutschland
Germany
Tel: +49 40 4135 7935
Fax: +49 40 4135 7945
E-mail: mabtech.de@mabtech.com

2013-01-24

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the following standards:



MABTECH AUSTRALIA Pty Ltd
resolvingIMAGES
Unit 22, 196 Settlement Road
Thomastown Victoria 3074
Australia
Tel: +61 3 9466 4007
Fax: +61 3 9466 4003
E-mail: mabtech.au@mabtech.com

MABTECH AB Bureau de liaison France
BP 255, 1300 route des Crêtes
06905 Sophia Antipolis
France
Tel: +33 (0)4 92 38 80 70
Fax: +33 (0)4 92 38 80 71
E-mail: mabtech.fr@mabtech.com