



TABLE OF CONTENTS

1. Kit Contents ..... 2

2. Materials Required but Not Supplied..... 2

3. Introduction ..... 3

    3.1 Principle of FRET and TR-FRET ..... 3

    3.2 Assay Overview..... 3

4. Instrument Settings ..... 4

    4.1 General Considerations ..... 4

    4.2 General Settings for Filter-based Instruments..... 5

    4.3 General Settings for Monochromator-based Instruments..... 6

    4.4 General Settings for CCD-based Instruments ..... 6

5. Guidelines and Recommendations ..... 6

    5.1 Reagent Handling..... 6

    5.2 Ligand Dilutions ..... 7

    5.3 Solvent Tolerance..... 7

    5.4 Reagent Order of Addition ..... 7

    5.5 Incubation Conditions ..... 7

6. Procedure ..... 8

    6.1 Assay Set-Up ..... 8

    6.2 Plate Layout..... 10

    6.3 Data Analysis ..... 10

    6.4 Reagent Volumes and Controls..... 11

    6.5 Ligand Assay—Representative Data ..... 12

7. Assay Pharmacology ..... 13

8. References..... 14

9. Limited Product Warranty ..... 14

10. Purchaser Notification ..... 14

For Research Use Only. Not for use in diagnostic procedures.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

## 1. KIT CONTENTS

The LanthaScreen® TR-FRET ER Alpha Competitive Binding Assay kit, catalog no. A15887, contains the following:

Component	Composition	Amount	Storage Temperature	Individual Part no.
Fluormone™ ES2 Green	1,800nM in 80% methanol/water	40 µL	-20°C	PV6045
ER Alpha-LBD (GST)	Human ER Alpha ligand-binding domain in buffer (pH 8.0) containing protein stabilizing reagents and 20% glycerol. See Certificate of Analysis for lot specific concentration.	25 µg	-80°C	A15677
LanthaScreen® Tb-anti-GST antibody	Terbium labeled anti-GST antibody in HEPES buffered saline (137 mM NaCl, 2.7 mM KCl, 10 mM HEPES pH 7.5). See Certificate of Analysis for lot specific concentration	25 µg	-20°C	PV3550
Nuclear Receptor Buffer K	Proprietary buffer (pH 7.4)	2 × 25 mL	Room Temp.	PV6043
DTT, 1 M	In water	1 mL	-20°C	P2325

## 2. MATERIALS REQUIRED BUT NOT SUPPLIED

The following materials are required but not supplied in the kit:

- A fluorescence plate reader with excitation capabilities at 340 nm and with the appropriate filter sets installed for detecting the fluorescent emission signals of terbium at 495 nm and fluorescein at 520 nm (see **Section 4**).
- Pipetting devices for 1–1000 µL volumes, suitable repeater pipettors, or multi-channel pipettor.
- Black, 384-well assay plates. We recommend black Corning® 384-well, low-volume, round-bottom, black (non-treated surface) assay plates, Cat. no. 4511. Other plate types may give satisfactory results as well but have not been fully tested.
- 96-well polypropylene plates that can accommodate a 300-µL volume per well for the ligand serial dilutions in 100% DMSO. Polypropylene plates are needed since they are tolerant to DMSO. We recommend Nalge Nunc, Cat. no. 249944.
- A known ER ligand, such as estradiol to serve as a positive control. We recommend Sigma, Cat. no. E1024.
- DMSO to perform serial dilutions. We recommend Sigma, Cat. no. 41647.
- Recommended for instrument testing, a fluorescein TR-FRET acceptor such as Fluorescein-Poly GT, Cat. no. PV3610.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

### 3. INTRODUCTION

---

The LanthaScreen® TR-FRET ER Alpha Competitive Binding Assay provides a sensitive and robust method for high-throughput screening of ligands for ER Alpha. The kit uses a terbium-labeled anti-GST antibody, a fluorescent small molecule ER Alpha ligand (Fluormone™ ES2 Green, also referred to as “tracer”), and a human ER Alpha ligand-binding domain (LBD) that is tagged with glutathione-S-transferase (GST) in a homogenous mix-and-read assay format. This kit contains enough reagents for 800 assays of 20 µL each.

#### 3.1 Principle of FRET and TR-FRET

For screening libraries of compounds, time-resolved FRET (TR-FRET) is a recognized method for overcoming interference from compound autofluorescence or light scatter from precipitated compounds. The premise of a TR-FRET assay is the same as that of a standard FRET assay: when a suitable pair of fluorophores is brought within close proximity of one another, excitation of the first fluorophore (the donor) can result in energy transfer to the second fluorophore (the acceptor). This energy transfer is detected by an increase in the fluorescence emission of the acceptor and a decrease in the fluorescence emission of the donor. In HTS assays, FRET is often expressed as a ratio of the intensities of the acceptor and donor fluorophores. The ratiometric nature of such a value corrects for differences in assay volumes between wells and corrects for quenching effects due to colored compounds.

In contrast to standard FRET assays, TR-FRET assays use a long-lifetime lanthanide chelate as the donor species. Lanthanide chelates are unique in that their excited-state lifetime (the average time that the molecule spends in the excited state after accepting a photon) can be on the order of a millisecond or longer. This is in sharp contrast to the lifetime of common fluorophores used in standard FRET assays, which are typically in the nanosecond range. Because interference from autofluorescent compounds or scattered light from precipitated compounds is also on the nanosecond timescale, these factors can negatively impact standard FRET assays. To overcome these interferences, TR-FRET assays are performed by measuring FRET after a suitable delay, typically 100 microseconds after excitation by a flashlamp excitation source in a microtiter plate reader. This delay not only overcomes interference from background fluorescence or light scatter, but also avoids interference from direct excitation due to the non-instantaneous nature of the flashlamp excitation source.

The most common lanthanides used in TR-FRET assays for HTS are terbium (Tb) and europium (Eu). Terbium offers unique advantages over europium when used as the donor species in a TR-FRET assay. In contrast to europium-based systems that employ the relatively large protein APC as the acceptor, terbium-based TR-FRET assays can use common small molecule fluorophores such as fluorescein as the acceptor. It is straightforward (and inexpensive) to directly label a molecule such as a peptide with fluorescein for use in terbium-based TR-FRET assays. In addition, the assay format is simpler than using a biotinylated peptide that would then be indirectly labeled via streptavidin-mediated recruitment of APC. Therefore, the use of directly labeled molecules (in this case the fluorescent ligand) in a terbium-based TR-FRET assay simplifies assay development by reducing the number of components to optimize, avoids problems due to steric interactions involving large APC conjugates, reduces costs, and improves kinetics.

#### 3.2 Assay Overview

A terbium-labeled anti-GST antibody is used to indirectly label a nuclear (NR) by binding to its GST tag. Competitive binding to the NR is detected by a test compound's ability to displace a fluorescent ligand (tracer) from the NR, which results in a loss of FRET signal between the Tb-anti-GST antibody and the tracer. This type of binding assay is analogous to radioligand-based assays, except that it eliminates handling of radioactivity and enables a homogenous, “addition-only” format.

When running the LanthaScreen® TR-FRET ER Alpha Competitive Binding Assay, Fluormone™ ES2 Green tracer is added to ligand test compounds or solvent controls followed by addition of a mixture of the ER Alpha-LBD (GST) and terbium anti-GST antibody. After an incubation period at room temperature, the TR-FRET ratio of 520:495 emissions is calculated and can be used to determine the IC<sub>50</sub> from a dose response curve of the compound.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266


For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

The recommended concentration of ER Alpha-LBD (GST) and Tb anti-GST antibody have been optimized to produce a satisfactory assay window ( $Z' > 0.5$ ) while using the lowest receptor concentration to achieve the best sensitivity to tight binding ligands. The assay window is highly dependent on the sensitivity of the fluorescent plate reader. Different plate readers may provide different assay windows and  $Z'$  factors. These assays were optimized using a BMG LABTECH® PHERAstar® fluorescent plate reader, which is considered very sensitive for these assays. If the assay window is not acceptable on your instrument, additional receptor or receptor and antibody may be added. Please note that this will decrease the sensitivity of the assay to tight binders.

**Figure 1. Principle of a LanthaScreen® TR-FRET NR Competitive Binding assay:** Tb-anti-GST antibody indirectly labels the nuclear receptor by binding to the GST tag. Binding of the fluorescently labeled tracer results in close proximity of the tracer to the terbium-labeled antibody, resulting in a high TR-FRET signal. When competitor displaces the fluorescently labeled tracer, there is a decrease in the TR-FRET signal.

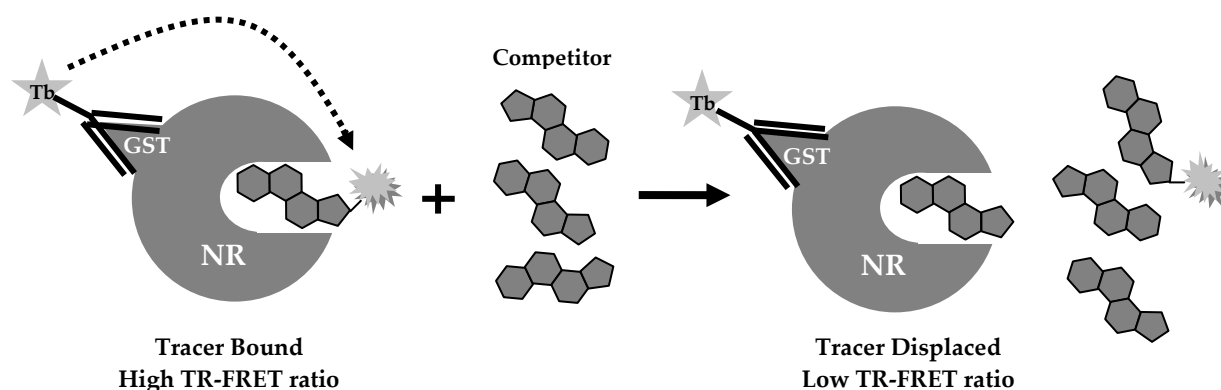

## 4. INSTRUMENT SETTINGS

To facilitate instrument set-up we recommend downloading the Terbium Assay Application Note which can be found at [www.lifetechnologies.com/instrumentsetup](http://www.lifetechnologies.com/instrumentsetup) and following the instructions. You will need Fluorescein-Poly GT, PV3610 to run this test. For additional help in setting up your instrument to perform LanthaScreen® TR-FRET assays, please contact Technical Support at 800-955-6288 ext. 40266, or email [drugdiscoverytech@lifetech.com](mailto:drugdiscoverytech@lifetech.com). Please see [www.lifetechnologies.com/instrumentsetup](http://www.lifetechnologies.com/instrumentsetup) for instrument specific information, including step-by-step instrument specific guides for optimizing LanthaScreen® Terbium assays on your particular instrument.

The LanthaScreen® technology has been tested on a variety of microplate readers. These include filter-based instruments, Tecan® ULTRA, GENios™ pro, and Infinite® F500, BMG LABTECH® PHERAstar®; BioTek® Synergy™ 2 and Synergy™ 4; PerkinElmer® EnVision® and VICTOR™, and Molecular Devices® Analyst®; CCD-based imagers (PerkinElmer® ViewLux®) and monochromator-based instruments (Tecan® Safire™ 2 and Infinite® M1000).

### 4.1 General Considerations

The excitation and emission spectra of terbium and fluorescein are shown in Figure 2, next page. To read a LanthaScreen® TR-FRET assay, the instrument is configured to excite the terbium donor around 340 nm, and to separately read the terbium emission peak that is centered at approximately 490 nm, and the fluorescein emission that is centered at approximately 520 nm. Separation of the terbium emission signal from the fluorescein emission signal is critical for assay success and is achieved with the proper selection of filter bandwidths which cannot be compromised. For this reason a standard fluorescein filter cannot typically be used. After taking the measurements, the signal from the fluorescein emission is divided (or "ratioed") by the terbium signal to provide a TR-FRET emission ratio and improve the  $Z'$  for the assay.


For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

Figure 2. Excitation and emission spectra of fluorescein and terbium.

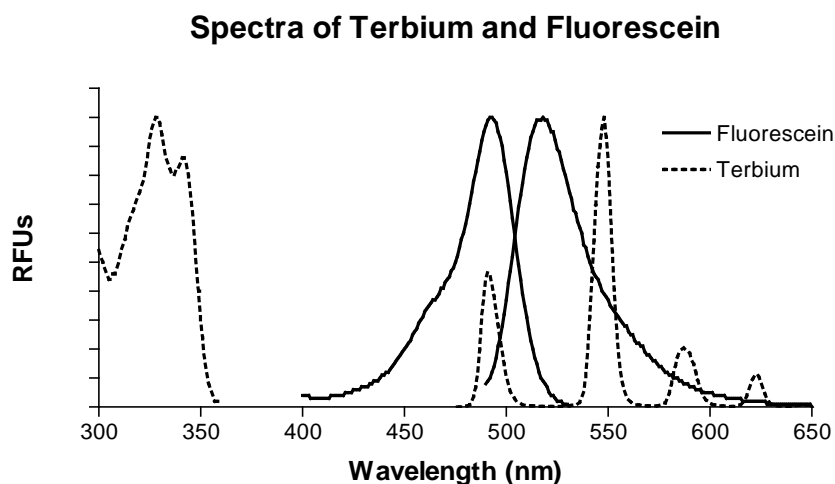

Aside from filter choices, other instrument settings are similar to the settings used with europium-based technologies. A delay time of 100  $\mu$ s followed by a 200- $\mu$ s integration time is typical for a LanthaScreen® TR-FRET assay. The number of flashes or measurements per well is highly instrument dependent and should be set as advised by your instrument manufacturer.

## 4.2 General Settings for Filter-based Instruments

Excitation	340 nm filter (30 nm bandwidth)
Fluorescein Emission	520 nm filter (25 nm bandwidth)
Terbium Emission	490 or 495 nm filter (10 nm bandwidth)
Dichroic Mirror	Fluorescein (Tecan® ULTRA, GENios™, GENios™ pro, Infinite®) PerkinElmer® LANCE®/TRF (EnVision®, VICTOR™) 380 nm (preferred) or 400 nm (Molecular Devices® Analyst®) Built In (BMG LABTECH® PHERAstar®)
Delay Time	100 $\mu$ s
Integration Time	200 $\mu$ s

Note that excitation filters with similar bandwidths will give satisfactory performance as long as the center wavelength falls at or between 330 nm and 340 nm. Excitation and emission filters (from Chroma Technologies) for most microplate readers are available directly from Life Technologies. Visit [www.lifetechnologies.com/instrumentsetup](http://www.lifetechnologies.com/instrumentsetup) to purchase filters for your instrument. A LanthaScreen® Tb filter optical module is available directly from BMG LABTECH® for use on the BMG LABTECH® PHERAstar® instrument. PerkinElmer® EnVision® instrument users may require a specific filter holder which is available from PerkinElmer® as part no. 2100-8110 for a package of 10.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

### 4.3 General Settings for Monochromator-based Instruments

Excitation	332 nm (20 nm bandwidth)
Fluorescein Emission	515 nm (20 nm bandwidth)
Terbium Emission	486 nm (20 nm bandwidth)
Delay Time	100 $\mu$ s
Integration Time	200 $\mu$ s

Some monochromator-based instruments, such as the Tecan® Infinite® M1000, allow each emission wavelength to be individually set. In this case, match the settings to those listed for the filter-based instruments. When the bandwidths of the two emission wavelengths cannot be individually set, use the general settings for monochromator-based instruments listed in the table above as a guide.

Note that not all monochromator-based instruments are capable of a 20 nm bandwidth. Smaller bandwidth settings may be used, but with a decrease in assay performance. Additionally, we have found that while some monochromator based instruments (Tecan® Safire<sup>2</sup>™, Tecan® Infinite® M1000) give satisfactory performance in LanthaScreen® TR-FRET assays, other monochromator-based instruments (e.g. Tecan® Safire™ and Molecular Devices® Gemini™, and SpectraMax® M2 and M5) may not be optimal for detection, although some change in signal and a ratio may still be observed. We have also found that with some assays, white plates give better assay performance when using monochromator-based instruments. In filter-based instruments the difference is typically negligible.

### 4.4 General Settings for CCD-based Instruments

Please see [www.lifetechnologies.com/instrumentsetup](http://www.lifetechnologies.com/instrumentsetup) and click on "ViewLux®" for a step-by-step guide to optimizing LanthaScreen® Terbium assays on the PerkinElmer® ViewLux® instrument.

## 5. GUIDELINES AND RECOMMENDATIONS

### 5.1 Reagent Handling

#### *ER Alpha-LBD (GST)*

Store ER Alpha-LBD (GST) at –80°C. Thaw on ice before use and perform all dilutions while on ice. Mix by gentle pipetting or inversion of the tube. Never vortex the ER Alpha-LBD (GST) stock or dilutions.

#### *Nuclear Receptor Buffer K*

Thaw Nuclear Receptor Buffer K at room temperature upon receipt. Mix well before first use, as the buffer may not have thawed evenly. The buffer should be stored at room temperature.

#### *LanthaScreen® Tb anti-GST Antibody*

Store Tb anti-GST antibody at –20°C and thaw before use.

**Important:** To help minimize the potential effect of spurious (random) donor emission spikes in your assay readout, centrifuge the stock vial of Tb anti-GST Antibody (~10,000 rpm for 10 minutes) prior to use. After centrifugation, pipet the quantity of antibody needed for your assays from the top of the liquid, thereby minimizing the potential mixing of any precipitate that has been spun to the bottom of the tube.

#### *Fluormone™ ES2 Green*

Store Fluormone™ ES2 Green at –20 °C. Place on ice for at least 15 minutes and then vortex before use. Keep container closed and on ice when not in use to prevent evaporation.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

## 5.2 Ligand Dilutions

The procedure described below involves diluting the ligand to 100X in DMSO followed by transfer into complete assay buffer, resulting in a 2X ligand, 2% DMSO dilution. This may be done in a DMSO-tolerant assay plate or tubes such as those made from polypropylene or HDPE. Dilution in DMSO reduces the occurrence of compound precipitation while performing the dilution series, reduces sticking to plastics, improves IC<sub>50</sub> reproducibility and maintains constant DMSO in all wells. After addition of all reagents to the assay, the final concentration will be 1X ligand and 1% DMSO. We suggest using this method for compound dilutions.

**Note:** Handling of some ligands can be problematic due to their tendency to stick to various plastics. These ligands may show varying degrees of stickiness, causing differences in the actual concentration of the ligand, depending on the type of pipette tips and plates used in the set up of the assay. Therefore, you may observe different IC<sub>50</sub> values than reported here.

## 5.3 Solvent Tolerance

The assay was validated in the presence of 1% DMSO. However, the assay has been performed with up to 4% DMSO, 8% ethanol, 8% acetonitrile and 8% methanol (in addition to the 1% DMSO present from the ligand dilution) with good results. The use of higher concentrations of DMSO may still yield a robust assay, however the observed assay window .

## 5.4 Reagent Order of Addition

The LanthaScreen® ER Alpha Competitive Binding Assay has been developed using separate tracer and receptor/antibody additions. Use of a pre-mixture containing all three components is not recommended. If the receptor is first allowed to bind to the tracer in the absence of competitor, more time will be required for binding of the competitor to reach equilibrium, and the assay window will be reduced due to higher a higher background signal.

## 5.5 Incubation Conditions

### *Incubation Time*

The incubation time can be set by the user. As a guide, results for various time points using estradiol as the ligand are shown in Table 1 below. The IC<sub>50</sub> and Z' are stable from 1 to 6 hours (Table 1). We suggest maintaining the same incubation time and temperature during testing and screening for consistent results.

**Table 1. Effect of Incubation Time on Assay Performance.** Sample data represents mean values from 4 separate experiments (n = 24). IC<sub>50</sub> values were determined by fitting the data to a sigmoidal dose response (variable slope) equation in IDBS XLfit® 4.2 (data not shown). Z'-factors were calculated using the method of Zhang *et al.* (Zhang *et al.*, 1999) on the 24 replicates of maximum ligand and no ligand. Z'-factor is an indication of the robustness of the assay, where values ≥ 0.5 indicate an excellent assay, while a value of 1 indicates a theoretically ideal assay with no variability.

Incubation Time (hours)	IC <sub>50</sub> Estradiol (nM)	Z'-Factor
1	0.56	0.62
2	0.55	0.61
4	0.57	0.61
6	0.57	0.63

### *Temperature*

We recommend that assays be conducted at room temperature (20–23°C).

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

## 6. PROCEDURE

The procedure in this section describes a method for determining the  $IC_{50}$  of a ligand and the  $Z'$  factor for maximum ligand and no ligand controls ( $n = 24$ ). The only variable is the ligand concentration. All other assay components (ER Alpha-LBD (GST), tracer, Tb anti-GST antibody) are fixed at concentrations optimized to produce a satisfactory assay window ( $Z'$ -factor  $>0.5$ ) while using the lowest receptor concentration to achieve the best sensitivity to tight binding ligands. Assays were measured on a BMG LABTECH® PHERAstar® fluorescent plate reader. Different plate readers may provide very different assay windows because assay window is highly dependent on fluorescent plate reader sensitivity. Higher concentrations of nuclear receptor may give a larger TR-FRET signal window, however it will decrease the sensitivity of the assay with regard to differentiating tight binding ligands. The recommended final concentrations are listed in the following table. If a component concentration is changed, the concentrations of the other components may need to be reoptimized.

Component	Recommended Final Assay Concentrations
Fluormone™ ES2 Green	3 nM
Tb anti-GST antibody	2 nM
ER Alpha-LBD (GST)	See Certificate of Analysis for the recommended molar concentration for this kit

### 6.1 Assay Set-Up

#### *Prepare Complete Nuclear Receptor Buffer K and Ligand Controls*

**Note:** The ER Alpha-LBD (GST) should be thawed on ice just prior to use. During the assay set-up, keep the ER Alpha-LBD (GST) and Fluormone™ ES2 Green on ice. Equilibrate all other assay components to room temperature.

1. Prepare Complete Nuclear Receptor Buffer K by adding 1 M DTT to Nuclear Receptor Buffer K for a final concentration of 5 mM DTT. Complete Nuclear Receptor Buffer K must be prepared fresh daily.

*For example:* Add 30  $\mu$ L of 1 M DTT to 5.97 mL of Nuclear Receptor Buffer K.

**Note:** Buffer K is a specially formulated buffer that has been optimized to work with this kit. To ensure performance of this assay, we highly recommend using Buffer K as provided with the kit.

2. For the “no ligand” controls, add DMSO to Complete Nuclear Receptor Buffer K for a final concentration of 2% DMSO (2X solvent). Add 10  $\mu$ L of this buffer containing DMSO to row C, columns 1–24 of a 384-well assay plate (see the plate layout in **Section 6.2**).

*For example:* Add 10  $\mu$ L of DMSO to 490  $\mu$ L of Complete Nuclear Receptor Buffer K.

3. Prepare a solution of control ligand (we recommend estradiol) at 100X of the final desired maximum starting concentration using 100% DMSO.

*For example:* If the final desired maximum starting concentration of ligand is 1  $\mu$ M, prepare a solution of 100  $\mu$ M ligand in 100% DMSO.

4. For the “maximum ligand” controls, dilute the 100X ligand solution from step 3 to 2X using Complete Nuclear Receptor Buffer K. Add 10  $\mu$ L of this solution to row D, columns 1–24 in the 384-well assay plate (see the plate layout in **Section 6.2**).

*For example:* Add 10  $\mu$ L of 100X ligand solution to 490  $\mu$ L of Complete Nuclear Receptor Buffer K.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)



### Prepare 2X Ligand Dilution Series

**Note:** Although steps 5 and 6 below require more pipetting than other methods of preparing a serial dilution of ligand, we have found that this approach provides a robust method for preparing the dilution series without problems due to ligand solubility. Dilution with 100% DMSO facilitates compound solubility during serial dilutions, which is important for obtaining consistent IC<sub>50</sub> values.

5. Prepare a 12-point 100X dilution series of ligand in a 96-well polypropylene plate (DMSO tolerant) by serially diluting the 100X ligand solution from step 3 using 100% DMSO. We recommend a 3-fold dilution series.

*For example:* Add 20 µL of 100% DMSO to wells A2–A12 in a 96-well polypropylene plate. To well A1, add 30 µL of the 100X ligand solution prepared in step 3. Perform a three-fold serial dilution by transferring 10 µL of the 100X ligand solution from well A1 to the 20 µL of DMSO in well A2. Mix by pipetting up and down. Repeat for wells A2–A12.

6. Dilute each 100X agonist serial dilution from step 5 to 2X using Complete TR-FRET Coregulator Buffer K.

*For example:* Transfer 5 µL of each of the 100X agonist serial dilutions from row A of the 96-well plate (wells A1–A12) to row B (wells B1–B12). Add 245 µL of Complete TR-FRET Coregulator Buffer K to each well in row B of the 96-well plate. Mix by pipetting up and down.

7. In order to assay 4 replicates of each ligand concentration, transfer 10 µL of each of the 2X ligand serial dilutions to the 384 well assay plate according to the plate layout in **Section 6.2**

*For example:* Using a 12 channel pipette transfer 10 µL aliquots from row B of the 96 well plate to alternate columns across row A of the 384 well plate, A1, A3, A5, A7, etc. Continue this process for Row A, this time pipetting into A2, A4, A6, etc. Repeat this process for row B of the 384-well plate.

Alternatively, use a 16 channel pipette where two tips fit into one well of the 96-well plate, then transfer to individual wells of the 384-well assay plate in rows A and B. Repeat to complete all of the columns (1–24).

With either method, wells A1, A2, B1, and B2 all contain the highest concentration of ligand from B1 of the 96-well plate. Wells A3, A4, B3, and B4 all contain the next concentration of ligand from B2 of the 96-well plate, etc for 4 replicates of each concentration with constant percent DMSO.

### Prepare 4X Fluormone™ ES2 Green

8. Prepare 4X Fluormone™ ES2 Green (12 nM) using Complete Nuclear Buffer K.

*For example:* The Fluormone™ ES2 Green has a stock concentration of 1800 nM and the recommended concentration for this kit is 3.0 nM, so prepare a 4X solution at 12 nM by adding 6.7 µL of Fluormone™ ES2 Green stock to 993 µL of Complete Nuclear Receptor Buffer K.

9. Add 5 µL of 4X Fluormone™ ES2 Green to rows A–D, columns 1–24 of the 384-well assay plate. We recommend adding from low to high concentration of ligand (right to left) to prevent ligand carry over.

**Note:** The 4X Fluormone™ ES2 Green solution exhibits some sticking to certain plastics so we recommend coating the pipet tips several times (~ 6 times) with the solution before transferring to the assay plate. This will decrease assay noise and improve Z' values.

**Prepare 4X ER Alpha-LBD (GST)/4X Tb anti-GST Antibody**

10. Prepare a solution containing 4X ER Alpha-LBD (GST) and 4X Tb anti-GST antibody (8 nM) using Complete Nuclear Receptor Buffer K. The recommended assay concentration of ER Alpha for this kit is listed on the Certificate of Analysis. *Never vortex the ER Alpha-LBD (GST) stock or dilutions.* Mix by pipetting or gentle inversion. Keep this solution on ice until needed for use in the assay. The stock concentration of Tb anti-GST antibody is indicated on both the vial label and the Certificate of Analysis (0.5 mg/mL = ~3.4 µM antibody)

*For example:* If the ER Alpha-LBD (GST) has a stock concentration of 4000 nM and the recommended concentration for this kit is 2.0 nM, prepare a 4X solution at 8 nM by adding 2 µL of ER Alpha-LBD (GST) stock and 2.4 µL of the Tb anti-GST antibody to 995 µL of Complete Nuclear Receptor Buffer K.

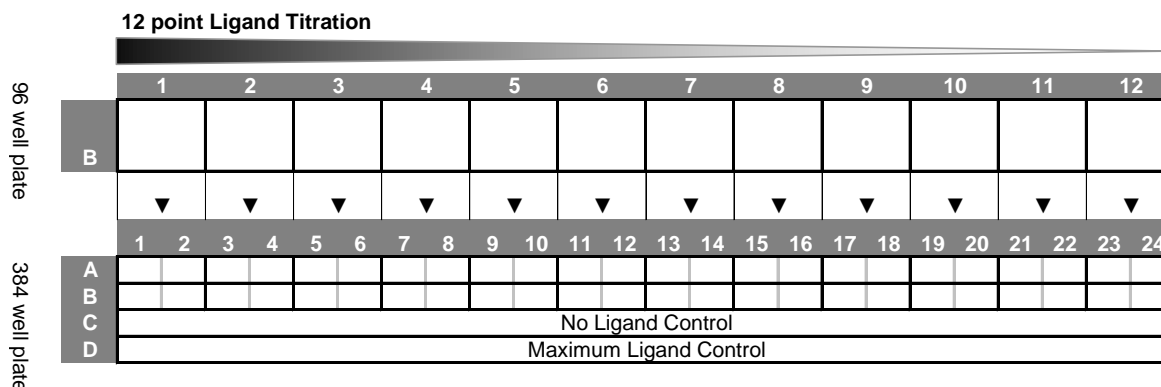
11. Add 5 µL of 4X ER Alpha-LBD (GST)/4X antibody solution to rows A–D, columns 1–24 of the 384-well assay plate (see plate layout in **Section 6.2**).

**Plate Incubation and Reads**

12. Briefly and gently mix the 384-well plate on a plate shaker and incubate at room temperature protected from light for the chosen equilibration time. The plate may be sealed with a cover to minimize evaporation.

**Note:** When using assay plates that are not coated such as Corning® 3677, drops may cling to the side of the wells during additions to the plate. It is critical that the plates be gently tapped or centrifuged to ensure that all reagents reach the bottom of the well.

13. Read the plate between 1 and 6 hours at your chosen equilibration time at wavelengths of 520 nm and 495 nm, using the instrument settings described in **Section 4**. If using a plate seal, spin the sealed plate in a centrifuge with an appropriate balance to spin down any condensation on the bottom of the seal, shake the plate gently, and read. Then proceed to data analysis as described in the next section.

**6.2 Plate Layout****6.3 Data Analysis**

Calculate the TR-FRET ratio by dividing the emission signal at 520 nm by the emission signal at 495 nm. Generate a binding curve by plotting the emission ratio vs. the log [ligand]. To determine the IC<sub>50</sub> value, fit the data using an equation for a sigmoidal dose response (variable slope), as provided by GraphPad Prism™ 4.0 or another comparable graphing program.

The “maximum ligand” and “no ligand” control data can be used to calculate Z'-factor based on the equation of Zhang *et al* (Zhang *et al.*, 1999).

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

## 6.4 Reagent Volumes and Controls

The following table summarizes the reagent volumes, order of addition, and recommended controls for the assay. The protocol in **Section 6.1** may be adapted to screen test compounds using the following summary table as a guide.

Assay Condition	Reagent Additions	Purpose
Test Compound	<ol style="list-style-type: none"> <li>1. 10 µL 2X Test Compound</li> <li>2. 5 µL 4X Fluormone™ ES2 Green</li> <li>3. 5 µL 4X ER Alpha-LBD (GST)/Tb-anti-GST Ab</li> </ol>	Assess competition by test compound of interest using a single concentration or dilution series.
Positive Control	<ol style="list-style-type: none"> <li>1. 10 µL 2X Control Competitor</li> <li>2. 5 µL 4X Fluormone™ ES2 Green</li> <li>3. 5 µL 4X ER Alpha-LBD (GST)/Tb-anti-GST Ab</li> </ol>	Represents 100% tracer displacement (minimum TR-FRET ratio) by a known ER Alpha ligand. We recommend using 1 µM estradiol as the Control Competitor.
Negative Control	<ol style="list-style-type: none"> <li>1. 10 µL 2X Test Compound Solvent</li> <li>2. 5 µL 4X Fluormone™ ES2 Green</li> <li>3. 5 µL 4X ER Alpha-LBD (GST)/Tb-anti-GST Ab</li> </ol>	Represents 0% tracer displacement (maximum TR-FRET ratio) and accounts for possible interference from a compound's vehicle solvent
No Receptor Control (Optional)	<ol style="list-style-type: none"> <li>1. 10 µL 2X Test Compound Solvent</li> <li>2. 5 µL 4X Fluormone™ ES2 Green</li> <li>3. 5 µL 4X Tb-anti-GST Ab</li> </ol>	Provides absolute minimum TR-FRET ratio possible, where only non-specific, diffusion-enhanced FRET between the antibody and Fluormone™ is detected.*

**Note:** All controls should contain the same percentage of solvent as the wells containing test compound.

\* Diffusion enhanced TR-FRET occurs when the donor (Tb-anti-GST Antibody) passes by the acceptor (ES2 Green) in solution during the excited state lifetime of the donor and is not related to a binding event. At higher concentrations of donor or acceptor, the probability of this occurrence increases and results in a larger background signal due to diffusion enhanced TR-FRET. The concentrations of ES2 Green used in this assay would result in minimal diffusion enhanced TR-FRET.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266


For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

## 6.5 Ligand Assay—Representative Data

**Figure 3. Representative experiment of LanthaScreen® TR-FRET ER Alpha Competitive Binding assay.** Assay composition: serial dilution of estradiol (1% final DMSO concentration), 2 nM ER Alpha-LBD (GST), 3 nM Fluormone™ ES2 Green, and 2 nM Tb anti-GST antibody. Results for 1-hour, 2-hour, 4-hour, and 6-hour incubations are shown with the corresponding  $IC_{50}$  values. The curves were generated using a sigmoidal dose response (variable slope) equation in IDBS XLfit® 4.2 software.

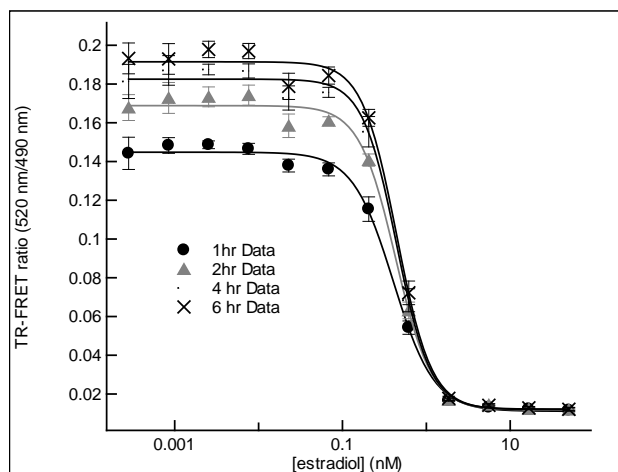

Time Point	$IC_{50}$ (nM)
1 hour	0.40
2 hour	0.43
4 hour	0.44
6 hour	0.44

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)


Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)


## 7. ASSAY PHARMACOLOGY

**Figure 4: Relative IC<sub>50</sub> Values of Selected Ligands for ER Alpha-LBD (GST) in the LanthaScreen® TR-FRET ER Alpha Competitive Binding Assay.** Serial dilutions of various test compounds (1% final DMSO concentration) were assayed (n = 4). Curves were fit using a sigmoidal dose-response equation (variable slope) in GraphPad Prism™ 4.0 software. Figures 4A and 4B are the same data plotted differently. The values plotted in 4B have been normalized by the minimum TR-FRET ratio obtained for estradiol so that the emission ratio for estradiol starts at 1.0. There is no effect on the Z' or calculated IC<sub>50</sub>'s.

A.

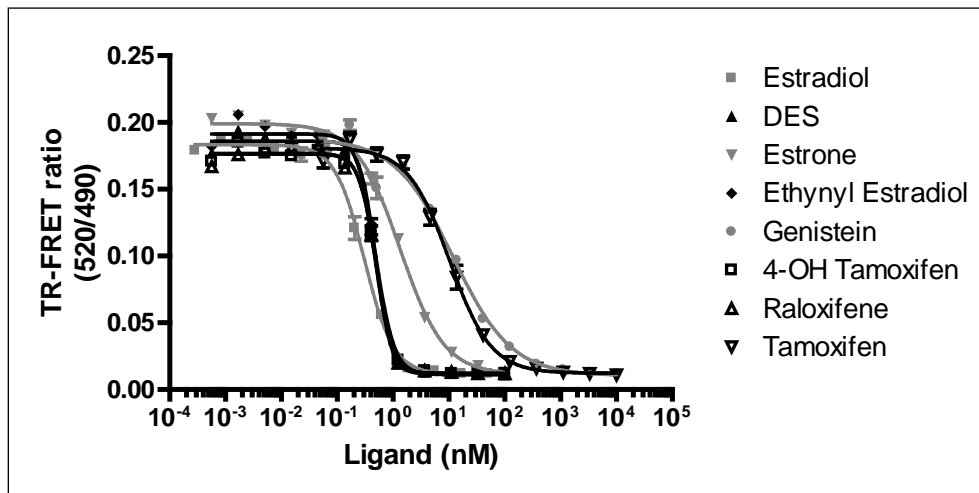

B.

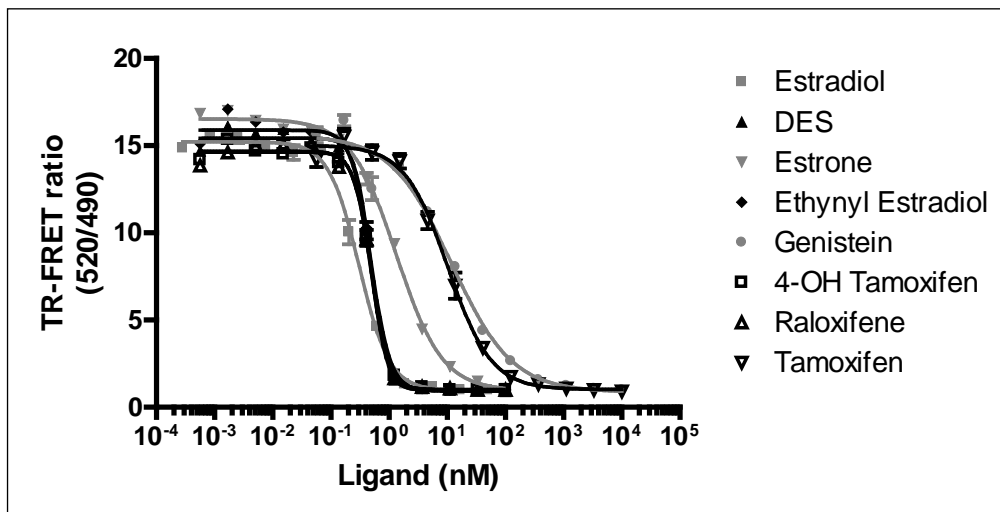

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)

## 8. REFERENCES

---

Zhang, J. H., Chung, T. D., and Oldenburg, K. R. (1999) A simple statistical parameter for use in evaluation and validation of high-throughput screening assays. *J. Biomol. Screen.*, 4, 67-73

## 9. LIMITED PRODUCT WARRANTY

---

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at [www.lifetechnologies.com/termsandconditions](http://www.lifetechnologies.com/termsandconditions). If you have any questions, please contact Life Technologies at [www.lifetechnologies.com/support](http://www.lifetechnologies.com/support).

## 10. PURCHASER NOTIFICATION

---

### Important Licensing Information

This product may be covered by one or more Limited Use Label Licenses. By use of this product, you accept the terms and conditions of all applicable Limited Use Label Licenses.

### DISCLAIMER

LIFE TECHNOLOGIES CORPORATION AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

### TRADEMARKS

All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.

GraphPad Prism is a trademark of GraphPad Software, Inc.

Corning is a registered trademark of Corning, Inc.

Tecan is a registered trademark of, and Infinite, GENious, Safire, and Safire2 are trademarks of Tecan Group AG.

BMG LABTECH and PHERAstar are registered trademarks of BMG LABTECH, GmbH.

Molecular Devices, Analyst, and SpectraMax are registered trademarks of, and Gemini is a trademark of Molecular Devices, LLC.

PerkinElmer and EnVision are registered trademarks of, and Victor is a trademark of, PerkinElmer, Inc.

BioTek is a registered trademark of, and Synergy is a trademark of, BTI Holdings, Inc.

LANCE is a registered trademark of Wallac Oy.

ViewLux is a registered trademark of PerkinElmer Singapore Pte Ltd.

XLfit is a registered trademark of ID Business Solutions Ltd.

© 2014 Thermo Fisher Scientific Inc. All rights reserved.

For Technical Support for this or other Drug Discovery Solutions Products, dial 760-603-7200 extension 40266

For information on frequently asked questions regarding the LanthaScreen® technology, please go to [www.lifetechnologies.com/lanthascreen](http://www.lifetechnologies.com/lanthascreen)

Thermo Fisher Scientific Inc. • 5791 Van Allen Way • Carlsbad, CA 92008 • Phone: 760 603 7200 • FAX: 760 602 6500

[www.lifetechnologies.com](http://www.lifetechnologies.com)