

ELISA for Mouse IL-17A

Product Code: 3521-1A-6

CONTENTS, development kit for 6 plates:

Vial 1 (yellow top)

Monoclonal antibody IL17-I (150 μ l)

Concentration: 0.5 mg/ml

Vial 2 (red top)

Biotinylated monoclonal antibody IL17-II (80 μ l)

Concentration: 0.5 mg/ml

Vial 3 (white top)

Streptavidin-Alkaline Phosphatase (80 μ l)

Vial 4

Recombinant mouse IL-17A standard (1 μ g)

To ensure total recovery of stated quantity, vials have been overfilled.

STORAGE:

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.


General

Intended use: For quantitative determination of native and recombinant mouse IL-17A in solution, e.g. cell culture supernatant.

Reagents: Antibodies are supplied in sterile-filtered (0.2 μm) PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.15% Kathon CG.

Standard range: 5-500 pg/ml

Standard calibration: No international standard exists for calibration.

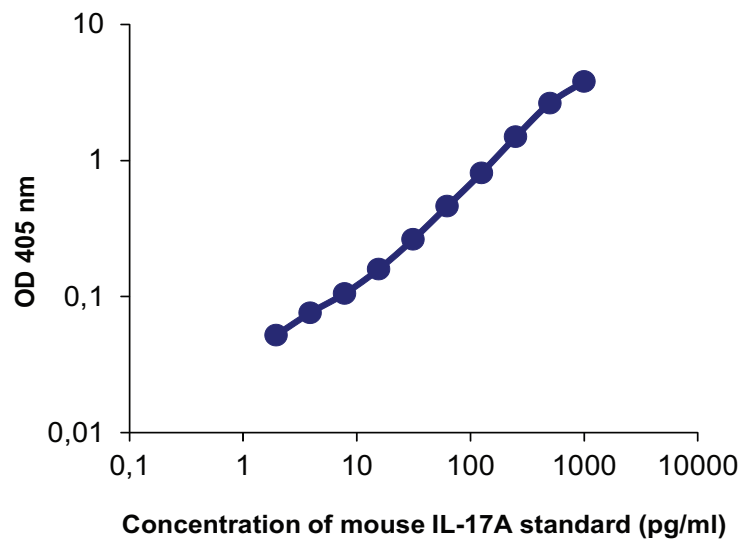

Guidelines for Mouse IL-17A ELISA

- Day 1**
1. Coat a high protein binding ELISA plate with mAb IL17-I, diluted to 1 $\mu\text{g/ml}$ in PBS, pH 7.4, by adding 100 $\mu\text{l/well}$. Incubate overnight at 4-8°C.
- Day 2**
2. Wash twice with PBS (200 $\mu\text{l/well}$).
 3. Block plate by adding 200 $\mu\text{l/well}$ of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
 4. Wash five times with PBS containing 0.05% Tween.
 5. Prepare mouse IL-17A standard by reconstituting contents of vial 4 in 1 ml PBS to make a stock solution of 1 $\mu\text{g/ml}$. Allow the standard to dissolve for 5 minutes, mix thoroughly and aliquot. Store at -20°C and avoid repeated freeze-thaw cycles of the standard aliquots. For the test, prepare dilutions of the stock using the standard range as a guideline.
 6. Add 100 $\mu\text{l/well}$ of samples or standards diluted in incubation buffer and incubate for 2 hours at room temperature.
 7. Wash as in step 4.
 8. Add 100 $\mu\text{l/well}$ of mAb IL17-II-biotin at 0.5 $\mu\text{g/ml}$ in incubation buffer or ELISA diluent for serum/plasma samples. Incubate for 1 hour at room temperature.
 9. Wash as in step 4.
 10. Add 100 $\mu\text{l/well}$ of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
 11. Wash as in step 4.
 12. Add 100 $\mu\text{l/well}$ of appropriate substrate solution e.g. p-nitrophenyl-phosphate (pNPP).
 13. Measure the optical density (405 nm for pNPP) in an ELISA reader after suitable developing time.

NOTE; for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



MABTECH AB
Box 1233
SE-131 28 Nacka Strand
Sweden
Tel: +46 8 716 27 00
Fax: +46 8 716 27 01
E-mail: mabtech@mabtech.com
www.mabtech.com

MABTECH Inc
M.E.B. 220
3814 West Street
Cincinnati, OH 45227
USA
Tel: +1 513 871 4500
Fax: +1 513 871 7353
E-mail: mabtech.usa@mabtech.com

MABTECH AB Büro Deutschland
Germany
Tel: +49 40 4135 7935
Fax: +49 40 4135 7945
E-mail: mabtech.de@mabtech.com

2013-01-24

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the following standards:



MABTECH AUSTRALIA Pty Ltd
resolvingIMAGES
Unit 22, 196 Settlement Road
Thomastown Victoria 3074
Australia
Tel: +61 3 9466 4007
Fax: +61 3 9466 4003
E-mail: mabtech.au@mabtech.com

MABTECH AB Bureau de liaison France
BP 255, 1300 route des Crêtes
06905 Sophia Antipolis
France
Tel: +33 (0)4 92 38 80 70
Fax: +33 (0)4 92 38 80 71
E-mail: mabtech.fr@mabtech.com