


Anti-Human CD8a APC

Catalog Number: 17-0086

Also Known As: CD8 alpha, leu-2a

RUO: For Research Use Only

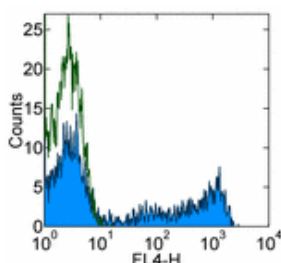

Staining of normal human peripheral blood cells with Mouse IgG2a κ Isotype Control APC (cat. 17-4724) (open histogram) or Anti-Human CD8a APC (filled histogram). Cells in the lymphocyte gate were used for analysis.

Product Information

Contents: Anti-Human CD8a APC


REF Catalog Number: 17-0086

Clone: OKT8 (OKT-8)

Concentration: 5 µl (0.06 µg)/test


Host/Isotype: Mouse IgG2a

Formulation: Phosphate buffer pH 7.2,
150 mM NaCl, 0.09% NaN₃, 0.2% BSA

 Temperature Limitation: Store at 2-8°C. Do not freeze. Light sensitive material.

LOT Batch Code: Refer to Vial

 Use By: Refer to Vial

 Caution, contains Azide

Description

The OKT8 monoclonal antibody reacts with the human CD8a molecule, an approximately 32-34 kDa cell surface receptor expressed either as a heterodimer with the CD8 β chain (CD8 αβ) or as a homodimer (CD8 αα). A majority of thymocytes and a subpopulation of mature T cells and NK cells express CD8a. CD8 binds to MHC class I and through its association with protein tyrosine kinase p56lck plays a role in T-cell development and activation of mature T cells. Preliminary testing indicates that OKT8 and two other mouse anti-human CD8 antibodies (clone RPA-T8, Cat. No.14-0088 and clone HIT8a, Cat. No.14-0089) do not compete with each other for binding to human peripheral blood leukocytes by flow cytometric analysis, suggesting that they do not bind to similar epitopes or block each other by steric hindrance.

Applications Reported

For research use only, not for diagnostic or therapeutic use. This OKT8 (OKT-8) antibody has been reported for use in flow cytometric analysis.

Applications Tested

This OKT8 (OKT-8) antibody has been pre-titrated and tested by flow cytometric analysis of human peripheral blood leukocytes. This can be used at 5 µl (0.06 µg)/per test. A test is defined as the amount (µg)/test of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

References

Kay HD, Horwitz DA. 1980. Evidence by reactivity with hybridoma antibodies for a probable myeloid origin of peripheral blood cells active in natural cytotoxicity and antibody-dependent cell-mediated cytotoxicity. *J Clin Invest.* 66(4):847-51.

Sayos J, Wu C, Morra M, Wang N, Zhang X, Allen D, van Schaik S, Notarangelo L, Geha R, Roncarolo MG, Oettgen H, De Vries JE, Aversa G, Terhorst C. The X-linked lymphoproliferative-disease gene product SAP regulates signals induced through the co-receptor SLAM. *Nature.* 1998 Oct 1;395(6701):462-9. (OKT8, IP, PubMed)

Thomas Y, Sosman J, Irigoyen O, Friedman SM, Kung PC, Goldstein G, Chess L. 1980. Functional analysis of human T cell subsets defined by monoclonal antibodies. I. Collaborative T-T interactions in the immunoregulation of B cell differentiation. *J Immunol.* 125(6):2402-8.

Campanelli R, Palermo B, Garbelli S, Mantovani S, Lucchi P, Necker A, Lantelme E, Giachino C. 2002. Human CD8 co-receptor is strictly involved in MHC-peptide tetramer-TCR binding and T cell activation. *Int Immunol.* 14(1):39-44.

Related Products

17-4724 Mouse IgG2a K Isotype Control APC

Not for further distribution without written consent.

Copyright © 2000-2010 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.eBioscience.com • info@eBioscience.com