


## Anti-Mouse CD24 PerCP-eFluor® 710

**Catalog Number:** 46-0242

**Also Known As:** Heat Stable Antigen, HSA

**RUO: For Research Use Only. Not for use in diagnostic procedures.**

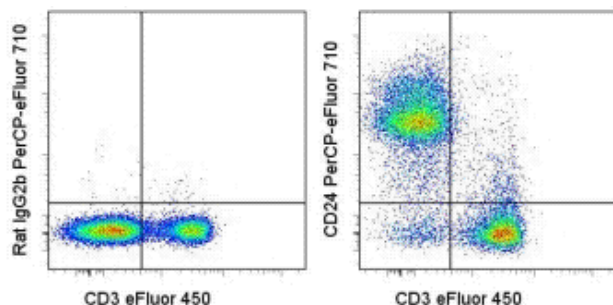

Staining of C57Bl/6 splenocytes with Anti-Mouse CD3 eFluor® 450 (cat. 48-0032) and 0.06 ug of Rat IgG2b K Isotype Control PerCP-eFluor® 710 (cat. 46-4031) (left) or 0.06 ug of Anti-Mouse CD24 PerCP-eFluor® 710 (right). Cells in the lymphocyte gate were used for analysis.

### Product Information

**Contents:** Anti-Mouse CD24 PerCP-eFluor® 710

**REF** **Catalog Number:** 46-0242

**Clone:** M1/69

**Concentration:** 0.2 mg/mL

**Host/Isotype:** Rat IgG2b, kappa

**Formulation:** aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

**Temperature Limitation:** Store at 2-8°C. Do not freeze. Light sensitive material.

**LOT** **Batch Code:** Refer to Vial

**Use By:** Refer to Vial

**Contains sodium azide**

### Description

The M1/69 monoclonal antibody reacts with the mouse CD24 molecule, also known as Heat Stable Antigen (HSA). This 35-50 kDa molecule is anchored in the plasma membrane via phosphatidylinositol and is expressed by erythrocytes, thymocytes, peripheral lymphocytes and myeloid lineage. CD24 is a variably glycosylated molecule resulting in heterogeneity of molecular mass of this antigen on cells of different lineages and antibodies to CD24 exhibit subtle differences in staining level on lymphocyte populations. The expression of CD24 has been used to resolve stages of B lymphopoiesis in mouse bone marrow. It has been reported that P-selectin (CD62P) binds to CD24.

### Applications Reported

This M1/69 antibody has been reported for use in flow cytometric analysis.

### Applications Tested

This M1/69 antibody has been tested by flow cytometric analysis of mouse spleen cells. This can be used at less than or equal to 0.125 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

**PerCP-eFluor® 710 can be used in place of PE-Cy5, PE-Cy5.5 or PerCP-Cy5.5. PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm). Please make sure that your instrument is capable of detecting this fluorochrome. For a filter configuration, we recommend using the 685 LP dichroic mirror and 710/40 band pass filter, however the 695/40 band pass filter is an acceptable alternative.**

**Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 uL cell sample + 100 uL IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.**

### References

Hunte BE, Capone M, Zlotnik A, Rennick D, Moore TA. 1998. Acquisition of CD24 expression by Lin-CD43+B220(low)ckit(hi) cells coincides with commitment to the B cell lineage. *Eur J Immunol.* 28(11):3850-6.

Wilson, A., L. M. Day, et al. 1988. Subpopulations of mature murine thymocytes: properties of CD4-CD8+ and CD4+CD8- thymocytes lacking the heat-stable antigen. *Cell Immunol* 117(2): 312-26.

Alterman, L. A., I. N. Crispe, et al. 1990. Characterization of the murine heat-stable antigen: an hematolymphoid differentiation antigen defined by the J11d, M1/69 and B2A2 antibodies. *Eur J Immunol* 20(7): 1597-602.

Springer T, Galfre G, Secher DS, Milstein C. 1978. Monoclonal xenogeneic antibodies to murine cell surface antigens: identification of novel leukocyte differentiation antigens. Eur J Immunol. 8(8):539-51.

**Related Products**

17-0031 Anti-Mouse CD3e APC (145-2C11)

46-4031 Rat IgG2b K Isotype Control PerCP-eFluor® 710

---

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • [www.eBioscience.com](http://www.eBioscience.com) • [info@eBioscience.com](mailto:info@eBioscience.com)