


Anti-Human CD62L (L-Selectin) PerCP-eFluor® 710

Catalog Number: 46-0629

Also Known As: LECAM-1, Leu-8

RUO: For Research Use Only. Not for use in diagnostic procedures.

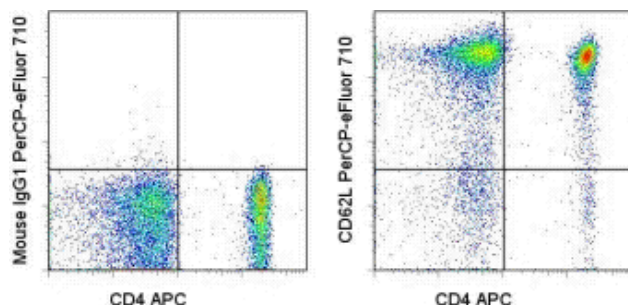

Staining of normal human peripheral blood cells with Anti-Human CD4 APC (cat. 17-0048) and Mouse IgG1 K Isotype Control PerCP-eFluor® 710 (cat. 46-4714) (left) or Anti-Human CD62L (L-Selectin) PerCP-eFluor® 710 (right). Cells in the lymphocyte gate were used for analysis.

Product Information

Contents: Anti-Human CD62L (L-Selectin) PerCP-eFluor® 710

REF **Catalog Number:** 46-0629

Clone: DREG-56 (DREG56)

Concentration: 5 µL (0.125 µg)/test

Host/Isotype: Mouse IgG1, kappa

HLDA Workshop: V S056

Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

Temperature Limitation: Store at 2-8°C. Do not freeze. Light sensitive material.

LOT **Batch Code:** Refer to Vial

Use By: Refer to Vial

Contains sodium azide

Description

The DREG-56 monoclonal antibody reacts with human CD62L, a 76 kDa member of the selectin family. CD62L is expressed by neutrophils, monocytes, and subsets of T, B, and NK cells and binds a number of glycosylated, fucosylated, sulfated sialylated glycoproteins including CD34, glycamin-1 and MAdCAM-1. These interactions mediate rolling of lymphocytes on activated endothelium at the sites of inflammation and homing of cells to the high endothelial venules (HEV) of peripheral lymphoid tissues.

Applications Reported

This DREG-56 (DREG56) antibody has been reported for use in flow cytometric analysis.

Applications Tested

This DREG-56 (DREG56) antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.125 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

PerCP-eFluor® 710 can be used in place of PE-Cy5, PE-Cy5.5 or PerCP-Cy5.5. PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm). Please make sure that your instrument is capable of detecting this fluorochrome. For a filter configuration, we recommend using the 685 LP dichroic mirror and 710/40 band pass filter, however the 695/40 band pass filter is an acceptable alternative.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

References

Jutila MA, Kurk S, Jackiw L, Knibbs RN, Stoolman LM. J Immunol. L-selectin serves as an E-selectin ligand on cultured human T lymphoblasts. 2002 Aug 15;169(4):1768-73. (**Dreg-56**, WB, PubMed)

Schlossman, S., L. Bloumsell, et al. eds. 1995. Leucocyte Typing V: White Cell Differentiation Antigens. Oxford University Press. New York.

Kishimoto TK, Jutila MA, Butcher EC. Identification of a human peripheral lymph node homing receptor: a rapidly down-regulated adhesion molecule. Proc Natl Acad Sci U S A. 1990 Mar;87(6):2244-8. (**Dreg-56**, mAb generation, WB, FC, PubMed)

Related Products

11-0048 Anti-Human CD4 FITC (OKT4 (OKT-4))

17-0048 Anti-Human CD4 APC (OKT4 (OKT-4))

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.eBioscience.com • info@eBioscience.com