


An Affymetrix Company

# Anti-Mouse CD199 (CCR9) PerCP-eFluor® 710

Catalog Number: 46-1991 Also known as: CCR-9

RUO: For Research Use Only. Not for use in diagnostic procedures.

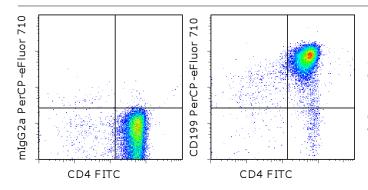

Staining of C57Bl/6 thymocytes with Anti-Mouse CD4 APC (cat. 17-0041) and 0.25 ug of Mouse IgG2a kappa Isotype Control PerCP-eFluor® 710 (cat. 46-4724) (left) or 0.25 ug of Anti-Mouse CD199 (CCR9) PerCP-eFluor® 710 (right). Total viable cells were used for analysis.

## **Product Information**

**Contents:** Anti-Mouse CD199 (CCR9) PerCP-eFluor® 710

Clone: eBioCW-1.2 (CW-1.2)
Concentration: 0.2 mg/mL
Host/Isotype: Mouse IgG2a



**Formulation:** aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer **Temperature Limitation:** Store at 2-8°C. Do not

freeze. Light-sensitive material. **Batch Code:** Refer to vial **Use By:** Refer to vial



The eBioCW-1.2 monoclonal antibody reacts with mouse CCR9 (CD199), which is the receptor for thymus-expressed chemokine (TECK). CCR9 is a member of the G protein coupled receptor (GPCR) supergene family, and is involved in trafficking of T cell progenitors within the thymus. CCR9 expression during thymocyte development commences at the double-negative (DN) 3 stage (CD4-CD8-CD25+CD44-), peaks in the double-positive (DP) stage (CD4+CD8+CD25-CD44-), and is down-regulated in committed CD4+ or CD8+ single-positive (SP) thymoyctes. CCR9-deficient mice show a mild impairment in thymoycte development. In the periphery, CCR9 is thought to be expressed on naïve CD8+ T cells, but not on naïve CD4+ T cells.

### **Applications Reported**

This eBioCW-1.2 (CW-1.2) antibody has been reported for use in flow cytometric analysis.

## **Applications Tested**

This eBioCW-1.2 (CW-1.2) antibody has been tested by flow cytometric analysis of mouse thymocytes. This can be used at less than or equal to 0.5  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

PerCP-eFluor® 710 can be used in place of PE-Cy5, PE-Cy5.5 or PerCP-Cy5.5. PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm). Please make sure that your instrument is capable of detecting this fluorochrome. For a filter configuration, we recommend using the 685 LP dichroic mirror and 710/40 band pass filter, however the 695/40 band pass filter is an acceptable alternative.

Our testing indicates that PerCP-eFluor® 710 conjugated antibodies are stable when stained samples are exposed to freshly prepared 2% formaldehyde overnight at 4°C, but please evaluate for alternative fixation protocols.

Click here or contact eBioscience Technical Support for more information on eFluor™ Organic Dyes including PerCP-eFluor® 710.



An Affymetrix Company

# Anti-Mouse CD199 (CCR9) PerCP-eFluor® 710

Catalog Number: 46-1991 Also known as: CCR-9

RUO: For Research Use Only. Not for use in diagnostic procedures.

#### References

Liu C, Saito F, Liu Z, Lei Y, Uehara S, Love P, Lipp M, Kondo S, Manley N, Takahama Y. Coordination between CCR7- and CCR9-mediated chemokine signals in pre-vascular fetal thymus colonization. Blood. 2006 Jun 29.

Wurbel MA, Malissen B, Campbell JJ. Complex regulation of CCR9 at multiple discrete stages of T cell development. Eur J Immunol. 2006 Jan;36(1):73-81. (CW-1.2, FC, Development of mAb, PubMed)

Zaballos A, Gutierrez J, Varona R, Ardavin C, Marquez G. Cutting edge: identification of the orphan chemokine receptor GPR-9-6 as CCR9, the receptor for the chemokine TECK. J Immunol. 1999 May 15;162(10):5671-5.

### **Related Products**

17-0041 Anti-Mouse CD4 APC (GK1.5) 46-4724 Mouse IgG2a K Isotype Control PerCP-eFluor® 710