


Anti-Human CD274 (B7-H1) PerCP-eFluor® 710

Catalog Number: 46-5983

Also known as: PD-L1

RUO: For Research Use Only. Not for use in diagnostic procedures.

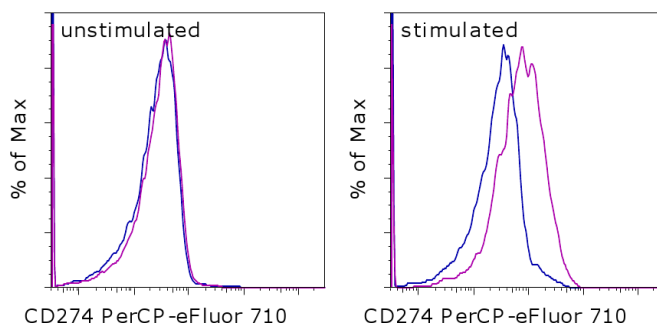

Staining of unstimulated (left) or PHA-stimulated (right) normal human peripheral blood cells with Mouse IgG1 K Isotype Control PerCP-eFluor® 710 (cat. 46-4714) (blue histogram) or Anti-Human CD274 (B7-H1) PerCP-eFluor® 710 (purple histogram). Total viable cells were used for analysis.

Product Information

Contents: Anti-Human CD274 (B7-H1)
PerCP-eFluor® 710



Catalog Number: 46-5983

Clone: MIH1

Concentration: 5 µL (1 µg)/test

Host/Isotype: Mouse IgG1



Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

Temperature Limitation: Store at 2-8°C. Do not freeze. Light-sensitive material.



Batch Code: Refer to vial



Use By: Refer to vial

Description

The MIH1 monoclonal antibody reacts with human B7-H1, also known as PD-L1. B7-H1, a member of the B7 family, has a predicted molecular weight of approximately 40 kDa and belongs to the Ig superfamily. B7-H1 is expressed on a majority of leukocytes. B7-H1 is a ligand for PD-1. Interaction of PD-1 with either PD-L1 (B7-H1) or PD-L2 (B7-DC) results in inhibition of T and B cell responses. MIH1 is reported to be a blocking antibody.

Applications Reported

This MIH1 antibody has been reported for use in flow cytometric analysis.

Applications Tested

This MIH1 antibody has been pre-titrated and tested by flow cytometric analysis of stimulated human peripheral blood cells. This can be used at 5 µL (1 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

PerCP-eFluor® 710 can be used in place of PE-Cy5, PE-Cy5.5 or PerCP-Cy5.5. PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm). Please make sure that your instrument is capable of detecting this fluorochrome. For a filter configuration, we recommend using the 685 LP dichroic mirror and 710/40 band pass filter, however the 695/40 band pass filter is an acceptable alternative.

Our testing indicates that PerCP-eFluor® 710 conjugated antibodies are stable when stained samples are exposed to freshly prepared 2% formaldehyde overnight at 4°C, but please evaluate for alternative fixation protocols.

Click here or contact eBioscience Technical Support for more information on eFluor™ Organic Dyes including PerCP-eFluor® 710.

References

Chen Y, Zhang J, Li J, Zou L, Zhao T, Tang Y, Wu Y. Expression of B7-H1 in inflammatory renal tubular epithelial

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.ebioscience.com •
info@ebioscience.com

Anti-Human CD274 (B7-H1) PerCP-eFluor® 710

Catalog Number: 46-5983

Also known as: PD-L1

RUO: For Research Use Only. Not for use in diagnostic procedures.

cells.Nephron Exp Nephrol. 2006; 102(3-4):e81-92. (**MIH1**, IH/F, FA, PubMed)

Saudemont, A., Jouy, N., Hetuin, D., and Quesnel, B. 2005. NK cells that are activated by CXCL10 can kill dormant tumor cells that resist CTL-mediated lysis and can express B7-H1 that stimulates T cells. *Blood*. 105:2428-2435.

Ohigashi Y, Sho M, Yamada Y, Tsurui Y, Hamada K, Ikeda N, Mizuno T, Yoriki R, Kashizuka H, Yane K, Tsushima F, Otsuki N, Yagita H, Azuma M, Nakajima Y. Clinical significance of programmed death-1 ligand-1 and programmed death-1 ligand-2 expression in human esophageal cancer. *Clin Cancer Res*. 2005 Apr 15;11(8):2947-53. (**MIH1**, IH/F, PubMed)

Kanai T, Totsuka T, Uraushihara K, Makita S, Nakamura T, Koganei K, Fukushima T, Akiba H, Yagita H, Okumura K, Machida U, Iwai H, Azuma M, Chen L, Watanabe M. 2003. Blockade of B7-H1 suppresses the development of chronic intestinal inflammation. *J Immunol*. 171(8):4156-63.

Freeman, G. J., A. J. Long, et al. 2000. Engagement of the PD-1 immunoinhibitory receptor by a novel B7 family member leads to negative regulation of lymphocyte activation. *J Exp Med* 192(7): 1027-34.

Dong, H., G. Zhu, et al. 1999. B7-H1, a third member of the B7 family, co-stimulates T-cell proliferation and interleukin-10 secretion. *Nat Med* 5(12): 1365-9.

Related Products

46-4714 Mouse IgG1 K Isotype Control PerCP-eFluor® 710 (P3.6.2.8.1)

65-0865 Fixable Viability Dye eFluor® 780