


Anti-Mouse/Rat CD29 (Integrin beta 1) APC-eFluor® 780

Catalog Number: 47-0291

Also known as: Fibronectin Receptor Beta

RUO: For Research Use Only. Not for use in diagnostic procedures.

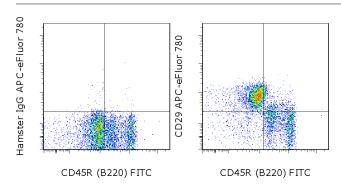

Staining of C57BI/6 bone marrow cells with Anti-Human/Mouse CD45R (B220) FITC (cat. 11-0452) and 0.25 ug of Armenian Hamster IgG Isotype Control APC-eFluor® 780 (cat. 47-4888) (left) or 0.25 ug of Anti-Mouse/Rat CD29 (Integrin beta 1) APC-eFluor® 780 (right). Total viable cells were used for analysis.

Product Information

Contents: Anti-Mouse/Rat CD29 (Integrin

beta 1) APC-eFluor® 780 REF Catalog Number: 47-0291 Clone: eBioHMb1-1 (HMb1-1) Concentration: 0.2 mg/mL

Host/Isotype: Armenian Hamster IgG



Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer Temperature Limitation: Store at 2-8°C. Do not freeze. Light-sensitive material. This tandem dye is sensitive to photo-induced oxidation. Protect this vial from light during storage, handling &



Batch Code: Refer to vial Use By: Refer to vial Contains sodium azide

experimental procedures.



Description

The eBioHMb1-1 monoclonal antibody reacts with mouse and rat CD29 (integrin beta 1), a 110-120 kDa member of the beta integrin family expressed by leukocytes, endothelial, smooth muscle and epithelial cells. CD29 binds noncovalently with the alpha integrins CD49a-f to form the VLA-1 through VLA-6 complexes, as well as with CD51. These alpha-beta integrin heterodimers are capable of mediating a variety of cellular responses including adhesion, trafficking, proliferation and differentiation. All integrins which include CD29 bind to extracellular matrix proteins including collagen, laminin, fibronectin and vitronectin, whereas some CD29-containing integrins can also interact with cellular receptors such as VCAM-1 and MadCAM-1.

Applications Reported

This eBioHMb1-1 (HMb1-1) antibody has been reported for use in flow cytometric analysis.

Applications Tested

This eBioHMb1-1 (HMb1-1) antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This can be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 μL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

APC-eFluor® emits at 780 nm and is excited with the Red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Light sensitivity: Tandem is sensitive photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 uL cell sample + 100 uL IC Fixation



Anti-Mouse/Rat CD29 (Integrin beta 1) APC-eFluor® 780

Catalog Number: 47-0291

Also known as: Fibronectin Receptor Beta

RUO: For Research Use Only. Not for use in diagnostic procedures.

Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

References

Noto K, Kato K, Okumura K, Yagita H. Identification and functional characterization of mouse CD29 with a mAb. Int Immunol. 1995 May;7(5):835-42. (**HMbeta1-1**, mAb development, IP, FA, PubMed)

Ridger VC, Wagner BE, Wallace WA, Hellewell PG. Differential effects of CD18, CD29, and CD49 integrin subunit inhibition on neutrophil migration in pulmonary inflammation. J Immunol. 2001 Mar 1;166(5):3484-90. (**HMbeta1-1**, FA, PubMed)

Related Products

11-0452 Anti-Human/Mouse CD45R (B220) FITC (RA3-6B2) 47-4888 Armenian Hamster IgG Isotype Control APC-eFluor® 780 (eBio299Arm)