


## Anti-Human CD56 (NCAM) APC-eFluor<sup>®</sup> 780

**Catalog Number:** 47-0567

**Also known as:** Neural cell adhesion molecule

**RUO: For Research Use Only. Not for use in diagnostic procedures.**

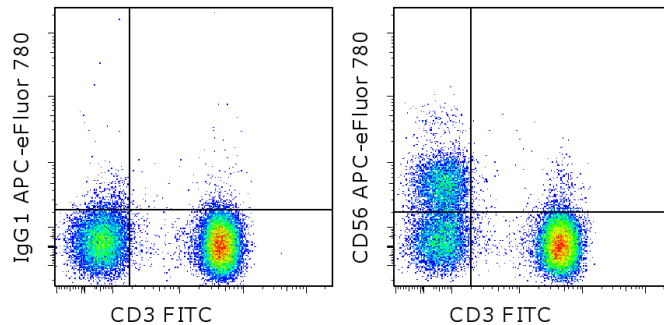

Staining of normal human peripheral blood cells with Anti-Human CD3 FITC (cat. 11-0038) and Mouse IgG1 K Isotype Control APC-eFluor<sup>®</sup> 780 (cat. 47-4714) (left) or Anti-Human CD56 (NCAM) APC-eFluor<sup>®</sup> 780 (right). Cells in the lymphocyte gate were used for analysis.

### Product Information

**Contents:** Anti-Human CD56 (NCAM) APC-eFluor<sup>®</sup> 780

**REF** **Catalog Number:** 47-0567

**Clone:** CMSSB

**Concentration:** 5  $\mu$ L (0.06  $\mu$ g)/test

**Host/Isotype:** Mouse IgG1, kappa

**Formulation:** aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

**Temperature Limitation:** Store at 2-8°C. Do not freeze. Light-sensitive material. This tandem dye is sensitive to photo-induced oxidation. Protect this vial from light during storage, handling & experimental procedures.

**Batch Code:** Refer to vial

**Use By:** Refer to vial

**Contains sodium azide**

### Description

This CMSSB monoclonal antibody reacts with human CD56, also known as Neural Cell Adhesion Molecule (NCAM). CD56 is a highly glycosylated transmembrane molecule expressed by neurons and plays a role in the homotypic adhesion of neural cells. In the hematopoietic system, CD56 is expressed on NK cells and a subset of T cells referred to as NKT cells. Staining with CMSSB does not block binding of MEM188 or CB56.

### Applications Reported

This CMSSB antibody has been reported for use in flow cytometric analysis.

### Applications Tested

This CMSSB antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5  $\mu$ L (0.06  $\mu$ g) per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from  $10^5$  to  $10^8$  cells/test.

**APC-eFluor<sup>®</sup> emits at 780 nm and is excited with the Red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochrome.**

**Light sensitivity:** Tandem is sensitive photo-induced oxidation. Please protect this vial and stained samples from light.

**Fixation:** Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100  $\mu$ L cell sample + 100  $\mu$ L IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • [www.ebioscience.com](http://www.ebioscience.com) •  
[info@ebioscience.com](mailto:info@ebioscience.com)

---

## Anti-Human CD56 (NCAM) APC-eFluor® 780

**Catalog Number:** 47-0567

**Also known as:** Neural cell adhesion molecule

**RUO: For Research Use Only. Not for use in diagnostic procedures.**

---

### References

Sonnenberg GF, Monticelli LA, Alenghat T, Fung TC, Hutnick NA, Kunisawa J, Shibata N, Grunberg S, Sinha R, Zahm AM, Tardif MR, Sathaliyawala T, Kubota M, Farber DL, Collman RG, Shaked A, Fouser LA, Weiner DB, Tessier PA, Friedman JR, Kiyono H, Bushman FD, Chang KM, Artis D. Innate lymphoid cells promote anatomical containment of lymphoid-resident commensal bacteria. *Science*. 2012 Jun 8;336(6086):1321-5. (CMSSB, FC, PubMed)

Hayakawa Y, Huntington ND, Nutt SL, Smyth MJ. Functional subsets of mouse natural killer cells. *Immunol Rev*. 2006 Dec;214:47-55.

Kishimoto, T., A.E.G., von dem Borne, et al. eds. 1998. *Leucocyte Typing VI: White Cell Differentiation Antigens*. Garland Publishing, Inc. London.

### Related Products

11-0038 Anti-Human CD3 FITC (UCHT1)

12-0577 Anti-Human CD57 PE (TB01 (TBO1))

46-0087 Anti-Human CD8a PerCP-eFluor® 710 (SK1 (SK-1))

47-4714 Mouse IgG1 K Isotype Control APC-eFluor® 780 (P3.6.2.8.1)

---

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • [www.ebioscience.com](http://www.ebioscience.com) •  
[info@ebioscience.com](mailto:info@ebioscience.com)