


Anti-Human CD14 eFluor® 450

Catalog Number: 48-0149

RUO: For Research Use Only. Not for use in diagnostic procedures.

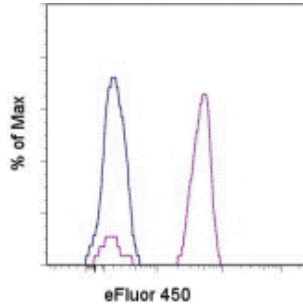

Staining of normal human peripheral blood cells with Mouse IgG1 K Isotype Control eFluor® 450 (cat. 48-4714) (blue histogram) or Anti-Human CD14 eFluor® 450 (purple histogram). Cells in the monocyte gate were used for analysis.

Product Information

Contents: Anti-Human CD14 eFluor® 450

REF **Catalog Number:** 48-0149


Clone: 61D3

Concentration: 5 µL (0.5 µg)/test

Host/Isotype: Mouse IgG1, kappa


HLDA Workshop: V MA085

Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

 **Temperature Limitation:** Store at 2-8°C. Do not freeze. Light sensitive material.

 **Batch Code:** Refer to Vial

 **Use By:** Refer to Vial

 **Contains sodium azide**

Description

The 61D3 monoclonal antibody reacts with human CD14, a 53-55 kDa GPI-linked glycoprotein. CD14 is expressed on monocytes, interfollicular macrophages and some dendritic cells. Complexes of LPS and LBP (LPS-Binding Protein) bind with high affinity to monocytes through the surface CD14.

Applications Reported

This 61D3 antibody has been reported for use in flow cytometric analysis.

eFluor® 450 is a replacement for Pacific Blue®. eFluor® 450 emits at 456 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochoime.

Applications Tested

This 61D3 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.5 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

eFluor® 450 is a replacement for Pacific Blue®. eFluor® 450 emits at 456 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochoime.

References

Fadok VA, Warner ML, Bratton DL, Henson PM. CD36 is required for phagocytosis of apoptotic cells by human macrophages that use either a phosphatidylserine receptor or the vitronectin receptor (alpha v beta 3). J Immunol 1998 Dec 1;161(11):6250-7.

Kishimoto, T., A.E.G., von dem Borne, et al. eds. 1998 Leucocyte Typing VI: White Cell Differentiation Antigens. Garland Publishing Inc. London.

Schlossman, S., L. Bloumsell, et al. eds 1995. Leucocyte Typing V: White Cell Differentiation Antigens. Oxford University Press. New York.

Related Products

48-4714 Mouse IgG1 K Isotype Control eFluor® 450 (P3.6.2.8.1)

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.eBioscience.com • info@eBioscience.com