Applications	Reactivity	Sensitivity	MW (kDa)	Source
W	H M R Mk	Endogenous	53	Rabbit

Applications Key: W=Western Blotting

Reactivity Key: H=Human M=Mouse R=Rat Mk=Monkey

Species cross-reactivity is determined by western blot. Species enclosed in parentheses are predicted to react based on 100% sequence homology.

Protocols

Western Immunoblotting Protocol (Primary Ab Incubation In Milk)

For Western blots, incubate membrane with diluted antibody in 5% w/v nonfat dry milk, 1X TBS, 0.1% Tween-20 at 4 $^{\circ}$ C with gentle shaking, overnight.

Products available from Cell Signaling Technology are linked by their respective catalog numbers.

A. Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1. 1X Phosphate Buffered Saline (PBS).
- 2. **1X SDS Sample Buffer:** (#7722, #7723) 62.5 mM Tris-HCl (pH 6.8 at 25 °C), 2% w/v SDS, 10% glycerol, 50 mM DTT, 0.01% w/v bromophenol blue or phenol red.
- 3. **Transfer Buffer:** 25 mM Tris base, 0.2 M glycine, 20% methanol (pH 8.5).
- 4. **10X Tris Buffered Saline (TBS):** (#9997) To prepare 1 liter of 10X TBS: 24.2 g Tris base, 80 g NaCl; adjust pH to 7.6 with HCl (use at 1X).
- 5. **Nonfat Dry Milk:** (#9999) (weight to volume [w/v]).
- 6. **Blocking Buffer:** 1X TBS, 0.1% Tween-20 with 5% w/v nonfat dry milk; for 150 ml, add 15 ml 10X TBS to 135 ml water, mix. Add 7.5 g nonfat dry milk and mix well. While stirring, add 0.15 ml Tween-20 (100%).
- 7. Wash Buffer: 1X TBS, 0.1% Tween-20 (TBS/T).
- 8. **Primary Antibody Dilution Buffer:** 1X TBS, 0.1% Tween-20 with 5% nonfat dry milk; for 20 ml, add 2 ml 10X TBS to 18 ml water, mix. Add 1.0 g nonfat dry milk and mix well. While stirring, add 20 μl Tween-20 (100%).
- 9. Phototope®-HRP Western Blot Detection System: (#7071 anti-rabbit) or (#7072 anti-mouse) Includes biotinylated protein ladder, secondary (#7074 anti-rabbit) or (#7076 anti-mouse) antibody conjugated to horseradish peroxidase (HRP), anti-biotin antibody conjugated to HRP, LumiGLO® chemiluminescent reagent and peroxide.
- 10. Prestained Protein Marker, Broad Range (Premixed Format): (#7720).
- 11. Biotinylated Protein Ladder Detection Pack: (#7727).
- Blotting Membrane: This protocol has been optimized for nitrocellulose membranes, which CST recommends. PVDF
 membranes may also be used.

B. Protein Blotting

A general protocol for sample preparation is described below.

- 1. Treat cells by adding fresh media containing regulator for desired time.
- 2. Aspirate media from cultures; wash cells with 1X PBS; aspirate.

- 3. Lyse cells by adding 1X SDS sample buffer (100 μl per well of 6-well plate or 500 μl per plate of 10 cm diameter plate). Immediately scrape the cells off the plate and transfer the extract to a microcentrifuge tube. Keep on ice.
- 4. Sonicate for 10-15 seconds for complete cell lysis and to shear DNA (to reduce sample viscosity).
- 5. Heat a 20 μ l sample to 95–100 °C for 5 minutes; cool on ice.
- 6. Microcentrifuge for 5 minutes.
- 7. Load 20 μl onto SDS-PAGE gel (10 cm x 10 cm).

NOTE: CST recommends loading prestained molecular weight marker (#7720, 10 µl/lane) to verify electrotransfer and biotinylated protein ladder (#7727, 10 µl/lane) to determine molecular weights.

1. Electrotransfer to nitrocellulose or PVDF membrane.

C. Membrane Blocking and Antibody Incubations

NOTE: Volumes are for 10 cm x 10 cm (100 cm²) of membrane; for different sized membranes, adjust volumes accordingly.

- 1. (Optional) After transfer, wash nitrocellulose membrane with 25 ml TBS for 5 minutes at room temperature.
- 2. Incubate membrane in 25 ml of blocking buffer for 1 hour at room temperature.
- 3. Wash three times for 5 minutes each with 15 ml of TBS/T.
- 4. Incubate membrane and primary antibody (at the appropriate dilution) in 10 ml primary antibody dilution buffer with gentle agitation overnight at 4 °C.
- 5. Wash three times for 5 minutes each with 15 ml of TBS/T.
- 6. Incubate membrane with appropriate HRP-conjugated secondary antibody (1:2000) and HRP-conjugated anti-biotin antibody (1:1000) to detect biotinylated protein markers in 10 ml of blocking buffer with gentle agitation for 1 hour at room temperature.
- 7. Wash three times for 5 minutes each with 15 ml of TBS/T.

D. Detection of Proteins

Incubate membrane with 10 ml LumiGLO® (0.5 ml 20X LumiGLO®, 0.5 ml 20X Peroxide and 9.0 ml Milli-Q water)
 with gentle agitation for 1 minute at room temperature.

NOTE: LumiGLO® substrate can be further diluted if signal response is too fast.

1. Drain membrane of excess developing solution (do not let dry), wrap in plastic wrap and expose to x-ray film. An initial 10-second exposure should indicate the proper exposure time.

NOTE: Due to the kinetics of the detection reaction, signal is most intense immediately following LumiGLO[®]incubation and declines over the following 2 hours.


Specificity / Sensitivity

XIAP Antibody detects endogenous levels of total XIAP protein. The antibody does not cross-react with other inhibitors of apoptosis.

Source / Purification


Polyclonal antibodies are produced by immunizing animals with a synthetic peptide corresponding to the amino terminus of human XIAP. Antibodies are purified by protein A and peptide affinity chromatography.

Western Blotting

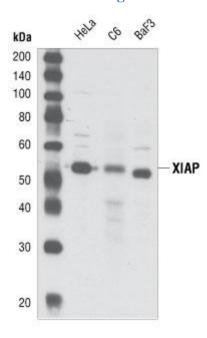


Western blot analysis of extracts from HeLa (human), C6 (rat) and BaF3 (mouse) cell lines, using XIAP Antibody.

Western Blotting

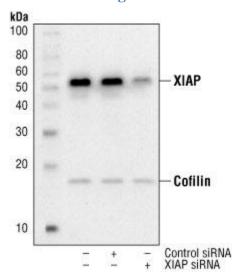

Western blot analysis of extracts from A549 cells 48 hours following mock transfection (-), transfection with non-targeted (control) siRNA or transfection with XIAP siRNA. XIAP was detected using XIAP Antibody #2042, and cofilin was detected using Cofilin Antibody #3312. The XIAP Antibody confirms silencing of XIAP expression, and the cofilin Antibody is used to control for protein loading and siRNA specificity.

Western Blotting

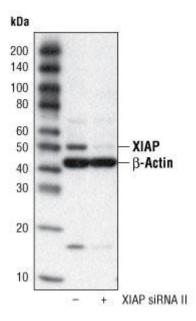

Western blot analysis of extracts from HeLa cells, transfected with 100 nM SignalSilence® Control siRNA (Fluorescein Conjugate) #6201 (-) or SignalSilence® XIAP siRNA II (+), using XIAP Antibody #2042 and β -Actin Antibody #4967. The XIAP antibody confirms silencing of XIAP expression, while the β -Actin Antibody is used to control for loading and specificity of XIAP siRNA.