

Applications	Reactivity	Sensitivity	MW (kDa)	Source
W IP	H M R Mk	Endogenous	115	Rabbit

Applications Key: W=Western Blotting IP=Immunoprecipitation

Reactivity Key: H=Human M=Mouse R=Rat Mk=Monkey

Species cross-reactivity is determined by western blot. Species enclosed in parentheses are predicted to react based on 100% sequence homology.

Protocols

Immunoprecipitation Protocol / (For Analysis By Western Immunoblotting)

For **shorter assay times** please try our [Immunoprecipitation Protocol Utilizing Magnetic Separation / \(For Analysis By Western Immunoblotting\)](#).

A. Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1X Phosphate Buffered Saline (PBS)
- 1X Cell Lysis Buffer:** ([#9803](#)) 20 mM Tris (pH 7.5), 150 mM NaCl, 1 mM EDTA, 1 mM EGTA, 1% Triton X-100, 2.5 mM Sodium pyrophosphate, 1 mM β -glycerophosphate, 1 mM Na_3VO_4 , 1 $\mu\text{g/ml}$ Leupeptin
NOTE: Add 1 mM PMSF immediately prior to use.
- Protein A or G Agarose Beads:** (Protein A [#9863](#)) Please prepare according to manufacturer's instructions. Use Protein A for rabbit IgG pull down and Protein G for mouse IgG pull down.
- 3X SDS Sample Buffer:** ([#7722](#)) 187.5 mM Tris-HCl (pH 6.8 at 25 °C), 6% w/v SDS, 30% glycerol, 150 mM DTT, 0.03% w/v bromophenol blue

B. Preparing Cell Lysates

- Aspirate media. Treat cells by adding fresh media containing regulator for desired time.
- To harvest cells under nondenaturing conditions, remove media and rinse cells once with ice-cold PBS.
- Remove PBS and add 0.5 ml ice-cold 1X cell lysis buffer to each plate (10 cm) and incubate the plates on ice for 5 minutes.
- Scrape cells off the plates and transfer to microcentrifuge tubes. Keep on ice.
- Sonicate samples on ice three times for 5 seconds each.
- Microcentrifuge for 10 minutes at 14,000 X g, 4 °C, and transfer the supernatant to a new tube. If necessary, lysate can be stored at -80 °C.

C. Immunoprecipitation

Optional: It may be necessary to perform a lysate pre-clearing step to reduce non-specific binding to the Protein A/G agarose beads (See section below).

- Take 200 μl cell lysate and add primary antibody. Incubate with gentle rocking overnight at 4 °C.
- Add either protein A or G agarose beads (20 μl of 50% bead slurry). Incubate with gentle rocking for 1–3 hours at 4 °C.
- Microcentrifuge for 30 seconds at 4 °C. Wash pellet five times with 500 μl of 1X cell lysis buffer. Keep on ice during washes.

4. Resuspend the pellet with 20 μ l 3X SDS sample buffer. Vortex, then microcentrifuge for 30 seconds.
5. Heat the sample to 95–100 $^{\circ}$ C for 2–5 minutes and microcentrifuge for 1 minute at 14,000 X g.
6. Load the sample (15–30 μ l) on SDS-PAGE gel (12–15%).
7. Analyze sample by Western blotting (see Western Immunoblotting Protocol: [Western BSA](#), [Western Milk](#)).

Cell Lysate Pre-Clearing (Optional)

1. Take 200 μ l cell lysate and add to either Protein A or G agarose beads (20 μ l of 50% bead slurry).
2. Incubate at 4 $^{\circ}$ C for 30 – 60 minutes.
3. Spin for 10 minutes at 4 $^{\circ}$ C. Transfer the supernatant to a fresh tube.
4. Proceed to step 1 of Immunoprecipitation.

NOTE: For proteins with molecular weights of 50 kDa, we recommend using [Mouse Anti-Rabbit IgG \(Light-Chain Specific\) \(L57A3\) mAb #3677](#) or [Mouse Anti-Rabbit IgG \(Conformation Specific\) \(L27A9\) mAb #3678](#) as a secondary antibody to minimize masking produced by denatured heavy chains. For proteins with molecular weights of 25 kDa, [Mouse Anti-Rabbit IgG \(Conformation Specific\) \(L27A9\) mAb #3678](#) is recommended.

Western Immunoblotting Protocol (Primary Ab Incubation In BSA)

For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4 $^{\circ}$ C with gentle shaking, overnight.

Products available from Cell Signaling Technology are linked by their respective catalog numbers.

A. Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

1. **1X Phosphate Buffered Saline (PBS).**
2. **1X SDS Sample Buffer:** ([#7722](#), [#7723](#)) 62.5 mM Tris-HCl (pH 6.8 at 25 $^{\circ}$ C), 2% w/v SDS, 10% glycerol, 50 mM DTT, 0.01% w/v bromophenol blue or phenol red.
3. **Transfer Buffer:** 25 mM Tris base, 0.2 M glycine, 20% methanol (pH 8.5).
4. **10X Tris Buffered Saline (TBS):** ([#9997](#)) To prepare 1 liter of 10X TBS: 24.2 g Tris base, 80 g NaCl; adjust pH to 7.6 with HCl (use at 1X).
5. **Nonfat Dry Milk:** ([#9999](#)) (weight to volume [w/v]).
6. **Blocking Buffer:** 1X TBS, 0.1% Tween-20 with 5% w/v nonfat dry milk; for 150 ml, add 15 ml 10X TBS to 135 ml water, mix. Add 7.5 g nonfat dry milk and mix well. While stirring, add 0.15 ml Tween-20 (100%).
7. **Wash Buffer:** 1X TBS, 0.1% Tween-20 (TBS/T).
8. **Bovine Serum Albumin (BSA):** ([#9998](#)).
9. **Primary Antibody Dilution Buffer:** 1X TBS, 0.1% Tween-20 with 5% BSA; for 20 ml, add 2 ml 10X TBS to 18 ml water, mix. Add 1.0 g BSA and mix well. While stirring, add 20 μ l Tween-20 (100%).
10. **Phototope[®]-HRP Western Blot Detection System:** ([#7071 anti-rabbit](#)) or ([#7072 anti-mouse](#)) Includes biotinylated protein ladder, secondary ([#7074 anti-rabbit](#)) or ([#7076 anti-mouse](#)) antibody conjugated to horseradish peroxidase (HRP), anti-biotin antibody conjugated to HRP, LumiGLO[®] chemiluminescent reagent and peroxide.
11. **Prestained Protein Marker, Broad Range (Premixed Format):** ([#7720](#)).

12. **Biotinylated Protein Ladder Detection Pack:** ([#7727](#)).
13. **Blotting Membrane:** This protocol has been optimized for nitrocellulose membranes, which CST recommends. PVDF membranes may also be used.

B. Protein Blotting

A general protocol for sample preparation is described below.

1. Treat cells by adding fresh media containing regulator for desired time.
2. Aspirate media from cultures; wash cells with 1X PBS; aspirate.
3. Lyse cells by adding 1X SDS sample buffer (100 µl per well of 6-well plate or 500 µl per plate of 10 cm diameter plate).
Immediately scrape the cells off the plate and transfer the extract to a microcentrifuge tube. Keep on ice.
4. Sonicate for 10–15 seconds for complete cell lysis and to shear DNA (to reduce sample viscosity).
5. Heat a 20 µl sample to 95–100 °C for 5 minutes; cool on ice.
6. Microcentrifuge for 5 minutes.
7. Load 20 µl onto SDS-PAGE gel (10 cm x 10 cm). **NOTE:** CST recommends loading prestained molecular weight markers ([#7720](#), 10 µl/lane) to verify electrotransfer and biotinylated protein ladder ([#7727](#), 10 µl/lane) to determine molecular weights.
8. Electrotransfer to nitrocellulose or PVDF membrane.

C. Membrane Blocking and Antibody Incubations

NOTE: Volumes are for 10 cm x 10 cm (100 cm²) of membrane; for different sized membranes, adjust volumes accordingly.

1. (Optional) After transfer, wash nitrocellulose membrane with 25 ml TBS for 5 minutes at room temperature.
2. Incubate membrane in 25 ml of blocking buffer for 1 hour at room temperature.
3. Wash three times for 5 minutes each with 15 ml of TBS/T.
4. Incubate membrane and primary antibody (at the appropriate dilution) in 10 ml primary antibody dilution buffer with gentle agitation overnight at 4 °C.
5. Wash three times for 5 minutes each with 15 ml of TBS/T.

I. For Unconjugated Primary Antibodies

1. Incubate membrane with **appropriate** HRP-conjugated secondary antibody (1:2000) and HRP-conjugated anti-biotin antibody (1:1000) to detect biotinylated protein markers in 10 ml of blocking buffer with gentle agitation for 1 hour at room temperature.
2. Wash three times for 5 minutes each with 15 ml of TBS/T.

II. For HRP Conjugated Primary Antibodies

Skip to Detection of Proteins (Step D).

III. For Biotinylated Primary Antibodies

1. Incubate membrane with HRP-Streptavidin (at the appropriate dilution) in milk for one hour with gentle agitation at room temperature.
2. Wash three times for 5 minutes each with 15 ml of TBS/T.

D. Detection of Proteins

1. Incubate membrane with 10 ml LumiGLO[®] (0.5 ml 20X LumiGLO[®], 0.5 ml 20X Peroxide and 9.0 ml Milli-Q water) with gentle agitation for 1 minute at room temperature. **NOTE:** LumiGLO[®] substrate can be further diluted if signal response is too fast.
2. Drain membrane of excess developing solution (do not let dry), wrap in plastic wrap and expose to x-ray film. An initial 10-second exposure should indicate the proper exposure time. **NOTE:** Due to the kinetics of the detection reaction, signal is most intense immediately following LumiGLO[®] incubation and declines over the following 2 hours.


Specificity / Sensitivity

Phospho-PKD/PKC μ (Ser916) Antibody detects endogenous levels of PKD1/PKC μ only when phosphorylated at serine 916. This antibody may also cross-react with isoform PKD2, in some species.

Source / Purification

Polyclonal antibodies are produced by immunizing animals with a synthetic phosphopeptide corresponding to residues surrounding Ser916 of mouse PKD. Antibodies are purified by protein A and peptide affinity chromatography.

Western Blotting

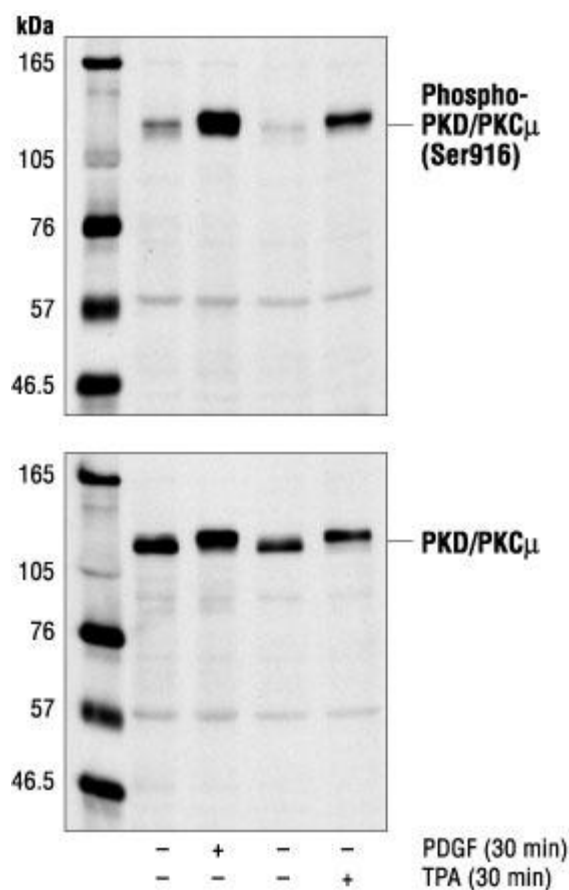

Western blot analysis of extracts from NIH/3T3 cells, untreated, PDGF-treated (50 ng/ml) or TPA-treated (0.2 μ M), using Phospho-PKD/PKC μ (Ser916) Antibody (upper) or PKD/PKC μ Antibody #2052 (lower).