


## Anti-Mouse TER-119 eFluor 450

Catalog Number: 48-5921

Also Known As: TER119, Erythroid cell marker, Ly-76, Ly76

RUO: For Research Use Only. Not for use in diagnostic procedures.

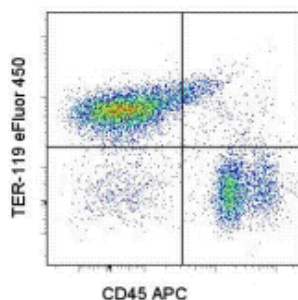

Staining of C57BL/6 bone marrow cells with Anti-Mouse CD45 APC (cat. 17-0451) and 0.25  $\mu$ g of Anti-Mouse TER-119 eFluor® 450. Total viable cells were used for analysis.

### Product Information

Contents: Anti-Mouse TER-119 eFluor 450


**REF** Catalog Number: 48-5921

Clone: TER-119


Concentration: 0.2 mg/mL

Host/Isotype: Rat IgG2b, kappa

Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

 Temperature Limitation: Store at 2-8°C. Do not freeze. Light sensitive material.

**LOT** Batch Code: Refer to Vial

 Use By: Refer to Vial

 Contains sodium azide

### Description

The TER-119 monoclonal antibody reacts with mouse erythroid cells from early proerythroblast to mature erythrocyte stages. The TER-119 antigen is present in yolk sac, fetal and newborn liver, but is not expressed by cells carrying BFU-E and CFU-E activities. Several erythroleukemia cell lines tested so far are negative for expression of TER-119 antigen even after dimethylsulfoxide stimulation. Biochemical and molecular analysis of the TER-119 antigen indicate that this molecule is associated with the surface glycoprotein A, but is not a typical glycoprotein.

### Applications Reported

This TER-119 antibody has been reported for use in flow cytometric analysis.

eFluor® 450 is a replacement for Pacific Blue®. eFluor® 450 emits at 456 nm and is excited with the Violet laser. Please make sure that your instrument is capable of detecting this fluorochrome.

### Applications Tested

This TER-119 antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This can be used at less than or equal to 0.5  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from  $10^5$  to  $10^8$  cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

eFluor™ 450 is a replacement for Pacific Blue®. eFluor™ 450 emits at 456 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochrome.

### References

Kina, T., K. Ikuta, et al. (2000). The monoclonal antibody TER-119 recognizes a molecule associated with glycoprotein A and specifically marks the late stages of murine erythroid lineage. *Br J Haematol* 109(2): 280-87.

Vannucchi, A. M., F. Paoletti, et al. (2000). Identification and characterization of a bipotent (erythroid and megakaryocytic) cell precursor from the spleen of phenylhydrazine-treated mice. *Blood* 95(8): 2559-68.

### Related Products

17-0451 Anti-Mouse CD45 APC (30-F11)

48-4031 Rat IgG2b K Isotype Control eFluor® 450

---

Not for further distribution without written consent.

Copyright © 2000-2010 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • [www.eBioscience.com](http://www.eBioscience.com) • [info@eBioscience.com](mailto:info@eBioscience.com)