

Applications	Reactivity	Sensitivity	MW (kDa)	Isotype
--------------	------------	-------------	----------	---------

W IHC-P F	H	Endogenous	18	Rabbit
-----------	---	------------	----	--------

Applications Key: W=Western Blotting IHC-P=Immunohistochemistry (Paraffin) F=Flow Cytometry

Reactivity Key: H=Human

Species cross-reactivity is determined by western blot. Species enclosed in parentheses are predicted to react based on 100% sequence homology.

Protocols

Flow Cytometry Protocol

A. Solutions and Reagents

NOTE: Prepare solutions with purified water.

1. **1X Phosphate Buffered Saline (PBS):** Dissolve 8 g NaCl, 0.2 g KCl, 1.44 g Na₂HPO₄ and 0.24 g KH₂PO₄ in 800 ml dH₂O. Adjust the pH to 7.4 with HCl and the volume to 1 L. Store at room temperature.
2. **Formaldehyde** (methanol free).
3. **100% Methanol**
4. **Incubation Buffer:** Dissolve 0.5 g bovine serum albumin (BSA) in 100 ml 1X PBS. Store at 4 °C.

B. Fixation

1. Collect cells by centrifugation and aspirate supernatant.
2. Resuspend cells briefly in 0.5–1 ml PBS. Add formaldehyde to a final concentration of 2–4% formaldehyde.
3. Fix for 10 min at 37 °C.
4. Chill tubes on ice for 1 min.
5. For extracellular staining with antibodies that do not require permeabilization, proceed to Section D, Step 1 or store cells in PBS with 0.1% sodium azide at 4 °C; for intracellular staining, proceed to permeabilization (Section C, Step 1).

C. Permeabilization

1. Permeabilize cells by adding ice-cold 100% methanol slowly to pre-chilled cells, while gently vortexing, to a final concentration of 90% methanol. Alternatively, to remove fix prior to permeabilization, pellet cells by centrifugation and resuspend in 90% methanol.
2. Incubate 30 min on ice.
3. Proceed with immunostaining (Section D, Step 1) or store cells at –20 °C in 90% methanol.

D. Immunostaining

NOTE: Account for isotype matched controls for monoclonal antibodies or species matched IgG for polyclonal antibodies. Count cells using a hemocytometer or alternative method.

1. Aliquot 0.5–1x10⁶ cells into each assay tube (by volume).
2. Add 2–3 ml Incubation Buffer to each tube and rinse by centrifugation. Repeat.
3. Resuspend cells in 100 µl Incubation Buffer per assay tube.
4. Block in Incubation Buffer for 10 min at room temperature.

5. Add the unconjugated, biotinylated, or fluorochrome-conjugated primary antibody at the appropriate dilution to the assay tubes (see individual antibody datasheet for the appropriate dilution).
6. Incubate for 1 hr at room temperature.
7. Rinse as before in Incubation Buffer by centrifugation.
8. If using a fluorochrome-conjugated primary antibody, resuspend cells in 0.5 ml PBS and analyze on flow cytometer; for unconjugated or biotinylated primary antibodies, proceed to immunostaining (Section D, Step 9).
9. Resuspend cells in fluorochrome-conjugated secondary antibody or fluorochrome-conjugated avidin, diluted in Incubation Buffer at the recommended dilution.
10. Incubate for 30 min at room temperature.
11. Rinse as before in Incubation Buffer by centrifugation.
12. Resuspend cells in 0.5 ml PBS and analyze on flow cytometer; alternatively, for DNA staining, proceed to optional DNA stain (Section E, Step 1).

E. Optional DNA Stain

1. Resuspend cells in 0.5 ml of DNA dye (e.g. Propidium Iodide (PI)/RNase Staining Solution [#4087](#)).
2. Incubate for at least 5 min at room temperature.
3. Analyze cells in DNA stain on flow cytometer.

Immunohistochemistry Protocol (Paraffin)

***IMPORTANT:** See product data sheet for the appropriate antibody diluent and antigen unmasking procedure. **IHC**

Protocol: Unmasking buffer/antibody diluent.

A. Solutions and Reagents

1. Xylene
2. Ethanol, anhydrous denatured, histological grade (100% and 95%)
3. Deionized water (dH₂O)
4. Hematoxylin (optional)
5. **Wash Buffer:**

1X TBS/0.1% Tween-20 (1X TBST): To prepare 1 L add 100 ml 10X TBS to 900 ml dH₂O. Add 1 ml Tween-20 and mix.

10X Tris Buffered Saline (TBS): To prepare 1 L add 24.2 g Trizma[®] base (C₄H₁₁NO₃) and 80 g sodium chloride (NaCl) to 1 L dH₂O. Adjust pH to 7.6 with concentrated HCl.

6. ***Antibody Diluent:**

- a. **SignalStain[®] Antibody Diluent #8112**
- b. **TBST/5% normal goat serum ([#5425](#)):** To 5 ml 1X TBST add 250 µl normal goat serum.
- c. **PBST/5% normal goat serum ([#5425](#)):** To 5 ml 1X PBST add 250 µl normal goat serum.

1X PBS/0.1% Tween-20 (1X PBST): To prepare 1 L add 100 ml 10X PBS to 900 ml dH₂O. Add 1 ml Tween-20 and mix.

10X Phosphate Buffered Saline (PBS): To prepare 1 L add 80 g sodium chloride (NaCl), 2 g potassium chloride (KCl),

14.4 g sodium phosphate, dibasic (Na_2HPO_4) and 2.4 g potassium phosphate, monobasic (KH_2PO_4) to 1 L dH_2O . Adjust pH to 7.4.

7. ***Antigen Unmasking:**

- a. **Citrate:** 10 mM Sodium Citrate Buffer: To prepare 1 L add 2.94 g sodium citrate trisodium salt dihydrate ($\text{C}_6\text{H}_5\text{Na}_3\text{O}_7 \cdot 2\text{H}_2\text{O}$) to 1 L dH_2O . Adjust pH to 6.0.
- b. **EDTA:** 1 mM EDTA: To prepare 1 L add 0.372 g EDTA ($\text{C}_{10}\text{H}_{14}\text{N}_2\text{O}_8\text{Na}_2 \cdot 2\text{H}_2\text{O}$) to 1 L dH_2O . Adjust pH to 8.0.
- c. **TE:** 10 mM Tris/1 mM EDTA, pH 9.0: To prepare 1L add 1.21 g Trizma[®] base ($\text{C}_4\text{H}_{11}\text{NO}_3$) and 0.372 g EDTA ($\text{C}_{10}\text{H}_{14}\text{N}_2\text{O}_8\text{Na}_2 \cdot 2\text{H}_2\text{O}$) to 950 ml dH_2O . Adjust pH to 9.0, then adjust final volume to 1000 ml with dH_2O .
- d. **Pepsin:** 1 mg/ml in Tris-HCl pH 2.0.
8. **3% Hydrogen Peroxide:** To prepare, add 10 ml 30% H_2O_2 to 90 ml dH_2O .
9. **Blocking Solution:** TBST/5% normal goat serum ([#5425](#)): to 5 ml 1X TBST add 250 μl normal goat serum.
10. Biotinylated secondary antibody.
11. **ABC Reagent:** (Vectastain ABC Kit, Vector Laboratories, Inc., Burlingame, CA) Prepare according to manufacturer's instructions 30 minutes before use.
12. **DAB Reagent or suitable substrate:** Prepare according to manufacturer's recommendations.

B. Deparaffinization/Rehydration

NOTE: Do not allow slides to dry at any time during this procedure.

1. **Deparaffinize/hydrate sections:**
 - a. Incubate sections in three washes of xylene for 5 minutes each.
 - b. Incubate sections in two washes of 100% ethanol for 10 minutes each.
 - c. Incubate sections in two washes of 95% ethanol for 10 minutes each.
2. Wash sections twice in dH_2O for 5 minutes each.

C. *Antigen Unmasking

NOTE: Consult product data sheet for specific recommendation for the unmasking solution.

1. **For Citrate:** Bring slides to a boil in 10 mM sodium citrate buffer pH 6.0 then maintain at a sub-boiling temperature for 10 minutes. Cool slides on bench top for 30 minutes.
2. **For EDTA:** Bring slides to a boil in 1 mM EDTA pH 8.0 followed by 15 minutes at a sub-boiling temperature. No cooling is necessary.
3. **For TE:** Bring slides to a boil in 10 mM TE/1 mM EDTA, pH 9.0 then maintain at a sub-boiling temperature for 18 minutes. Cool on the bench for 30 minutes.
4. **For Pepsin:** Digest for 10 minutes at 37 °C.

D. Staining

1. Wash sections in dH_2O three times for 5 minutes each.
2. Incubate sections in 3% hydrogen peroxide for 10 minutes.
3. Wash sections in dH_2O twice for 5 minutes each.

NOTE: Consult product data sheet for recommended antibody diluent.

4. Wash sections in wash buffer for 5 minutes.
5. Block each section with 100-400 µl blocking solution for 1 hour at room temperature.
6. Remove blocking solution and add 100-400 µl primary antibody diluted in recommended antibody diluent to each section.
Incubate overnight at 4 °C.
7. Remove antibody solution and wash sections in wash buffer three times for 5 minutes each.
8. Add 100-400 µl biotinylated secondary antibody, diluted in TBST per manufacturer's recommendation, to each section.
Incubate 30 minutes at room temperature.
9. If using ABC avidin/biotin method, prepare ABC reagent according to the manufacturer's instructions and incubate solution for 30 minutes at room temperature.
10. Remove secondary antibody solution and wash sections three times with wash buffer for 5 minutes each.
11. Add 100-400 µl ABC reagent to each section and incubate for 30 minutes at room temperature.
12. Remove ABC reagent and wash sections three times in wash buffer for 5 minutes each.
13. Add 100-400 µl DAB or suitable substrate to each section and monitor staining closely.
14. As soon as the sections develop, immerse slides in dH₂O.
15. If desired, counterstain sections in hematoxylin per manufacturer's instructions.
16. Wash sections in dH₂O two times for 5 minutes each.
17. Dehydrate sections:
 - a. Incubate sections in 95% ethanol two times for 10 seconds each.
 - b. Repeat in 100% ethanol, incubating sections two times for 10 seconds each.
 - c. Repeat in xylene, incubating sections two times for 10 seconds each.
18. Mount coverslips.

Western Immunoblotting Protocol (Primary Ab Incubation In BSA)

For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4 °C with gentle shaking, overnight.

Products available from Cell Signaling Technology are linked by their respective catalog numbers.

A. Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

1. **1X Phosphate Buffered Saline (PBS).**
2. **1X SDS Sample Buffer:** ([#7722](#), [#7723](#)) 62.5 mM Tris-HCl (pH 6.8 at 25 °C), 2% w/v SDS, 10% glycerol, 50 mM DTT, 0.01% w/v bromophenol blue or phenol red.
3. **Transfer Buffer:** 25 mM Tris base, 0.2 M glycine, 20% methanol (pH 8.5).
4. **10X Tris Buffered Saline (TBS):** ([#9997](#)) To prepare 1 liter of 10X TBS: 24.2 g Tris base, 80 g NaCl; adjust pH to 7.6 with HCl (use at 1X).
5. **Nonfat Dry Milk:** ([#9999](#)) (weight to volume [w/v]).
6. **Blocking Buffer:** 1X TBS, 0.1% Tween-20 with 5% w/v nonfat dry milk; for 150 ml, add 15 ml 10X TBS to 135 ml water, mix. Add 7.5 g nonfat dry milk and mix well. While stirring, add 0.15 ml Tween-20 (100%).

7. **Wash Buffer:** 1X TBS, 0.1% Tween-20 (TBS/T).
8. **Bovine Serum Albumin (BSA):** ([#9998](#)).
9. **Primary Antibody Dilution Buffer:** 1X TBS, 0.1% Tween-20 with 5% BSA; for 20 ml, add 2 ml 10X TBS to 18 ml water, mix. Add 1.0 g BSA and mix well. While stirring, add 20 µl Tween-20 (100%).
10. **Phototope®-HRP Western Blot Detection System:** ([#7071 anti-rabbit](#)) or ([#7072 anti-mouse](#)) Includes biotinylated protein ladder, secondary ([#7074 anti-rabbit](#)) or ([#7076 anti-mouse](#)) antibody conjugated to horseradish peroxidase (HRP), anti-biotin antibody conjugated to HRP, LumiGLO® chemiluminescent reagent and peroxide.
11. **Prestained Protein Marker, Broad Range (Premixed Format):** ([#7720](#)).
12. **Biotinylated Protein Ladder Detection Pack:** ([#7727](#)).
13. **Blotting Membrane:** This protocol has been optimized for nitrocellulose membranes, which CST recommends. PVDF membranes may also be used.

B. Protein Blotting

A general protocol for sample preparation is described below.

1. Treat cells by adding fresh media containing regulator for desired time.
2. Aspirate media from cultures; wash cells with 1X PBS; aspirate.
3. Lyse cells by adding 1X SDS sample buffer (100 µl per well of 6-well plate or 500 µl per plate of 10 cm diameter plate).
Immediately scrape the cells off the plate and transfer the extract to a microcentrifuge tube. Keep on ice.
4. Sonicate for 10–15 seconds for complete cell lysis and to shear DNA (to reduce sample viscosity).
5. Heat a 20 µl sample to 95–100 °C for 5 minutes; cool on ice.
6. Microcentrifuge for 5 minutes.
7. Load 20 µl onto SDS-PAGE gel (10 cm x 10 cm). **NOTE:** CST recommends loading prestained molecular weight markers ([#7720](#), 10 µl/lane) to verify electrotransfer and biotinylated protein ladder ([#7727](#), 10 µl/lane) to determine molecular weights.
8. Electrotransfer to nitrocellulose or PVDF membrane.

C. Membrane Blocking and Antibody Incubations

NOTE: Volumes are for 10 cm x 10 cm (100 cm²) of membrane; for different sized membranes, adjust volumes accordingly.

1. (Optional) After transfer, wash nitrocellulose membrane with 25 ml TBS for 5 minutes at room temperature.
2. Incubate membrane in 25 ml of blocking buffer for 1 hour at room temperature.
3. Wash three times for 5 minutes each with 15 ml of TBS/T.
4. Incubate membrane and primary antibody (at the appropriate dilution) in 10 ml primary antibody dilution buffer with gentle agitation overnight at 4 °C.
5. Wash three times for 5 minutes each with 15 ml of TBS/T.

I. For Unconjugated Primary Antibodies

1. Incubate membrane with **appropriate** HRP-conjugated secondary antibody (1:2000) and HRP-conjugated anti-biotin antibody (1:1000) to detect biotinylated protein markers in 10 ml of blocking buffer with gentle agitation for 1 hour at room temperature.
2. Wash three times for 5 minutes each with 15 ml of TBS/T.

II. For HRP Conjugated Primary Antibodies

Skip to Detection of Proteins (Step D).

III. For Biotinylated Primary Antibodies

1. Incubate membrane with HRP-Streptavidin (at the appropriate dilution) in milk for one hour with gentle agitation at room temperature.
2. Wash three times for 5 minutes each with 15 ml of TBS/T.

D. Detection of Proteins

1. Incubate membrane with 10 ml LumiGLO® (0.5 ml 20X LumiGLO®, 0.5 ml 20X Peroxide and 9.0 ml Milli-Q water) with gentle agitation for 1 minute at room temperature. **NOTE:** LumiGLO® substrate can be further diluted if signal response is too fast.
2. Drain membrane of excess developing solution (do not let dry), wrap in plastic wrap and expose to x-ray film. An initial 10-second exposure should indicate the proper exposure time. **NOTE:** Due to the kinetics of the detection reaction, signal is most intense immediately following LumiGLO® incubation and declines over the following 2 hours.

Specificity / Sensitivity


SCF (C19H6) Rabbit mAb detects endogenous levels of SCF by immunohistochemical staining and by flow cytometry.

Western blotting detects recombinant SCF.

Source / Purification


Monoclonal antibody is produced by immunizing animals with a synthetic peptide corresponding to amino acids in the N-terminal sequence of human SCF.

Western Blotting

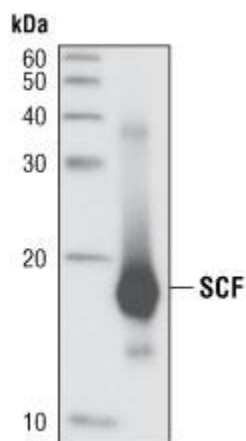


Western blot analysis of recombinant SCF (5 ng) using SCF (C19H6) Rabbit mAb.

IHC-P (paraffin)

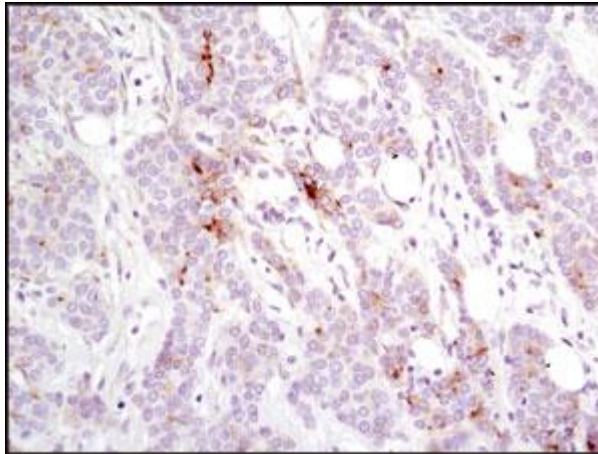

Immunohistochemical analysis of paraffin-embedded human breast carcinoma using SCF (C19H6) RmAb.

Flow Cytometry

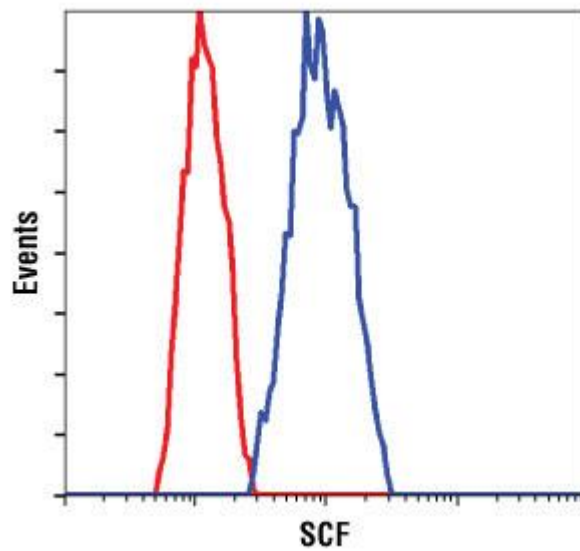

Flow cytometric analysis of U-87 MG cells using SCF (C19H6) Rabbit mAb (blue) compared to a nonspecific negative control antibody (red).