



An Affymetrix Company

# Anti-Mouse CD19 eFluor® 660 (Alexa Fluor® 647 Replacement)

Catalog Number: 50-0193

RUO: For Research Use Only. Not for use in diagnostic procedures.

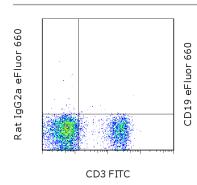

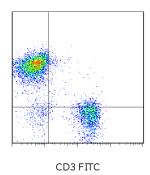

Staining of BALB/c splenocytes with Anti-Mouse CD3e FITC (cat. 11-0031) and 0.125 ug of Rat IgG2a K Isotype Control eFluor® 660 (cat. 50-4321) (left) or 0.125 ug of Anti-Mouse CD19 eFluor® 660 (right). Total viable cells were used for analysis.


#### **Product Information**

Contents: Anti-Mouse CD19 eFluor® 660 (Alexa Fluor® 647 Replacement)

**REF** Catalog Number: 50-0193 Clone: eBio1D3 (1D3) Concentration: 0.2 mg/mL Host/Isotype: Rat IgG2a, kappa



Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer Temperature Limitation: Store at 2-8°C. Do not

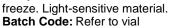



Use By: Refer to vial Contains sodium azide

## Description

The eBio1D3 (1D3) monoclonal antibody reacts with mouse CD19, a 95 kDa transmembrane glycoprotein. CD19 is expressed by B cells during all stages of development excluding the terminally differentiated plasma cells. Follicular dendritic cells also express CD19. Together CD21, CD81, MHC class II, and CD19 form a multimolecular complex that associates with the BCR. Signaling through CD19 induces tyrosine phosphorylation, calcium flux and proliferation of B cells.

### **Applications Reported**

This eBio1D3 (1D3) antibody has been reported for use in flow cytometric analysis.

## **Applications Tested**

This eBio1D3 (1D3) antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

eFluor® 660 is a replacement for Alexa Fluor® 647. eFluor® 660 emits at 659 nm and is excited with the red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochome.

### References

Cherukuri A, Cheng PC, Pierce SK. The role of the CD19/CD21 complex in B cell processing and presentation of complement-tagged antigens. J Immunol. 2001 Jul 1;167(1):163-72.

Inabe K, Kurosaki T. Tyrosine phosphorylation of B-cell adaptor for phosphoinositide 3-kinase is required for Akt activation in response to CD19 engagement. Blood. 2002 Jan 15;99(2):584-9.

Krop I, de Fougerolles AR, Hardy RR, Allison M, Schlissel MS, Fearon DT. Self-renewal of B-1 lymphocytes is dependent on CD19. Eur J Immunol. 1996 Jan;26(1):238-42. (1D3, FA, PubMed)



An Affymetrix Company

## **Anti-Mouse CD19 eFluor® 660 (Alexa Fluor® 647 Replacement)**

Catalog Number: 50-0193

RUO: For Research Use Only. Not for use in diagnostic procedures.

Shoham T, Rajapaksa R, Boucheix C, Rubinstein E, Poe JC, Tedder TF, Levy S. The tetraspanin CD81 regulates the expression of CD19 during B cell development in a postendoplasmic reticulum compartment. J Immunol. 2003 Oct 15;171(8):4062-72.

Krop I, Shaffer AL, Fearon DT, Schlissel MS. The signaling activity of murine CD19 is regulated during cell development. J Immunol. 1996 Jul 1;157(1):48-56. (1D3, FC, FA, PubMed)

### **Related Products**

11-0031 Anti-Mouse CD3e FITC (145-2C11) 12-0452 Anti-Human/Mouse CD45R (B220) PE (RA3-6B2) 25-0232 Anti-Mouse CD23 PE-Cy7 (B3B4)

50-4321 Rat IgG2a K Isotype Control eFluor® 660 (eBR2a)