


Anti-Mouse CD144 (VE-Cadherin) eFluor[®] 660 (Alexa Fluor[®] 647 Replacement)

Catalog Number: 50-1441

Also known as: Cdh5, Cadherin 5, Vascular endothelial cadherin

RUO: For Research Use Only. Not for use in diagnostic procedures.

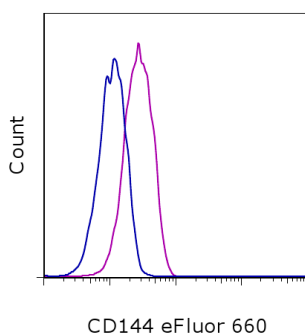

Staining of bEnd.3 cells with 0.5 ug of Rat IgG1 K Isotype Control eFluor[®] 660 (cat. 50-4301) (blue histogram) or 0.5 ug of Anti-Mouse CD144 (VE-Cadherin) eFluor[®] 660 (Alexa Fluor[®] 647 Replacement) (purple histogram). Total viable cells were used for analysis.

Product Information

Contents: Anti-Mouse CD144 (VE-Cadherin) eFluor[®] 660 (Alexa Fluor[®] 647 Replacement)

REF **Catalog Number:** 50-1441

Clone: eBioBV13 (BV13)

Concentration: 0.2 mg/mL

Host/Isotype: Rat IgG1

Formulation: aqueous buffer, 0.09% sodium azide, may contain carrier protein/stabilizer

Temperature Limitation: Store at 2-8°C. Do not freeze. Light-sensitive material.

Batch Code: Refer to vial

Use By: Refer to vial

Contains sodium azide



Description

The eBioBV13 (BV13) monoclonal antibody reacts with mouse VE-Cadherin (CD144). VE-Cadherin is a 120 kDa member of the type II Cadherin family, characterized by the presence of 5 extracellular cadherin domains (ECD), and anchored to the actin cytoskeleton through their cytoplasmic tail. VE-Cadherin mediates homophilic adhesion between neighbouring endothelial cells and is localized within specialized structures at cell-cell contacts, called adherens junctions. VE-Cadherin is expressed constitutively throughout the entire vasculature, and is required for numerous endothelial cell functions including migration, survival, contact-dependent growth inhibition and endothelial cell assembly into tubular structures. Furthermore, it is thought that VE-Cadherin+CD45- cells from the yolk sac or aorta-gonad-mesonephros (AGM) have the potential to give rise to hematopoietic cells.

Applications Reported

This eBioBV13 (BV13) antibody has been reported for use in flow cytometric analysis.

Applications Tested

This eBioBV13 (BV13) antibody has been tested by flow cytometric analysis of bEnd.3 cells. This can be used at less than or equal to 1 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

eFluor[®] 660 is a replacement for Alexa Fluor[®] 647. eFluor[®] 660 emits at 659 nm and is excited with the red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochrome.

References

Corada M, Mariotti M, Thurston G, Smith K, Kunkel R, Brockhaus M, Lampugnani MG, Martin-Padura I, Stoppacciaro A, Ruco L, McDonald DM, Ward PA, Dejana E. Vascular endothelial-cadherin is an important determinant of microvascular integrity in vivo. Proc Natl Acad Sci U S A. 1999 Aug 17;96(17):9815-20. (BV13, FA, IHC, PubMed)

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.ebioscience.com •
info@ebioscience.com

Anti-Mouse CD144 (VE-Cadherin) eFluor® 660 (Alexa Fluor® 647 Replacement)

Catalog Number: 50-1441

Also known as: Cdh5, Cadherin 5, Vascular endothelial cadherin

RUO: For Research Use Only. Not for use in diagnostic procedures.

Liao F, Li Y, O'Connor W, Zanetta L, Bassi R, Santiago A, Overholser J, Hooper A, Mignatti P, Dejana E, Hicklin DJ, Bohlen P. Monoclonal antibody to vascular endothelial-cadherin is a potent inhibitor of angiogenesis, tumor growth, and metastasis. *Cancer Res.* 2000 Dec 15;60(24):6805-10. (**BV13**, FA, PubMed)

Crosby CV, Fleming PA, Argraves WS, Corada M, Zanetta L, Dejana E, Drake CJ. VE-cadherin is not required for the formation of nascent blood vessels but acts to prevent their disassembly. *Blood.* 2005 Apr 1;105(7):2771-6. Epub 2004 Dec 16. (**BV13**, FA, PubMed)

Related Products

00-4222 Flow Cytometry Staining Buffer

50-4301 Rat IgG1 K Isotype Control eFluor® 660

Not for further distribution without written consent.

Copyright © 2000-2012 eBioscience, Inc.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • www.ebioscience.com •
info@ebioscience.com