

PRODUCT INFORMATION & MANUAL

Human sICAM-3 Platinum ELISA

BMS218 / BMS218TEN

Enzyme-linked Immunosorbent Assay for
quantitative detection of human sICAM-3.

For research use only.

Not for diagnostic or therapeutic procedures.



Human sICAM-3 Platinum ELISA

North America

Technical Support:

Research Products:
888.810.6168
858.642.2058
tech@eBioscience.com

Clinical Products:
877.726.8559
858.642.2058
tech@eBioscience.com

Customer Service:

888.999.1371
858.642.2058
info@eBioscience.com

Fax:

858.642.2046

Europe/International*

Technical Support:

+43 1 796 40 40-120
tech@eBioscience.com

Customer Service:

+43 1 796 40 40-304
info@eBioscience.com

Fax:

+43 1 796 40 40-400



Bender MedSystems GmbH
Campus Vienna Biocenter 2
1030 Vienna, Austria
www.eBioscience.com

* Customers outside North America and Europe may contact their eBioscience distributor listed on our website at www.eBioscience.com/distributors.

TABLE OF CONTENTS

1	Intended Use	3
2	Summary	3
3	Principles of the Test	5
4	Reagents Provided	6
5	Storage Instructions – ELISA Kit	8
6	Specimen Collection and Storage Instructions	8
7	Materials Required But Not Provided	9
8	Precautions for Use	10
9	Preparation of Reagents	12
10	Test Protocol	15
11	Calculation of Results	20
12	Limitations	23
13	Performance Characteristics	24
14	Ordering Information	29
15	Reagent Preparation Summary	29
16	Test Protocol Summary	31

1 Intended Use

The human sICAM-3 ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of human sICAM-3. **The human sICAM-3 ELISA is for research use only. Not for diagnostic or therapeutic procedures.**

2 Summary

Intercellular Adhesion Molecule-3 (ICAM-3) is a member of the immunoglobulin supergene family and functions as a ligand for the Lymphocyte Function-Associated Antigen-1 (LFA-1). Three counter-receptors have been described for LFA-1, intercellular adhesion molecule 1 (ICAM-1), ICAM-2 and ICAM-3. LFA-1, an alpha-beta complex, is a member of the leukocyte integrin family which mediate lymphocyte adhesion.

ICAM-3 is a heavily glycosylated protein of 124 kDa with a polypeptide core of 57 kDa. The integral membrane protein with five immunoglobulin-like domains shares high homology to ICAM-1 and ICAM-2 in the extracellular region. In contrast to ICAM-1 and ICAM-2, ICAM-3 is absent on endothelia. ICAM-3 is expressed on resting lymphocytes, monocytes and neutrophils, representing the major LFA-1 ligand on these cells. The finding that adhesion of resting T lymphocytes to LFA-1 occurs primarily via ICAM-3 combined with the fact that ICAM-3 is much better expressed than other LFA-1 ligands on monocytes and resting lymphocytes implies an important role for ICAM-3 in the initiation of immune responses.

ICAM-3 was found to be involved in the regulation of LFA-1/ICAM-1 dependent leukocyte intercellular interactions. The initial interaction of ICAM-3 with LFA-1 might increase LFA-1-mediated cell binding to ICAM-1.

Furthermore, ICAM-3 expression has been shown for dendritic epidermal Langerhans cells, whereas it is absent on other dendritic cells from different lymphoid organs. Thus potential function of ICAM-3 at the initiation phase of LC-leukocyte interactions taking place during skin localized immune reactions can be postulated.

Recent data suggest that ICAM-3 expression can be induced on endothelial cells in lymphoid neoplasms as shown for Hodgkin's and non-Hodgkin's disease.

ICAM-3 is a very interesting molecule involved in the initial immune response thus suggesting an important role as a disease marker for a number of different indications and pathological situations.

For literature update refer to **www.eBioscience.com**

3 Principles of the Test

An anti-human sICAM-3 coating antibody is adsorbed onto microwells.

Human sICAM-3 present in the sample or standard binds to antibodies adsorbed to the microwells. A HRP-conjugated anti-human sICAM-3 antibody is added and binds to human sICAM-3 captured by the first antibody.

Following incubation unbound HRP-conjugated anti-human sICAM-3 is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of human sICAM-3 present in the sample or standard. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from 7 human sICAM-3 standard dilutions and human sICAM-3 concentration determined.

Figure 1

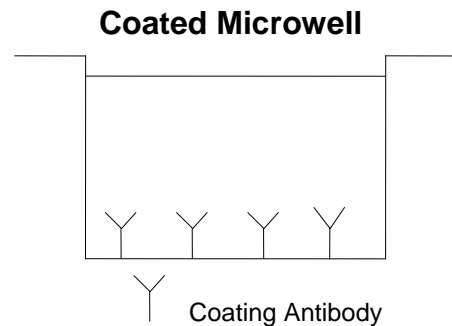


Figure 2

First Incubation

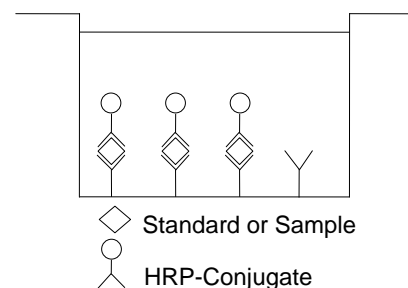


Figure 3

Second Incubation


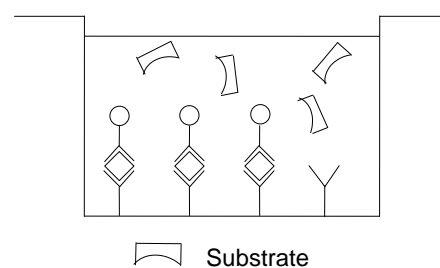


Figure 4

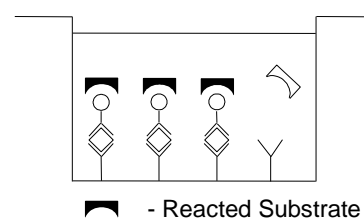

4 Reagents Provided

4.1 Reagents for human sICAM-3 ELISA BMS218 (96 tests)

- 1 aluminium pouch with a **Microwell Plate coated** with monoclonal antibody to human sICAM-3
- 1 vial (6 ml) **HRP-Conjugate** anti-human sICAM-3 monoclonal antibody
- 2 vials (500 µl) human sICAM-3 **Standard**, 50 ng/ml
- 1 vial (12 ml) **Sample Diluent**
- 1 bottle (50 ml) **Wash Buffer Concentrate 20x** (PBS with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 1 vial (0.4 ml) **Blue-Dye**
- 1 vial (0.4 ml) **Green-Dye**
- 2 **Adhesive Films**

4.2 Reagents for human sICAM-3 ELISA BMS218TEN (10x96 tests)

- 10 aluminium pouches with a **Microwell Plate coated** with monoclonal antibody to human sICAM-3
- 10 vials (6 ml) **HRP-Conjugate** anti-human sICAM-3 monoclonal antibody
- 10 vials (500 µl) human sICAM-3 **Standard**, 50 ng/ml
- 10 vials (12 ml) **Sample Diluent**
- 3 bottles (50 ml) **Wash Buffer Concentrate 20x**
(PBS with 1% Tween 20)
- 10 vials (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 10 vials (15 ml) **Stop Solution** (1M Phosphoric acid)
- 6 vials (0.4 ml) **Blue-Dye**
- 6 vials (0.4 ml) **Green-Dye**
- 10 **Adhesive Films**

5 Storage Instructions – ELISA Kit

Store kit reagents between 2° and 8°C. Immediately after use remaining reagents should be returned to cold storage (2° to 8°C). Expiry of the kit and reagents is stated on labels.

Expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

6 Specimen Collection and Storage Instructions

Cell culture supernatant, serum, plasma (EDTA, heparin), amniotic fluid and urine were tested with this assay. Other biological samples might be suitable for use in the assay.

Remove serum or plasma from the clot or cells as soon as possible after clotting and separation.

Pay attention to a possible “**Hook Effect**” due to high sample concentrations (see chapter 11).

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20°C to avoid loss of bioactive human sICAM-3. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13.5).

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- 50 μ l to 300 μ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.

- Exposure to acid inactivates the conjugate.
- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents

Buffer Concentrates should be brought to room temperature and should be diluted before starting the test procedure.

If crystals have formed in the **Buffer Concentrates**, warm them gently until they have completely dissolved.

9.1 Wash Buffer (1x)

Pour entire contents (50 ml) of the **Wash Buffer Concentrate** (20x) into a clean 1000 ml graduated cylinder. Bring to final volume of 1000 ml with glass-distilled or deionized water.

Mix gently to avoid foaming.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer (1x) is stable for 30 days.

Wash Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Wash Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

9.2 Human sICAM-3 Standard

Standard dilutions can be prepared directly on the microwell plate (see 10.c) or alternatively in tubes (see 9.2.1).

9.2.1 External Standard Dilution

Label 6 tubes, one for each standard point.

S2, S3, S4, S5, S6, S7

Then prepare 1:2 serial dilutions for the standard curve as follows:

Pipette 225 µl of Sample Diluent into tubes S2 – S7.

Pipette 225 µl of undiluted standard (serves as the highest standard S1, concentration of standard 1= 50 ng/ml) into the first tube, labelled S2,


and mix (concentration of standard 2 = 25 ng/ml).

Pipette 225 μ l of this dilution into the second tube, labelled S3, and mix thoroughly before the next transfer.

Repeat serial dilutions 4 more times thus creating the points of the standard curve (see Figure 5).

Sample Diluent serves as blank.

Figure 5

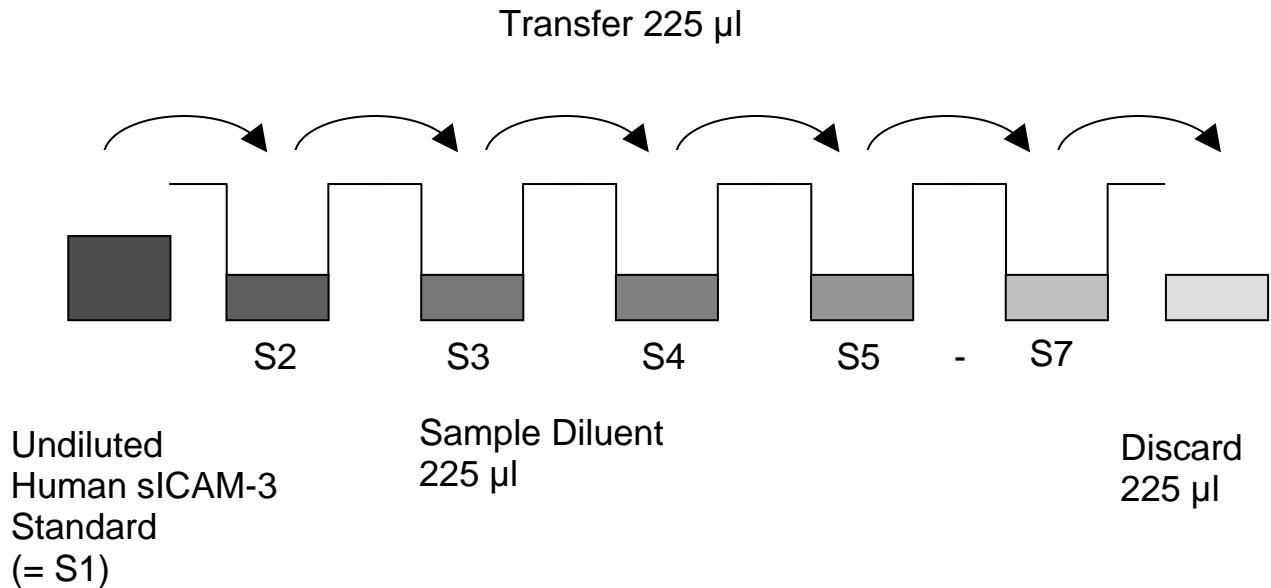

9.3 Addition of Colour-giving Reagents: Blue-Dye, Green-Dye

In order to help our customers to avoid any mistakes in pipetting the Platinum ELISAs, eBioscience offers a tool that helps to monitor the addition of even very small volumes of a solution to the reaction well by giving distinctive colours to each step of the ELISA procedure.

This procedure is optional, does not in any way interfere with the test results, and is designed to help the customer with the performance of the test, but can also be omitted, just following the instruction booklet.

Alternatively, the dye solutions from the stocks provided (**Blue-Dye, Green-Dye**) can be added to the reagents according to the following guidelines:

1. Diluent:

Before standard and sample dilution add the **Blue-Dye** at a dilution of 1:250 (see table below) to the appropriate diluent (1x) according to the test protocol. After addition of **Blue-Dye**, proceed according to the instruction booklet.

5 ml Sample Diluent	20 µl Blue-Dye
12 ml Sample Diluent	48 µl Blue-Dye
50 ml Sample Diluent	200 µl Blue-Dye

2. HRP-Conjugate:


Before dilution of the concentrated HRP-Conjugate add the **Green-Dye** at a dilution of 1:100 (see table below) to the Assay Buffer (1x) used for the final conjugate dilution. Proceed after addition of **Green-Dye** according to the instruction booklet: Preparation of HRP-Conjugate.

3 ml Assay Buffer (1x)	30 µl Green-Dye
6 ml Assay Buffer (1x)	60 µl Green-Dye

10 Test Protocol

- a. Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Each sample, standard, blank and optional control sample should be assayed in duplicate. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2°-8°C sealed tightly.
- b. Wash the microwell strips twice with approximately 400 µl **Wash Buffer** per well with thorough aspiration of microwell contents between washes. Allow the Wash Buffer to sit in the wells for about **10 – 15 seconds** before aspiration. Take care not to scratch the surface of the microwells.
After the last wash step, empty wells and tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing. Alternatively microwell strips can be placed upside down on a wet absorbent paper for not longer than 15 minutes. **Do not allow wells to dry.**
- c. **Standard dilution on the microwell plate** (Alternatively the standard dilution can be prepared in tubes - see 9.2.1):
Add 100 µl of Sample Diluent in duplicate to **standard wells** B1/2-G1/2, leaving A1/A2 empty. Pipette 200 µl of undiluted **standard** (concentration = 50.00 ng/ml) in duplicate into well A1 and A2 (see Table 1). Transfer 100 µl to wells B1 and B2. Mix the contents of wells B1 and B2 by repeated aspiration and ejection, and transfer 100 µl to wells C1 and C2, respectively. (see Figure 6). Take care not to scratch the inner surface of the microwells. Continue this procedure 4 times, creating two rows of human sICAM-3 standard dilutions ranging from 50.00 to 0.78 ng/ml. Discard 100 µl of the contents from the last microwells (G1, G2) used.

Figure 6

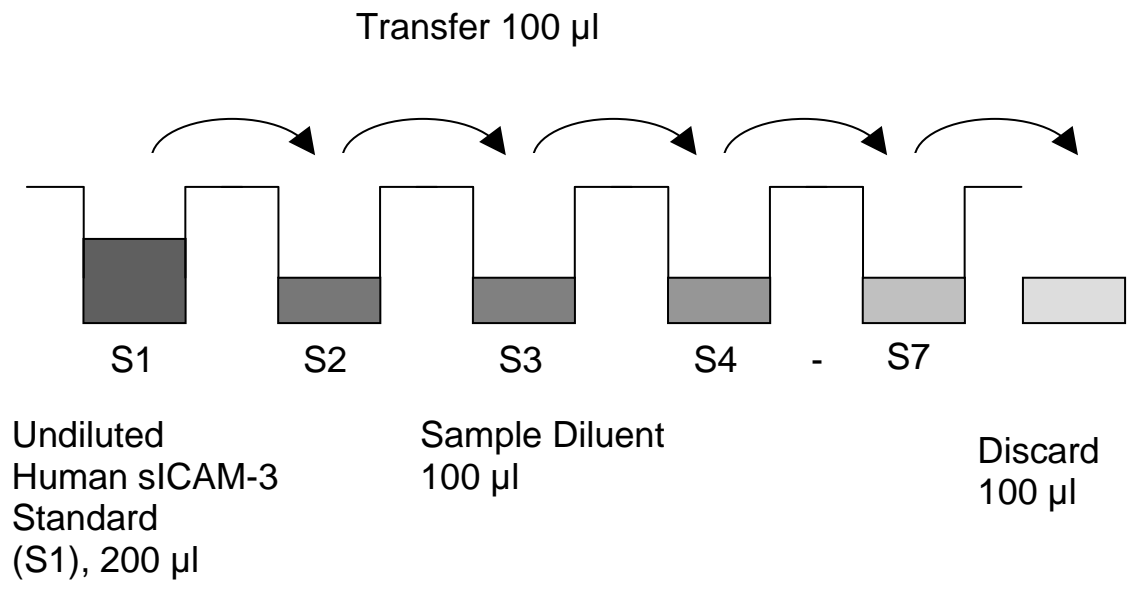

In case of an **external standard dilution** (see 9.2.1), pipette 100 µl of these standard dilutions (S1 - S7) in the standard wells according to Table 1.

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (50.00 ng/ml)	Standard 1 (50.00 ng/ml)	Sample 1	Sample 1
B	Standard 2 (25.00 ng/ml)	Standard 2 (25.00 ng/ml)	Sample 2	Sample 2
C	Standard 3 (12.50 ng/ml)	Standard 3 (12.50 ng/ml)	Sample 3	Sample 3
D	Standard 4 (6.25 ng/ml)	Standard 4 (6.25 ng/ml)	Sample 4	Sample 4
E	Standard 5 (3.13 ng/ml)	Standard 5 (3.13 ng/ml)	Sample 5	Sample 5
F	Standard 6 (1.56 ng/ml)	Standard 6 (1.56 ng/ml)	Sample 6	Sample 6
G	Standard 7 (0.78 ng/ml)	Standard 7 (0.78 ng/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- d. Add 100 µl of **Sample Diluent** in duplicate to the **blank wells**.
- e. Add 80 µl of **Sample Diluent** to the **sample wells**.
- f. Add 20 µl of each **sample** in duplicate to the **sample wells**.
- g. Add 50 µl of **HRP-Conjugate** to all wells.
- h. Cover with an adhesive film and incubate at room temperature (18 to 25°C) for 2 hours, if available on a microplate shaker set at 400 rpm.
- i. Remove adhesive film and empty wells. **Wash** microwell strips 3 times according to point b. of the test protocol. Proceed immediately to the next step.
- j. Pipette 100 µl of **TMB Substrate Solution** to all wells.
- k. Incubate the microwell strips at room temperature (18° to 25°C) for about 10 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see next point of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. Alternatively the colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.9 – 0.95.

- l. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.

m. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean value.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the human sICAM-3 concentration on the abscissa. Draw a best fit curve through the points of the graph (a 5-parameter curve fit is recommended).
- To determine the concentration of circulating human sICAM-3 for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding human sICAM-3 concentration.
- **If instructions in this protocol have been followed samples have been diluted 1:5 (20 µl sample + 80 µl Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 5).**
- **Calculation of samples with a concentration exceeding standard 1 will result in incorrect, low human sICAM-3 levels (Hook Effect). Such samples require further external predilution according to expected human sICAM-3 values with Sample Diluent in order to precisely quantitate the actual human sICAM-3 level.**
- It is suggested that each testing facility establishes a control sample of known human sICAM-3 concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 7. This curve cannot be used to derive test results. Each laboratory must prepare a standard curve for each group of microwell strips assayed.

Figure 7

Representative standard curve for human sICAM-3 ELISA. Human sICAM-3 was diluted in serial 2-fold steps in Sample Diluent. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

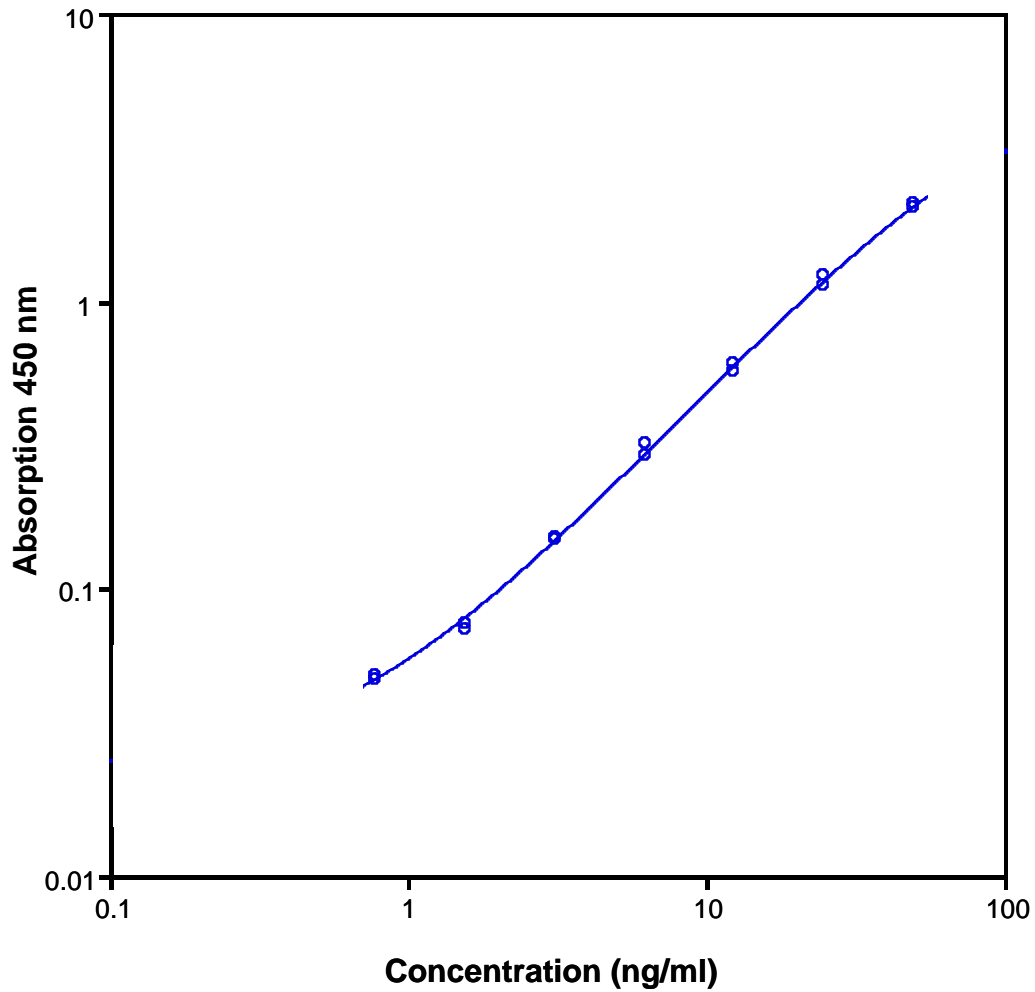


Table 2

Typical data using the human sICAM-3 ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	Human sICAM-3 Concentration (ng/ml)	O.D. at 450 nm	Mean O.D. at 450 nm	C.V. (%)
1	50.00	2.098 2.166	2.132	2.3
2	25.00	1.226 1.138	1.182	5.3
3	12.50	0.607 0.564	0.586	5.2
4	6.25	0.318 0.291	0.305	6.2
5	3.13	0.151 0.147	0.149	1.2
6	1.60	0.072 0.075	0.074	2.9
7	0.78	0.050 0.048	0.049	2.9
Blank	0	0.020 0.011	0.016	29.0

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of radioimmunotherapy has significantly increased the number of patients with human anti-mouse IgG antibodies (HAMA). HAMA may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analysed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the sample.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of human sICAM-3 defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 0.38 ng/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of human sICAM-3. 2 standard curves were run on each plate. Data below show the mean human sICAM-3 concentration and the coefficient of variation for each sample (see Table 3). The calculated overall intra-assay coefficient of variation was 2.5%.

Table 3

The mean human sICAM-3 concentration and the coefficient of variation for each sample

Sample	Experiment	Mean Human sICAM-3 Concentration (ng/ml)	Coefficient of Variation (%)
1	1	95.3	4.2
	2	98.2	2.2
	3	90.8	10.4
2	1	74.4	0.7
	2	74.1	0.6
	3	71.3	2.7
3	1	69.9	2.6
	2	69.1	2.4
	3	68.7	2.1
4	1	54.7	2.8
	2	56.9	4.3
	3	56.6	2.0
5	1	46.2	2.5
	2	45.2	2.2
	3	45.1	1.3
6	1	36.6	0.8
	2	37.5	4.3
	3	36.7	3.1
7	1	26.4	0.6
	2	27.3	1.5
	3	28.2	2.4
8	1	59.9	1.7
	2	59.5	1.1
	3	58.7	1.1

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of human sICAM-3. 2 standard curves were run on each plate. Data below show the mean human sICAM-3 concentration and the coefficient of variation calculated on 18 determinations of each sample (see Table 4). The calculated overall inter-assay coefficient of variation was 2.0%.

Table 4

The mean human sICAM-3 concentration and the coefficient of variation of each sample

Sample	Mean Human sICAM-3 Concentration (ng/ml)	Coefficient of Variation (%)
1	94.8	4.0
2	73.3	2.3
3	69.3	0.9
4	56.1	2.2
5	45.5	1.4
6	36.9	1.4
7	27.3	3.1
8	59.3	1.0

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of human sICAM-3 into different pooled normal human serum samples. Recoveries were determined in 3 independent experiments with 4 replicates each.

The amount of endogenous human sICAM-3 in unspiked serum was subtracted from the spike values.

The recovery ranged from 98% to 112% with an overall mean recovery of 106%.

13.4 Dilution Parallelism

4 serum samples with different levels of human sICAM-3 were analysed at serial 2 fold dilutions with 4 replicates each.

The recovery ranged from 87% to 104% with an overall recovery of 96% (see Table 5).

Table 5

Sample	Dilution	Expected Human sICAM-3 Concentration (ng/ml)	Observed Human sICAM-3 Concentration (ng/ml)	Recovery of Expected Human sICAM-3 Concentration (%)
1	1:5	--	59.6	--
	1:10	29.8	30.3	102
	1:20	14.9	15.5	104
2	1:5	--	60.6	--
	1:10	30.3	29.4	97
	1:20	15.2	14.1	93
3	1:5	--	78.8	--
	1:10	39.4	37.5	95
	1:20	19.7	18.1	92
4	1:5	--	102.7	--
	1:10	51.3	48.1	94
	1:20	25.7	22.3	87

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20°C and thawed 5 times, and the human sICAM-3 levels determined. There was no significant loss of human sICAM-3 immunoreactivity detected by freezing and thawing.

13.5.2 Storage Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20°C, 2-8°C, room temperature (RT) and at 37°C, and the human sICAM-3 level determined after 24, 48 and 96 h. There was no significant loss of human sICAM-3 immunoreactivity detected during storage under above conditions.

13.6 Comparison of Serum and Plasma

From 3 individuals, serum as well as EDTA and heparin plasma obtained at the same time point, was evaluated. Human sICAM-3 levels were not significantly different and therefore all these blood preparations are suitable for human sICAM-3 determinations.

13.7 Specificity

The interference of circulating factors of the immune system was evaluated by spiking these proteins at physiologically relevant concentrations into a human sICAM-3 positive serum. There was no crossreactivity detected, notably not with human sICAM-1.

13.8 Expected Values

Panels of 40 serum as well as EDTA and heparin plasma samples from randomly selected apparently healthy donors (males and females) were tested for human sICAM-3. The levels measured may vary with the sample collection used.

For detected human sICAM-3 levels see Table 6.

Table 6

Sample Matrix	Number of Samples Evaluated	Range (ng/ml)	Mean (ng/ml)	Standard Deviation (ng/ml)
Serum	40	26.6 – 112.7	71.0	21.4
Plasma (EDTA)	40	34.9 – 160.6	82.0	27.6
Plasma (Heparin)	40	22.9 – 114.2	47.8	17.6

14 Ordering Information

North America

Technical Support:

Research Products:
888.810.6168
858.642.2058
tech@eBioscience.com

Clinical Products:
877.726.8559
858.642.2058
tech@eBioscience.com

Customer Service:

888.999.1371
858.642.2058
info@eBioscience.com

Fax:

858.642.2046

Europe/International*

Technical Support:

+43 1 796 40 40-120
tech@eBioscience.com

Customer Service:

+43 1 796 40 40-304
info@eBioscience.com

Fax:

+43 1 796 40 40-400



Bender MedSystems GmbH
Campus Vienna Biocenter 2
1030 Vienna, Austria
www.eBioscience.com

* Customers outside North America and Europe may contact their eBioscience distributor listed on our website at www.eBioscience.com/distributors.

15 Reagent Preparation Summary

15.1 Wash Buffer (1x)

Add **Wash Buffer Concentrate** 20x (50 ml) to 950 ml distilled water.

Number of Strips	Wash Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

16 Test Protocol Summary

1. Determine the number of microwell strips required.
2. Wash microwell strips twice with Wash Buffer.
3. Standard dilution on the microwell plate: Add 100 µl Sample Diluent, in duplicate, to all standard wells leaving the first wells empty. Pipette 200 µl standard into the first wells and create standard dilutions by transferring 100 µl from well to well. Discard 100 µl from the last wells.
Alternatively external standard dilution in tubes (see 9.2.1): Pipette 100 µl of these standard dilutions in the microwell strips.
4. Add 100 µl Sample Diluent, in duplicate, to the blank wells.
5. Add 80 µl Sample Diluent to sample wells.
6. Add 20 µl sample in duplicate, to designated sample wells.
7. Add 50 µl HRP-Conjugate to all wells.
8. Cover microwell strips and incubate 2 hours at room temperature (18° to 25°C).
9. Empty and wash microwell strips 3 times with Wash Buffer.
10. Add 100 µl of TMB Substrate Solution to all wells.
11. Incubate the microwell strips for about 10 minutes at room temperature (18° to 25°C).
12. Add 100 µl Stop Solution to all wells.
13. Blank microwell reader and measure colour intensity at 450 nm.

Note: If instructions in this protocol have been followed samples have been diluted 1:5 (20 µl sample + 80 µl Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 5).