


Technical Data Sheet

Mouse IL-12p70 Flex Set

Product Information


Material Number:	558303
Size:	100 tests
Bead Position:	D7
Reactivity:	QC Testing: Mouse
Assay Range:	10 - 2,500 pg/mL

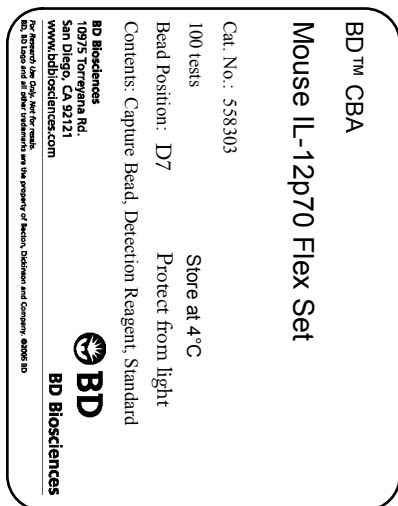

Component Description: Mouse IL-12p70 Standard
 Component Mat. No: 51-9003530
 Component Storage Buffer: Lyophilized in an aqueous buffered solution containing BSA and ProClin™ 150.

Component Description: Mouse IL-12p70 PE Detection Reagent
 Component Mat. No: 51-9004155
 Component Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Component Description: Mouse IL-12p70 Capture Bead D7
 Component Mat. No: 51-9005240
 Component Storage Buffer: Aqueous buffered solution containing fetal bovine serum and ≤0.09% sodium azide.

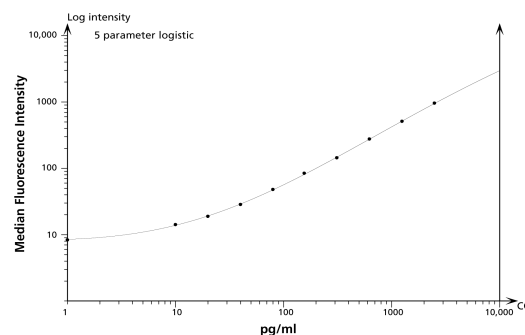

Example BD CBA Mouse IL-12p70 Flex Set standard curve. Data acquired on a BD FACSAArray bioanalyzer and analyzed using FCAP Array Software.

Description

The BD™ CBA Mouse IL-12p70 Flex Set is a bead-based immunoassay capable of measuring mouse interleukin-12p70 (IL-12p70) in serum and cell culture supernatant samples. Mouse reactivity was determined by testing samples with the BD CBA Mouse IL-12p70 Flex Set. The biology and function of IL-12p70 has been extensively reviewed in the literature. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266 or 558267).

Preparation and Storage

This BD™ CBA Flex Set contains one vial each of Capture Bead and PE Detection Reagent and two vials of Standard. The Capture Bead and PE Detection Reagent components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. The Standard component is lyophilized and should be transferred to a 15 mL polypropylene tube for reconstitution. When reconstituted in 4.0 mL Assay Diluent, the standard has a protein concentration of 2,500 pg/mL. Discard unused reconstituted standard, do not store or reuse. Store lyophilized standard and other components at 4°C. Protect Capture Beads and the PE Detection Reagent from prolonged exposure to light.

Application Notes

Recommended Assay Procedure: The BD CBA Mouse IL-12p70 Flex Set must be used in conjunction with a BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266, 100 tests, or 558267, 500 tests), a flow cytometer, and FCAP Array™ Software. Detailed instructions on the use of this product can be found in the manual for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Mouse IL-12p70 Flex Set will be 2,500 pg/mL. An example standard curve is shown in Figure 1.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



The BD CBA Mouse IL-12p70 Flex Set cannot be used in the same assay well with the following BD CBA Mouse Soluble Protein Flex Set reagents:

Flex Set	Bead Position	Catalog Number
BD CBA Mouse IL-12/IL-23p40 Flex Set	D7	560151

The BD CBA Mouse IL-12p70 Flex Set should not be used in the same assay well with any non-BD CBA Mouse Soluble Protein Flex Set reagents (such as BD CBA Human or Rat Soluble Protein or Cell Signaling Flex Sets). For an updated assay compatibility chart for the BD CBA Mouse Soluble Protein Flex Sets, please refer to the BD CBA Flex Set System homepage at <http://www.bdbiosciences.com/cbasetup>.

Performance

Limit of Detection: The theoretical limit of detection is 1.9 pg/mL and was determined by evaluating the estimated result of the average MFI of the negative control (0 pg/mL, $n=30$) + 2 standard deviations.

Specificity		Inter-Assay Reproducibility			Intra-Assay Reproducibility		
		Mean (pg/ml)	Standard Deviation	%CV	Mean (pg/ml)	Standard Deviation	%CV
Mouse IL-12p70	Sample 1	39.3	3.4	9%	38.8	2.7	7%
	Sample 2	158.2	5.3	3%	156.6	7.3	5%
	Sample 3	623.8	14.9	2%	613.8	16.7	3%

Reproducibility: The inter-assay and intra-assay reproducibility were determined for the BD CBA Mouse IL-12p70 Flex Set by evaluating ten replicates of three different sample levels (intra-assay) and two replicates of three different sample levels from four separate experiments (inter-assay).

Sample Dilution	Cell Culture Supernatant		Serum	
	Average % Recovery	Range	Average % Recovery	Range
Mouse IL-12p70	97%	95 - 100%	71%	65 - 74%

Recovery: Cell culture supernatant and pooled mouse serum were spiked with three different levels of protein. The spiked samples were assayed and the results were compared with expected values. Serum samples were diluted 1:4 before the protein was spiked into each. Serum used was commercially available pooled mouse serum.

Sample Dilution	Cell Culture Supernatant		Serum	
	Detected (pg/ml)	% of Expected	Detected (pg/ml)	% of Expected
Spiked sample	600.3	100%	503.2	100%
1 : 2	302.4	101%	279.4	111%
1 : 4	155.9	104%	144.5	115%

Linearity: Cell culture supernatant and 1:4 diluted pooled mouse serum were spiked with protein and serially diluted. The diluted samples were assayed and the results were compared with the original spiked sample.

Product Notices

1. ProClin is a trademark of Rohm and Haas Company.
2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 800.268.5430 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

