Technical Data Sheet

Alexa Fluor® 647 Mouse Anti-Human CD193

Product Information


Material Number: 561745 Alternate Name: CCR3 25 μg Size **Concentration:** 0.2 mg/ml Clone: 5E8

Mouse (C57BL/6) IgG2b, κ **Isotype:** Reactivity: QC Testing: Human

Aqueous buffered solution containing ≤0.09% sodium azide. Storage Buffer:

Description

Antibody 5E8 recognizes an epitope in CCR3. CCR3 is a G protein-linked, 7 transmembrane, chemokine receptor expressed on a variety of hematopoietic cells. Similar to CCR5 and CXCR4, CCR3 can be a co-receptor for HIV-1. It is primarily found on eosinophils and basophils during atopic conditions, dermatitis, allergic rhinitis, conjunctivitis and bronchial asthma. Chemokines like RANTES, Eotaxin, MCP-3, MIP1α have been reported to act as ligands for CCR3 and stimulate CCR3+ cells. Eotaxin stimulates Th2 cells expressing CCR3. Other studies describe HIV-1 specific T cell cytotoxicity can be mediated by RANTES and Eotaxin through CCR3. CCR3 expressed on dendritic cells may have a biological role on cell-cell interaction during antigen presentation. CCR3 has been clustered as CD193 in the HLDA VIIIth workshop.

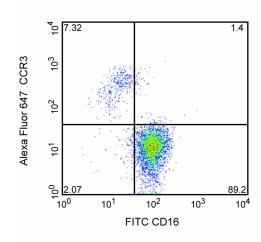

Profile of anti-CCR3 (5E8) reactivity on peripheral blood granulocytes (eosinophils) analyzed by flow cytometry

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Flow cytometry	Routinely Tested	

Suggested Companion Products

Catalog Number	Name	Size	Clone	
557903	Alexa Fluor® 647 Mouse IgG2b, κ Isotype Control	100 tests	27-35	
560996	FITC Mouse Anti-Human CD16	25 tests	3G8	
554656	Stain Buffer (FBS)	500 ml	(none)	

Product Notices

- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).

BD Biosciences

bdbiosciences.com


Asia Pacific Europe Japan 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD

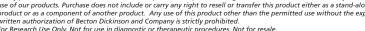



- 4. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
- 5. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- 6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 7. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 8. An isotype control should be used at the same concentration as the antibody of interest.

References

Daugherty BL, Siciliano SJ, DeMartino JA, Malkowitz L, Sirotina A, Springer MS. Cloning, expression, and characterization of the human eosinophil eotaxin receptor. *J Exp Med.* 1996; 83(5):2349-2354. (Biology)

Ghorpade A, Xia MQ, Hyman BT, et al. Role of the beta-chemokine receptors CCR3 and CCR5 in human immunodeficiency virus type 1 infection of monocytes and microglia. J Virol. 1998; 72(4):3351-3361. (Biology)

Hadida F, Vieillard V, Autran B, Clark-Lewis I, Baggiolini M, Debre P. HIV-specific T cell cytotoxicity mediated by RANTES via the chemokine receptor CCR3. *J Exp Med.* 1998; 188(3):609-614. (Biology)

Sallusto F, Mackay CR, Lanzavecchia A. Selective expression of the eotaxin receptor CCR3 by human T helper 2 cells. Science. 1997; 277(5334):2005-2007. (Biology)

Sato K, Kawasaki H, Nagayama H, et al. CC chemokine receptors, CCR-1 and CCR-3, are potentially involved in antigen-presenting cell function of human peripheral blood monocyte-derived dendritic cells. *Blood.* 1999; 93(1):34-42. (Biology)

Zimmermann N, Daugherty BL, Stark JM, Rothenberg ME. Molecular analysis of CCR-3 events in eosinophilic cells. J Immunol. 2000; 164(2):1055-1064. (Biology)

561745 Rev. 1 Page 2 of 2