Technical Data Sheet

APC Mouse Anti-Human CD90

Product Information

Material Number: 561971

Alternate Name: THY1; Thy-1 antigen; Thy-1 membrane glycoprotein

Size 0.2 mg/ml Concentration: Clone: 5E10 Isotype: Mouse IgG1, κ

Tested in Development: Baboon, Rhesus, Cynomolgus, Pig, and Dog.

Workshop:


Storage Buffer: Aqueous buffered solution containing ≤0.09% sodium azide.

Description

Reactivity:

The 5E10 monoclonal antibody specifically binds to human CD90. CD90 is a 25-35 kDa molecule expressed on 1-4% of human fetal liver cells, cord blood cells, and bone marrow cells. Anti-CD90 reacts with a subset of immature, CD34+ cells and a distinct subset of mature CD34- cells that are CD3+CD4+. The CD90+CD34+ population is highly enriched for cells capable of long-term culture. Anti-CD90 is useful for enriching high proliferative potential colony-forming cells (HIPP-CFC) which are primative progenitor cells.

QC Testing: Human

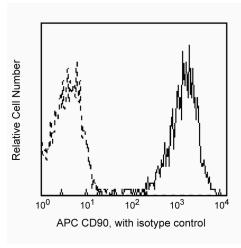

Profile of HEL cell line analyzed by flow cytometry

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number Name Size Clone 555751 APC Mouse IgG1, κ Isotype Control 100 tests MOPC-21

Product Notices

- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.

BD Biosciences

bdbiosciences.com

Asia Pacific Japar 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



561971 Rev. 1

6. An isotype control should be used at the same concentration as the antibody of interest.

References

Baum CM, Weissman IL, Tsukamoto AS, Buckle AM, Peault B. Isolation of a candidate human hematopoietic stem-cell population. *Proc Natl Acad Sci U S A*. 1992; 89(7):2804-2808. (Biology)

Craig W, Kay R, Cutler RL, Lansdorp PM. Expression of Thy-1 on human hematopoietic progenitor cells. J Exp Med. 1993; 177(5):1331-1342. (Biology)

Knapp W, Dorken B, Rieber EP, et al, ed. Leucocyte Typing IV. New York: Oxford University Press; 1989:1-1208. (Biology)

Lansdorp PM, Thomas TE. AP Gee, ed. Bone Marrow Processing and Purging. Boca Raton FL: CRC Press; 1991. (Biology)

Schlossman S, Boumell L, et al, ed. *Leucocyte Typing V*. New York: Oxford University Press; 1995. (Clone-specific)

561971 Rev. 1 Page 2 of 2