

Technical Data Sheet


FITC Rat Anti-Mouse Vα2 TCR

Product Information

Material Number:	562085
Size:	25 µg
Concentration:	0.5 mg/ml
Clone:	B20.1
Immunogen:	Soluble αβ TCR from mouse cytotoxic T-cell clone KB5-C20
Isotype:	Rat (LOU) IgG2a, λ
Reactivity:	QC Testing: Mouse
Storage Buffer:	Aqueous buffered solution containing ≤0.09% sodium azide.

Description

The B20.1 monoclonal antibody specifically binds to most members of the Vα2 T-cell Receptor (TCR) subfamily in mice having the *a*, *b*, and *c* haplotypes of the *Tcrb* gene complex. B20.1 antibody may crossreact with Vδ8 TCR, which shares >90% sequence homology with Vα2 TCR. Levels of B20.1+ T cells appear to be influenced by Vα haplotypes. Moreover, the frequencies of Vα2+ CD8+ and CD4+ T cells are influenced by H-2 haplotypes.

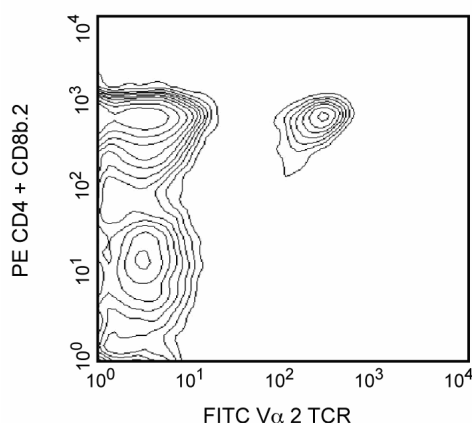

Two-color analysis of the expression of Vα2 TCR on peripheral lymphocytes. BALB/c lymph node cells were incubated simultaneously with FITC Rat anti-Mouse Vα2 TCR, PE Rat anti-Mouse CD4 (Cat. No. 561829/557308), and PE Rat anti-Mouse CD8b.2 (Cat. No. 553041) monoclonal antibodies. Flow cytometry was performed on a BD FACScan™ flow cytometer (BDIS, San Jose, CA).

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed.

Application Notes

Application

Flow cytometry

Routinely Tested

Recommended Assay Procedure:

For flow cytometry of cell suspensions from peripheral lymphoid tissues, it is recommended that multicolor staining be performed to distinguish T lymphocytes from non-T cells.

Suggested Companion Products

Catalog Number	Name	Size	Clone
557077	FITC Rat IgG2a, λ Isotype Control	0.25 mg	B39-4
554656	Stain Buffer (FBS)	500 ml	(none)
561829	PE Rat Anti-Mouse CD4	25 µg	GK1.5
557308	PE Rat Anti-Mouse CD4	0.1 mg	GK1.5
553041	PE Rat Anti-Mouse CD8b.2	0.1 mg	53-5.8

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.268.5430	32.53.720.550	0120.8555.90	65.6861.0633	0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
5. An isotype control should be used at the same concentration as the antibody of interest.

References

- Gregoire C, Rebai N, Schweisguth F, et al. Engineered secreted T-cell receptor alpha beta heterodimers. *Proc Natl Acad Sci U S A*. 1991; 88(18):8077-8081. (Clone-specific)
- Pircher H, Rebai N, Groettrup M, et al. Preferential positive selection of V alpha 2+ CD8+ T cells in mouse strains expressing both H-2k and T cell receptor V alpha a haplotypes: determination with a V alpha 2-specific monoclonal antibody. *Eur J Immunol*. 1992; 22(2):399-404. (Immunogen)
- Tomonari K, Fairchild S, Rosenwasser OA. Influence of viral superantigens on V beta- and V alpha-specific positive and negative selection. *Immunol Rev*. 1993; 131:131-168. (Biology)