## **Technical Data Sheet**

# V500 Rat anti-Mouse CD8a

## **Product Information**

Material Number:
Alternate Name:
Size:
Concentration:
Clone:
Immunogen:
Isotype:
Reactivity:
Storage Buffer:

560776 Ly-2, Lyt-2 0.1 mg 0.2 mg/ml 53-6.7 Mouse thymus / spleen Cells Rat (LOU) IgG2a, κ QC Testing: Mouse Aqueous buffered solution containing protein stabilizer, glycerol and ≤0.09% sodium azide.

## Description

The 53-6.7 antibody has been reported to react with the 38 kDa  $\alpha$  and 34 kDa  $\alpha'$  chains of the CD8 differentiation antigen (Ly-2 or Lyt-2) of all mouse strains tested. The CD8  $\alpha$  and  $\alpha'$  chains (CD8a) form heterodimers with the CD8  $\beta$  chain (CD8b, Ly-3, or Lyt-3) on the surface of most thymocytes. A subpopulation of mature T lymphocytes (i.e., MHC class I-restricted T cells, including most T suppressor/cytotoxic cells) expresses almost exclusively the CD8  $\alpha\beta$  heterodimer (the  $\alpha'$  chain is absent). Subsets of  $\gamma\delta$  TCR-bearing T cells, intestinal intrapithelial lymphocytes, and dendritic cells express CD8a without CD8b. It has been suggested that the expression of the CD8a/CD8b heterodimer is restricted to T lymphocytes which matured in the thymus or in an extrathymic environment that had been influenced by thymus-initiated neuroendocrine signals. CD8 is an antigen coreceptor on the T-cell surface which interacts with MHC class I molecules on antigen-presenting cells or epithelial cells. It participates in T-cell activation through its association with the T-cell receptor complex and protein tyrosine kinase lck (p56 [lck]). The CD8  $\alpha$  and  $\alpha'$  chains arise from alternatively spliced messengers of a single *CD8a* gene. The longer  $\alpha$  form associates with p56 [lck] via a CXCP motif in its cytoplasmic domain, which it shares with CD4, but not with CD8b. The truncated  $\alpha'$  chain is unable to associate with p56 [lck], and it may function to attenuate the CD8-mediated costimulatory signal during intrathymic T-cell maturation. In vivo and in vitro treatment with 53-6.7 mAb has reportedly been effective at depleting CD8+ peripheral T lymphocytes in the Egyptian toad, *Bufo regularis*.

The antibody is conjugated to BD Horizon<sup>TM</sup> V500, which has been developed for use in multicolor flow cytometry experiments and is available exclusively from BD Biosciences. It is excited by the Violet laser with an Ex max of 415 nm and Em Max at 500 nm. BD Horizon V500 conjugates emit at a similar wavelength to Amcyan yet exhibit reduced spillover into the FITC channel. For more information on BD Horizon V500, visit bdbiosciences.com/colors.

When compensating dyes in this spectral range (such as Horizon<sup>™</sup> V500 and AmCyan), the most accurate compensation can be obtained using single stained cellular controls. Due to spectral differences between cells and beads in this channel, using BD CompBeads can result in spillover errors for V500 and AmCyan reagents. Therefore, the use of BD CompBeads or BD CompBeads Plus to determine spillover values for these reagents is not recommended. Different V500 reagents (e.g. CD4 vs. CD45) can have slightly different fluorescence spillover therefore, it may also be necessary to use clone specific compensation controls when using these reagents.

#### **BD Biosciences**

560776 Rev. 1

 bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean


 877.232.8995
 888.268.5430
 32.53.720.550
 0120.8555.90
 65.8661.0633
 0800.771.7157

 For country-specific contact information, visit
 bdbiosciences.com/how\_to\_order/

 Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

 For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
 B), BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



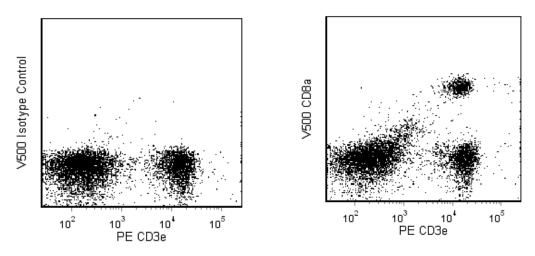

Analysis of CD8a on mouse splenocytes. Splenocytes from BALB/c mice were stained simultaneously with BD Horizon<sup>™</sup> V500 Rat anti-Mouse CD8a (right panel) or BD Horizon<sup>™</sup> V500 Rat IgG2a, κ Isotype Control (clone R35-95, Cat. No. 560786, left panel), and PE Hamster Anti-Mouse CD3e (clone 145-2C11, Cat. No. 553063/553064). The contour plots were derived from gated events based on light scattering characteristics of splenocytes. Flow cytometry was performed on a BD FACSCanto<sup>™</sup> II flow cytometry system.

### **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon<sup>TM</sup> V500 under optimum conditions, and unreacted BD Horizon<sup>TM</sup> V500 was removed.

## **Application Notes**

Application					
Flow cytometry		Routinely Tested			
Suggested Companion Products					
Catalog Number	Name		Size	Clone	
560786	V500 Rat IgG2a, к Isotype Control		0.1 mg	R35-95	
553063	PE Hamster Anti-Mouse CD3e		0.1 mg	145-2C11	

#### **Product Notices**

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.

- 2. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 3. BD Horizon<sup>™</sup> V500 has a maximum absorption of 415 nm and maximum emission of 500 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 5. An isotype control should be used at the same concentration as the antibody of interest.
- 6. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

#### References

Bierer BE, Sleckman BP, Ratnofsky SE, Burakoff SJ. The biologic roles of CD2, CD4, and CD8 in T-cell activation. *Annu Rev Immunol.* 1989; 7:579-599. (Biology) Fujiura Y, Kawaguchi M, Kondo Y, et al. Development of CD8 alpha alpha+ intestinal intraepithelial T cells in beta 2-microglobulin- and/or TAP1-deficient mice. *J Immunol.* 1996; 156(8):2710-2715. (Biology)

Hathcock KS. T cell depletion by cytotoxic elimination. In: Coligan JE, Kruisbeek AM, Margulies DH, Shevach EM, Strober W, ed. *Current Protocols in Immunology*. New York: John Wiley and Sons; 1991:3.4.1-3.4.3. (Biology)

Janeway CA Jr. The T cell receptor as a multicomponent signalling machine: CD4/CD8 coreceptors and CD45 in T cell activation. Annu Rev Immunol. 1992; 10:645-674. (Biology)

Kruisbeek AM, Shevach EM. Proliferative assays for T cell function. In: Coligan J, Kruisbeek AM, Margulies D, Shevach EM, Strober W, ed. Current Protocols in Immunology. New York: John Wiley and Sons; 1991:3.12.1-3.12.14. (Biology)

Ledbetter JA, Herzenberg LA. Xenogeneic monoclonal antibodies to mouse lymphoid differentiation antigens. *Immunol Rev.* 1979; 47:63-90. (Immunogen) Ledbetter JA, Rouse RV, Micklem HS, Herzenberg LA. T cell subsets defined by expression of Lyt-1,2,3 and Thy-1 antigens. Two-parameter immunofluorescence and cytotoxicity analysis with monoclonal antibodies modifies current views. *J Exp Med.* 1980; 152(2):280-295. (Biology) Ledbetter JA, Seama WE, Tsu TJ, Herzenberg LA, Lyt-2 and Ivt-3 antienes are on two different nolveneting slinked by disulfide bonds. Relationship of

Ledbetter JA, Seaman WE, Tsu TT, Herzenberg LA. Lyt-2 and lyt-3 antigens are on two different polypeptide subunits linked by disulfide bonds. Relationship of subunits to T cell cytolytic activity. J Exp Med. 1981; 153(6):1503-1516. (Biology)

LeFrancois L. Extrathymic differentiation of intraepithelial lymphocytes: generation of a separate and unequal T-cell repertoire. *Immunol Today.* 1991; 12(12):436-438. (Biology)

Leishman AJ, Naidenko OV, Attinger A, et al. T cell responses modulated through interaction between CD8alphaalpha and the nonclassical MHC class I molecule, TL. Science. 2001; 294(5548):1848-1849. (Biology)

MacDonald HR, Schreyer M, Howe RC, Bron C. Selective expression of CD8 alpha (Ly-2) subunit on activated thymic gamma/delta cells. Eur J Immunol. 1990; 20(4):927-930. (Biology)

Mitnacht R, Bischof A, Torres-Nagel N, Hunig T. Opposite CD4/CD8 lineage decisions of CD4+8+ mouse and rat thymocytes to equivalent triggering signals: correlation with thymic expression of a truncated CD8 alpha chain in mice but not rats. *J Immunol.* 1998; 160(2):700-707. (Biology)

Murosaki S, Yoshikai Y, Ishida A, et al. Failure of T cell receptor V beta negative selection in murine intestinal intra-epithelial lymphocytes. *Int Immunol.* 1991; 3(10):1005-1013. (Biology)

Nakayama K, Nakayama K, Negishi I, et al. Requirement for CD8 beta chain in positive selection of CD8-lineage T cells. Science. 1994; 263(5150):1131-1133. (Biology)

Negm HI, Mansour MH, Saad AH, Abdel Halim RS. Structural characterization of an Lyt-2/3 homolog expressed on Bufo regularis lymphocytes. Comp Biochem Physiol B Biochem Mol Biol. 1996; 113(1):79-87. (Biology)

O'Rourke AM, Mescher MF. The roles of CD8 in cytotoxic T lymphocyte function. Immunol Today. 1993; 14(4):183-188. (Biology)

Rocha B, Vassalli P, Guy-Grand D. The extrathymic T-cell development pathway. Immunol Today. 1992; 14(3):140-141. (Biology)

Sydora BC, Brossay L, Hagenbaugh A, Kronenberg M, Cheroutre H. TAP-independent selection of CD8+ intestinal intraepithelial lymphocytes. *J Immunol.* 1996; 156(11):4209-4216. (Biology)

Traver D, Akashi K, Manz M, et al. Development of CD8alpha-positive dendritic cells from a common myeloid progenitor. Science. 2000; 290(5499):2152-2154. (Biology)

van Ewijk W, van Soest PL, van den Engh GJ. Fluorescence analysis and anatomic distribution of mouse T lymphocyte subsets defined by monoclonal antibodies to the antigens Thy-1, Lyt-1, Lyt-2, and T-200. *J Immunol.* 1981; 127(6):2594-2604. (Biology)

Walker ID, Murray BJ, Hogarth PM, Kelso A, McKenzie IF. Comparison of thymic and peripheral T cell Ly-2/3 antigens. *Eur J Immunol.* 1984; 14(10):906-910. (Biology)

Wang J, Klein JR. Thymus-neuroendocrine interactions in extrathymic T cell development. Science. 1994; 265(5180):1860-1862. (Biology)

Zamoyska R. The CD8 coreceptor revisited: one chain good, two chains better. Immunity. 1994; 1(4):243-246. (Biology)

Zamoyska R, Derham P, Gorman SD, et al. Inability of CD8 alpha' polypeptides to associate with p56lck correlates with impaired function in vitro and lack of expression in vivo. *Nature*. 1989; 342(6247):278-281. (Biology)

Zamoyska R, Vollmer AC, Sizer KC, Liaw CW, Parnes JR. Two Lyt-2 polypeptides arise from a single gene by alternative splicing patterns of mRNA. *Cell.* 1985; 43(1):153-163. (Biology)