Technical Data Sheet

PE-Cy[™]5 Mouse Anti-Human CD137

Product Information


Material Number: 551137 Alternate Name: 4-1BB 100 tests Size $20~\mu l$ Vol. per Test: 4B4-1 Clone: Isotype: Mouse IgG1, κ Reactivity: QC Testing: Human

Workshop:

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Reacts with 4-1BB, a 30 kDa glycoprotein expressed on activated T cells. 4-1BB can be detected on Con A or PHA-stimulated peripheral blood T cells (CD4+ and CD8+) and on CEM cells (human T-cell leukemia) following 2 day PMA and ionomycin stimulation, but not on resting T cells. 4-1BB is reported to be a participant in T-cell activation, being responsible for rescuing T cells from activation-induced apoptosis, upregulation of Th1-type, T-helper cells, downregulation of Th2-type cytokine production, and induction of cell adhesion to fibronectin.

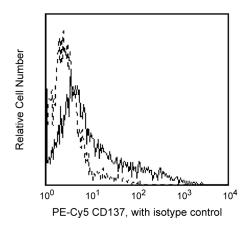

Profile of Con A-stimulated peripheral blood mononuclear cells (PBMC) analyzed by flow cytometry.

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with PE-Cy5 (formerly known as BD Cy-Chrome™) under optimum conditions, and unconjugated antibody and free PE-Cy5 were removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

	••	
Ap	plica	ation

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number Name Size Clone PE-CyTM5 Mouse IgG1 κ Isotype Control 555750 MOPC-21

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use $1 \times 10^{\circ}6$ cells in a 100- μ l experimental
- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

BD Biosciences

bdbiosciences.com

877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation cof any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



- 4. Cy is a trademark of Amersham Biosciences Limited. This conjugated product is sold under license to the following patents: US Patent Nos. 5,486,616; 5,569,587; 5,569,766; 5,627,027.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 7. This product is subject to proprietary rights of Amersham Biosciences Corp. and Carnegie Mellon University and made and sold under license from Amersham Biosciences Corp. This product is licensed for sale only for research. It is not licensed for any other use. If you require a commercial license to use this product and do not have one return this material, unopened to BD Biosciences, 10975 Torreyana Rd, San Diego, CA 92121 and any money paid for the material will be refunded.
- 8. PE-Cy5 is optimized for use with a single argon ion laser emitting 488-nm light. Because of the broad absorption spectrum of the PE-Cy5 tandem fluorochrome, extra care must be taken when using dual-laser cytometers which may directly excite both PE and Cy5TM.
- 9. PE-Cy5 is a tandem fluorochrome composed of R-phycoerythrin (PE), which is excited by the 488 nm light of an Argon ion laser and serves as an energy donor, coupled to the cyanine dye Cy5, which acts as an energy acceptor and fluoresces at 670 nm. BD Biosciences Pharmingen has maximized the fluorochrome energy transfer in PE-Cy5, thus maximizing its fluorescence emission intensity, minimizing residual emission from PE, and minimizing lot-to-lot variation.
- 10. PE-Cy5 tandem fluorochromes have been reported to bind some classes of human macrophages and granulocytes via Fc receptors, and PE has been reported to bind to mouse B lymphocytes via Fc receptors. Preincubation of mouse leukocytes with Mouse BD Fc Block™ purified anti-mouse CD16/CD32 mAb 2.4G2 can reduce the non-specific binding of PE-Cy5-conjugated reagents to mouse B cells. However, PE-Cy5 conjugated reagents should not be used to stain splenocytes of SJL, NOD, and MRL mice as B lymphocytes and/or other leukocytes have been reported to non-specifically stain regardless of the use of Mouse BD Fc Block™ (the CD72c complex has been implicated for PE-Cy5 binding in these strains). Reagents conjugated to PE, PerCP, PerCP-Cy5.5, APC, and APC-Cy7 tandem fluorochrome can be used on leukocytes from these mouse strains.

References

Garni-Wagner BA, Lee ZH, Kim YJ, Wilde C, Kang CY, Kwon BS. 4-1BB is expressed on CD45RAhiROhi transitional T cell in humans. *Cell Immunol.* 1996; 169(1):91-98. (Biology)

Hurtado JC, Kim SH, Pollok KE, Lee ZH, Kwon BS. Potential role of 4-1BB in T cell activation. Comparison with the costimulatory molecule CD28. *J Immunol.* 1995; 155(7):3360-3367. (Biology)

Hurtado JC, Kim YJ, Kwon BS. AAI Meeting, Abstract #841. San Francisco: 1997. (Biology)

Kim YJ, Kim SH, Kwon BS. AAI Meeting, Abstract #838. San Francisco: 1997. (Biology)

Zhou Z, Kim S, Hurtado J, et al. Characterization of human homologue of 4-1BB and its ligand. Immunol Lett. 1995; 45(1-2):67-73. (Biology)

551137 Rev. 9 Page 2 of 2