Technical Data Sheet

PE Mouse Anti-Human CD275

Product Information


Material Number: 552502 Alternate Name: B7-H2 100 tests Size Vol. per Test: 20 ul 2D3/B7-H2 Clone: Isotype: Mouse IgG2b, κ Reactivity: QC Testing: Human

Workshop:

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Reacts with a 302 amino acid cell surface protein termed B7-H2 (B-7 homologue 2). B7-H2 binds to a CD28-like receptor inducible costimulator (ICOS) and costimulates the proliferation and cytokine production of human T cells. It is weakly expressed on peripheral blood mononuclear cells, but it is expressed on monocyte-derived dendritic cells. This antibody could be useful in the study of dendritic cell development and/or costimulatory molecules.

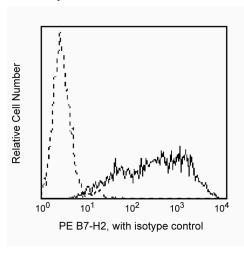

Profile of B7-H2 expression on CHO-B7-H2 cells analyzed by flow cytometry

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Flow cytometry	Routinely Tested	

Suggested Companion Products

Catalog Number	Name	Size	Clone
555743	PE Mouse IgG2b κ Isotype Control	100 tests	27-35

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-µl experimental
- 2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

BD Biosciences

bdbiosciences.com

877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

Bretscher PA. A two-step, two-signal model for the primary activation of precursor helper T cells. *Proc Natl Acad Sci U S A*. 1999; 96(1):185-190. (Biology)
Dong H, Zhu G, Tamada K, Chen L. B7-H1, a third member of the B7 family, co-stimulates T-cell proliferation and interleukin-10 secretion. *Nat Med.* 1999; 5(12):1365-1369. (Biology)
Wang S, Zhu G, Chapoval AI, et al. Costimulation of T cells by B7-H2, a B7-like molecule that binds ICOS. *Blood*. 2000; 96(8):2808-2813. (Biology)

552502 Rev. 4 Page 2 of 2