Technical Data Sheet

FITC Mouse Anti-Rat IgG2b

Product Information

 Material Number:
 553900

 Size:
 0.5 mg

 Concentration:
 0.5 mg/ml

 Clone:
 RG7/11.1

 Immunogen:
 Pooled rat IgG


 Isotype:
 Mouse (SJL) IgG2b, κ

 Reactivity:
 QC Testing: Rat

Storage Buffer: Aqueous buffered solution containing ≤0.09% sodium azide.

Description

The RG7/11.1 antibody reacts specifically with the Fc region of rat IgG2b. It does not react with other Ig isotypes.

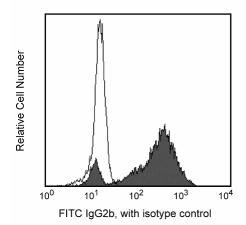

Detection of intracellular rat IgG2b in an antibody-secreting hybidoma cell line. Cells were fixed, permeabilized, and stained according to the method described in the recommended assay procedure using FITC-conjugated RG7/11.1 mAb (filled histogram) or the matched isotype control, FITC-conjugated 27-35 mAb (open histogram, Cat. No. 555057). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed. Store undiluted at 4° C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

-pppppppppppppppp.		
Flow cytometry	Routinely Tested	
Intracellular staining (flow cytometry)	Routinely Tested	

Recommended Assay Procedure:

RG7 antibody is effective for detection of cell-surface or intracellular Ig by immunofluorescent staining with flow cytometric analysis. FITC-conjugated RG7/11.1 mAb may be used as a primary or secondary reagent in immunofluorescent staining.

IMMUNOFLUORESCENT STAINING OF INTRACELLULAR IMMUNOGLOBULIN (Ig) PROTOCOL

- 1. Prepare a single-cell suspension and determine cell number.
- 2. Suspend cells in staining buffer (PBS + 2% FBS + 0.1% Sodium Azide) at $2 \times 10e7$ cells/ml and transfer to U-bottom microwell plates in 50 μ l/well for immunofluorescent staining.

Note: The BD Pharmingen™ Stain Buffer with FBS (Cat. No. 554656) is effective for use as a staining buffer in this protocol.

- 3. Block Fc γ receptors by adding 0.2 μ g of purified 2.4G2 antibody (Mouse BD Fc BlockTM purified anti-mouse CD16/CD32 mAb 2.4G2) (Cat. No. 553141/553142) in 50 μ l of staining buffer to each well.
- 4. Incubate 5 minutes on ice.

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 888.259.0187
 32.53.720.550
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



553900 Rev. 7 Page 1 of 2

- 5. Add 200 μl of staining buffer/well and resuspend cells. Centrifuge at 250×g for 5 minutes and aspirate supernatant.
- 6. Block surface Ig with purified RG7/11.1 mAb (Cat. No. 553897) by adding 1.0 µg per sample in 50 µl of staining buffer/well.

Note: Surface markers may be stained during this step as described in the "Immunofluorescent Staining of Mouse and Rat Leukocytes for Flow Cytometry" in the Technical Protocols section of our web site at

http://www.bdbiosciences.com/pharmingen/protocols/Mouse_and_Rat_Leukocytes.shtml

- 7. Incubate 15 minutes on ice.
- 8. Wash 2× as described in Step 5.
- 9. Resuspend cells in 100 μl of BD Cytofix/CytopermTM intracellular staining buffer (BD Cytofix/CytopermTM Kit, Cat. No. 554714) per well.
- 10. Incubate 30 minutes at room temperature.
- 11. Wash $2\times$ with $200~\mu$ l of $1\times$ Perm/Wash buffer (provided in the BD Cytofix/Cytoperm Kit) per well. Centrifuge at $250\times$ g for 5 minutes and aspirate supernatant between washes.
- 12. Stain intracellular Ig by adding ≤ 1 μg of FITC-conjugated RG7/11.1 mAb in 50 μl of 1 × Perm/Wash buffer/well.

Note: Other antibodies recommended for staining of intracellular markers may be added during this step as described in Step 12.

- 13. Incubate for 30 minutes at room temperature.
- 14. Wash 2× as described in Step 11.
- 15. Resuspend and transfer samples in 100 μl of staining buffer to tubes appropriate for analysis with a flow cytometer. Bring volume in each tube to 400 μl with staining buffer.
- 16. Analyze samples on a flow cytometer.

Suggested Companion Products

Catalog Number	Name	Size	Clone
555057	FITC Mouse IgG2b, κ Isotype Control	0.1 mg	27-35
553141	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block™)	0.1 mg	2.4G2
553897	Purified Mouse Anti-Rat IgG2b	0.5 mg	RG7/11.1
554656	Stain Buffer (FBS)	500 ml	(none)
554714	BD Cytofix/Cytoperm Fixation/Permeablization Kit	250 tests	(none)

Product Notices

- 1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

References

Springer TA, Bhattacharya A, Cardoza JT, Sanchez-Madrid F. Monoclonal antibodies specific for rat IgG1, IgG2a, and IgG2b subclasses, and kappa chain monotypic and allotypic determinants: reagents for use with rat monoclonal antibodies. *Hybridoma*. 1982; 1(3):257-273.(Immunogen)

553900 Rev. 7 Page 2 of 2