Technical Data Sheet FITC Mouse IgG3, κ Isotype Control

Product Information


Material Number:
Alternate Name:
Size:
Vol. per Test:
Clone:
Isotype:
Reactivity:
Storage Buffer:

555578 anti-Fructosan 100 tests 20 μl J606 Mouse (BALB/c) IgG3, κ QC Testing: Human Aqueous buffered solution containing BSA, protein stabilizer, and ≤0.09% sodium azide.

Description

Mouse $IgG3 \kappa$ isotype control antibody is specific for the β -2 - 6-linked Fructosan. The transplantable plasmacytoma J606 was induced by intraperitoneal injection of mineral oils into BALB/c mice.

This antibody is routinely tested by flow cytometric analysis. Other applications were tested at BD Biosciences Pharmingen during antibody development only or reported in the literature.

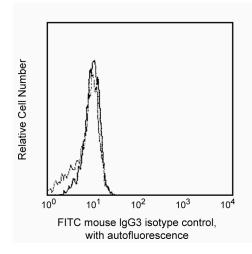

Profile of peripheral blood lymphocytes analyzed on a FACScan (BDIS, San Jose, CA)

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed. Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

A	Application	
	Flow cytometry	Routinely Tested
	Isotype control	Routinely Tested

Recommended Assay Procedure:

An isotype control should be used at the same concentration as the antibody of interest.

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^{6} cells in a 100-µl experimental sample (a test).

BD Biosciences

bdbiosciences.	com				
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
77.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995
or country-sp	ecific contact in	formation, visit	bdbiosciences.co	m/how_to_orde	r/
of any patents. Bl use of our produc product or as a co written authoriza For Research Use	D Biosciences will n ets. Purchase does n omponent of anoth tion of Becton Dick Only. Not for use ir	ot be held responsi not include or carry er product. Any us kinson and Compan n diagnostic or ther	ble for patent infrin any right to resell oi	gement or other vio r transfer this produc er than the permitte d. Not for resale.	e the above product in violation lations that may occur with the t either as a stand-alone ed use without the express 08 BD



- 2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.