

Technical Data Sheet


Biotin Mouse Anti-Human CD28

Product Information

| | |
|-------------------------|---|
| Material Number: | 555727 |
| Alternate Name: | CD28 antigen; T44; Tp44; TP44 |
| Size: | 100 Tests |
| Vol. per Test: | 20 µl |
| Clone: | CD28.2 |
| Immunogen: | Human CD28 Transfected Cell Line |
| Isotype: | Mouse (C3H x BALB/c) IgG1, κ |
| Reactivity: | QC Testing: Human |
| Workshop: | V 5T CD28.05 |
| Storage Buffer: | Aqueous buffered solution containing BSA and ≤0.09% sodium azide. |

Description

The CD28.2 monoclonal antibody specifically binds to CD28, a 44 kDa homodimeric transmembrane glycoprotein present on most mature T cells, thymocytes and plasma cells. CD28 is a costimulatory receptor that binds CD80 and CD86 as ligands and plays a very important role in T cell-B cell interactions. It has been suggested that CD28 initiates and regulates a separate and distinct signal transduction pathway from those stimulated by the TCR complex. Additionally, it has been reported that CD28 antibody clones vary in their ability to stimulate T cells to produce IL-2 and increase intracellular Ca²⁺ concentration. This finding suggests the existence of functionally distinct subregions on the CD28 molecule. CD28.2 has been demonstrated to bind to the same molecule as clone L293, another CD28 mAb, and has been reported to induce Ca²⁺ influx in Jurkat T cells.

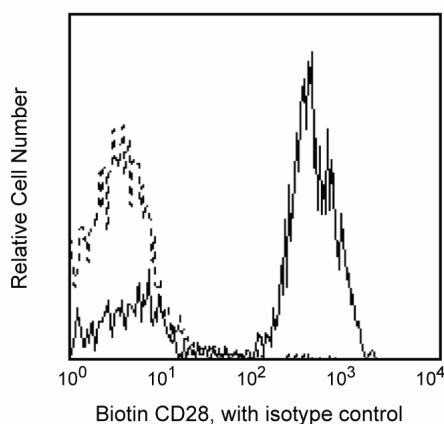

Profile of peripheral blood lymphocytes analyzed on a FACSscan (BDIS, San Jose, CA)

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.
The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.
The antibody was conjugated with biotin under optimum conditions, and unreacted biotin was removed.

Application Notes

Application

| | |
|----------------|------------------|
| Flow cytometry | Routinely Tested |
|----------------|------------------|

Suggested Companion Products

| Catalog Number | Name | Size | Clone |
|----------------|-------------------------------------|-----------|---------|
| 555747 | Biotin Mouse IgG1 κ Isotype Control | 100 Tests | MOPC-21 |
| 554061 | PE Streptavidin | 0.5 mg | (none) |
| 555899 | Lysing Buffer | 100 mL | (none) |
| 554656 | Stain Buffer (FBS) | 500 mL | (none) |

BD Biosciences

bdbiosciences.com

| | | | | | |
|---------------|--------------|----------------|--------------|--------------|-------------------------|
| United States | Canada | Europe | Japan | Asia Pacific | Latin America/Caribbean |
| 877.232.8995 | 866.979.9408 | 32.2.400.98.95 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995 |

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

555727 Rev. 8



Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Please refer to www.bdbiosciences.com/pharming/en/protocols for technical protocols.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
6. An isotype control should be used at the same concentration as the antibody of interest.
7. This product is sold under license. Purchase of this product does not include rights to (i) incorporate this product into the purchaser's own products for resale to end-users, or (ii) use this product to conduct for-profit research for or on behalf of another party. For information on obtaining a license to this product for such prohibited uses, contact INSERM, 7 rue Watt, 75013 Paris. Telephone: +33 1 55 03 01 60. Facsimile: +33 1 55 03 01 18. Email: techtransfert@inserm-transfert.fr

References

Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997. (Biology)

Guesdon JL, Ternynck T, Avrameas S. The use of avidin-biotin interaction in immunoenzymatic techniques. *J Histochem Cytochem*. 1979; 27(8):1131-1139. (Biology)

June CH, Bluestone JA, Nadler LM, Thompson CB. The B7 and CD28 receptor families. *Immunol Today*. 1994; 15(7):321-331. (Biology)

Kuiper H, Brouwer M, Vermeire S, van Lier R. Analysis of the Workshop CD28 Panel mAb: distinct signalling pathways coupled to CD28. In: Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. Oxford: Oxford University Press; 1995:373-374. (Clone-specific: Activation, Calcium Flux, (Co)-stimulation)

Nunes J, Klasen S, Franco MD, et al. Signalling through CD28 T-cell activation pathway involves an inositol phospholipid-specific phospholipase C activity. *Biochem J*. 1993; 293(3):835-842. (Biology)

Nunes J, Klasen S, Ragueneau M, et al. CD28 mAbs with distinct binding properties differ in their ability to induce T cell activation: analysis of early and late activation events. *Int Immunol*. 1993; 5(3):311-315. (Biology)

Olive D, Cerdan C, Costello R, Sielleur I, Ragueneau M, Pages F, Klasen S, Nunes J, Imbert J. CD28 and CTLA-4 cluster report. In: Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. Oxford: Oxford University Press; 1995:360-370. (Clone-specific: (Co)-stimulation, Flow cytometry, Functional assay, Inhibition, Stimulation)

Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. Oxford: Oxford University Press; 1995. (Clone-specific)

Verwilghen J, Vandenberghe P, Wallays G, et al. Simultaneous ligation of CD5 and CD28 on resting T lymphocytes induces T cell activation in the absence of T cell receptor/CD3 occupancy. *J Immunol*. 1993; 150(3):835-846. (Biology)

BD Biosciences

bdbiosciences.com

| United States | Canada | Europe | Japan | Asia Pacific | Latin America/Caribbean |
|---------------|--------------|----------------|--------------|--------------|-------------------------|
| 877.232.8995 | 866.979.9408 | 32.2.400.98.95 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995 |

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

555727 Rev. 8

