Technical Data Sheet

PE Mouse Anti-Human CD70

Product Information


555835 **Material Number:** 100 tests Size: 20 µl Vol. per Test: Ki-24 Clone: Mouse IgG3, κ Isotype: QC Testing: Human Reactivity:

IV A109 Workshop:

Aqueous buffered solution containing BSA and ≤0.09% sodium azide. Storage Buffer:

Description

Reacts with an antigen which immunoprecipitates as five bands (50, 70, 90, 100 and 160 kDa) under non-reducing conditions. CD70 is strongly expressed on Reed-Sternberg cells, some activated T or B cells and Epstein Barr Virus (EBV)+ lymphoblastoid cell lines. Reports suggest that CD70 plays a role in the regulation of special T-cell proliferative functions.

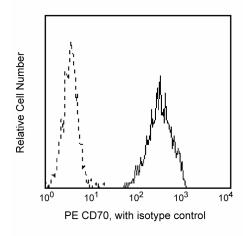

Profile of peripheral U266 cells (human plasma cell line) analyzed by flow cytometry

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed. Store undiluted at 4° C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

 FF	-	
Flow cytometry	Routinely Tested	

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 X 10e6 cells in a 100-µl experimental
- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/pharmingen/colors

BD Biosciences

www.bdbiosciences.com

Europe Asia Pacific 32.53.720.550 0120.8555.90 877.232.8995 888.259.0187 65.6861.0633 55.11.5185.9995 For country-specific contact information, visit www.bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation

Conditions: Ine information disclosed nerien is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2007 BD



- 5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

References

Schlossman SF, Boumsell L, Gilks W, et al, ed. Leukocyte Typing V: White Cell Differentiation Antigens. New York: Oxford University Press; 1995.(Biology) Knapp W, Dorken B, Rieber EP, et al, ed. Leucocyte Typing IV. New York: Oxford University Press; 1989.(Clone-specific)

555835 Rev. 8 Page 2 of 2