


Technical Data Sheet

Nile Red (High Intensity) Fluorescent Particles, 1.7 - 2.2 µm**Product Information**

Material Number: 556261
Size: 2.0 ml
Storage Buffer: Aqueous buffered solution containing ≤0.09% sodium azide.

Description

The vial contains a mixture of 1.7 - 2.2 µm Nile Red Fluorescent Particles that are dyed to have an excitation maximum of approximately 525 nm with a subsequent emission maximum of approximately 560 nm. This particle mixture is useful with flow cytometers equipped with a 488 nm-argon laser. This particle mixture has compatible excitation and emission spectra with R-Phycoerythrin and Propidium iodide.

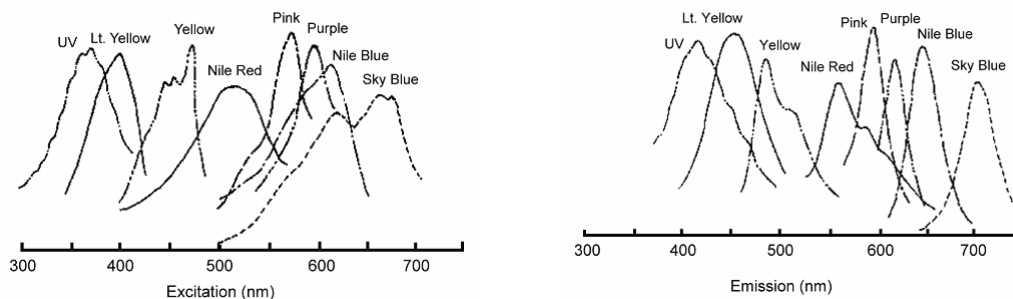

Excitation and emission spectra of SPHERO fluorescent particles in suspension with the exception of Nile Blue and Sky Blue, which were read from a solution of dye in styrene.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes**Recommended Assay Procedure:**

This particle mixture (~1.0% w/v) is useful for routine calibration of flow cytometers. Before use, resuspend the particles by vortexing. Dilution of 3-5 drops of particles to 1 ml of sheath fluid will provide an adequate number of particles for flow cytometric analysis.

Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD

