# Technical Data Sheet

# PE-Cy<sup>™</sup>5 Mouse Anti-Human CD18

#### **Product Information**

**Material Number:** 557528 Alternate Name: Integrin β2 chain 100 tests Size Vol. per Test: 20 µl Clone: 6.7

**Isotype:** Mouse IgG1, κ Reactivity: QC Testing: Human


Tested in Development: Rhesus, Cynomolgus, Pig, Rabbit

Workshop:

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

#### Description

Reacts with β2 integrin, a 90-95 kDa transmembrane protein, which associates non-covalently as a heterodimer with the a chains of CD11a, CD11b, or CD11c. CD18 is expressed on lymphocytes, monocytes, and more weakly on granulocytes. LFA-1 has been shown to play an important role in homotypic and heterotypic cellular adhesion in immune and inflammatory responses.

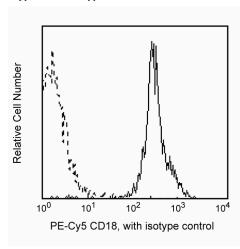

Profile of peripheral blood lymphocytes analyzed by flow cvtometry

## **Preparation and Storage**

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with PE-Cy5 (formerly known as BD Cy-Chrome<sup>TM</sup>) under optimum conditions, and unconjugated antibody and free PE-Cv5 were removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

#### **Application Notes**

Application

Flow cytometry Routinely Tested

# **Suggested Companion Products**

Catalog Number Clone Size 555750 PE-Cy<sup>TM</sup>5 Mouse IgG1 κ Isotype Control 100 tests MOPC-21

## **Product Notices**

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^{6}$  cells in a  $100 \mu l$  experimental sample (a test).
- 2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

## **BD Biosciences**

bdbiosciences.com

Asia Pacific Europe 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



557528 Rev. 9

- 5. PE-Cy5 is optimized for use with a single argon ion laser emitting 488-nm light. Because of the broad absorption spectrum of the PE-Cy5 tandem fluorochrome, extra care must be taken when using dual-laser cytometers which may directly excite both PE and Cy5<sup>TM</sup>.
- Cy is a trademark of Amersham Biosciences Limited. This conjugated product is sold under license to the following patents: US Patent Nos. 5,486,616; 5,569,587; 5,569,766; 5,627,027.
- 7. This product is subject to proprietary rights of Amersham Biosciences Corp. and Carnegie Mellon University and made and sold under license from Amersham Biosciences Corp. This product is licensed for sale only for research. It is not licensed for any other use. If you require a commercial license to use this product and do not have one return this material, unopened to BD Biosciences, 10975 Torreyana Rd, San Diego, CA 92121 and any money paid for the material will be refunded.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 9. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

#### References

Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997. (Biology) Schlossman S, Boumell L, et al, ed. *Leucocyte Typing V*. New York: Oxford University Press; 1995. (Clone-specific)

557528 Rev. 9 Page 2 of 2