Technical Data Sheet

Alexa Fluor®488 Mouse Anti-Human CD4

Product Information


557695 **Material Number:** 100 tests Size: Vol. per Test: 5 μl Clone: RPA-T4 Mouse IgG1, κ Isotype: QC Testing: Human Reactivity:

IV T114 Workshop:

Aqueous buffered solution containing BSA and ≤0.09% sodium azide. Storage Buffer:

Description

The RPA-T4 clone reacts with CD4, a 59 kDa single-chain transmembrane glycoprotein [receptor for human immunodeficiency virus (HIV)] present on T-helper/inducer cell populations. This antibody binds to the D1 domain (CDR1 and CDR3 epitopes) of the CD4 antigen and reacts with approximately 80% of thymocytes and 45% of peripheral blood lymphocytes. CD4 is also present in low density on peripheral blood monocytes. RPA-T4 is capable of blocking HIV-1, gp120, and inhibits syncytium formation.

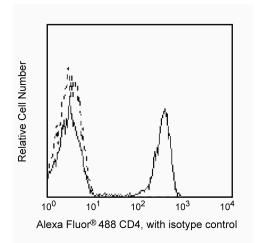

Profile of peripheral blood lymphocytes analyzed by flow cvtometry.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 488 under optimum conditions, and unreacted Alexa Fluor® 488 was removed.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone
557702	Alexa Fluor® 488 Mouse IgG1 κ Isotype Control	100 tests	MOPC-21
555899	Lysing Buffer	100 ml	(none)

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use $1 \times 10^{\circ}6$ cells in a 100- μ l experimental sample (a test).
- An isotype control should be used at the same concentration as the antibody of interest.
- Alexa Fluor® 488 fluorochrome emission is collected at the same instrument settings as for fluorescein isothiocyanate (FITC). 3.
- Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- The Alexa Fluor®, Pacific BlueTM, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.

BD Biosciences

bdbiosciences.com **United States** Canada Asia Pacific Latin America/Caribbean 877.232.8995 888.259.0187 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995 For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



- 6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 7. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 8. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- $9. \quad Please \ refer \ to \ www.bdbiosciences.com/pharmingen/protocols \ for \ technical \ protocols.$

References

Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989. (Biology)
Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V*: *White Cell Differentiation Antigens*. New York: Oxford University Press; 1995. (Biology)

557695 Rev. 4 Page 2 of 2