

Technical Data Sheet

Alexa Fluor® 647 Mouse anti-Btk (pY551)/Itk (pY511)**Product Information**


Material Number:	558134
Size:	50 tests
Vol. per Test:	20 µl
Clone:	24a/BTK (Y551)
Immunogen:	Phosphorylated Human Btk Peptide
Isotype:	Mouse IgG1, κ
Reactivity:	QC Tested: Human
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Bruton's tyrosine kinase (Btk) is a nonreceptor tyrosine kinase whose function is critical for proper B cell development and signaling. The activity of Btk is regulated by Src mediated phosphorylation of the kinase domain at tyrosine 551 (Y551). This event induces Btk kinase activity and subsequent autophosphorylation at Y223. Phosphorylated Btk then associates with the cell membrane via the interaction of the PH domain with phosphatidylinositol 3, 4, 5-triphosphate.

The Tec family kinase Itk plays a critical role in signal transduction downstream of the T cell antigen receptor and has been implicated in the activation of phospholipase C-γ1 (PLCγ1), a key regulator of calcium mobilization and extracellular signal-regulated kinase (ERK) activation. Itk is regulated by an activating transphosphorylation event in which Y511 in the kinase domain is phosphorylated by Lck.

The 24a/BTK (Y551) monoclonal antibody recognizes the Y551-phosphorylated form of human Btk and the Y511 phosphorylated form of human Itk.

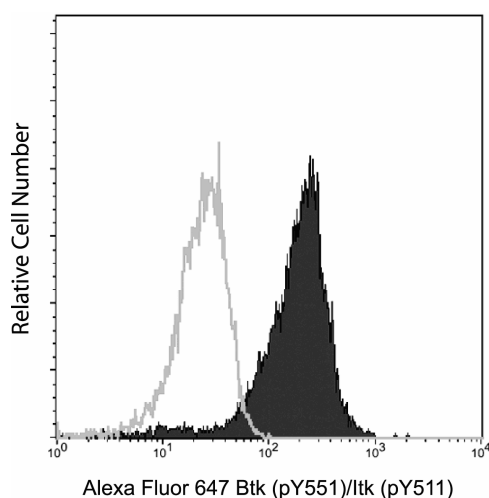

Analysis of Btk (pY551) in activated human B lymphoma cells. Ramos cells (ATCC CRL-1596) were serum starved overnight and then either stimulated with 5 mM hydrogen peroxide for 15 minutes (shaded histogram) or unstimulated (open histogram). The cells were fixed (BD Cytofix™ Fixation buffer, Cat. No. 554655) for 10 minutes, then permeabilized (BD Phosflow™ Perm Buffer III, Cat. No. 558050) on ice for at least 30 minutes, and then stained with Alexa Fluor® 647 anti-Btk (pY551)/Itk (pY511). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes**Application**

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Suggested Companion Products

Catalog Number	Name	Size	Clone
554655	Fixation Buffer	100 ml	(none)
558050	Perm Buffer III	125 ml	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.
3. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
4. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
8. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.

References

Mahajan S, Fargnoli J, Burkhardt AL, Kut SA, Saouaf SJ, Bolen JB. Src family protein tyrosine kinases induce autoactivation of Bruton's tyrosine kinase. *Mol Cell Biol.* 1995; 15:5304-5311. (Biology)

Rawlings DJ, Scharenberg AM, Park H, et al. Activation of BTK by a phosphorylation mechanism initiated by SRC family kinases. *Science.* 1996; 271:822-825. (Biology)

Wilcox HM, Berg LJ.. Itk phosphorylation sites are required for functional activity in primary T cells. *J Biol Chem.* 2003; 278:37112-37121. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

