BD™ Cytometric Bead Array (CBA)



Technical Data Sheet Human IL-10 Flex Set

Product Information

- Material Number: Size: Bead Position: Assay Range: Reactivity:
- **558274** 100 Tests B7 10-2,500 pg/mL QC Testing: Human

Component Description: Component Mat. No: Component Storage Buffer:

Component Description: Component Mat. No: Component Storage Buffer:

Component Description: Component Mat. No: Component Storage Buffer: Human IL-10 Standard 51-9003501 Lyophilized in an aqueous buffered solution containing BSA and ProClin[™] 150. Human IL-10 PE Detection Reagent

51-9004043 Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.

Human IL-10 Capture Bead B7 51-9005291 Aqueous buffered solution containing fetal bovine serum and ≤0.09% sodium azide.

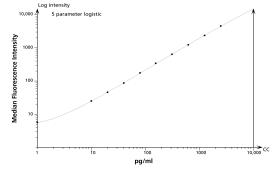


Figure 1. Example BD CBA Human IL-10 Flex Set standard curve. Data acquired on a BD FACSArray bioanalyzer and analyzed using FCAP Array Software.

Description

The BDTM CBA Human IL-10 Flex Set is a bead-based immunoassay capable of measuring human interleukin-10 (IL-10) in serum, plasma, and cell culture supernatant samples. Human reactivity was determined by testing samples with the BD CBA Human IL-10 Flex Set. The biology and function of IL-10 has been extensively reviewed in the literature. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Human Soluble Protein Master Buffer Kit (Cat. No. 558264 or 558265).

Preparation and Storage

This BDTM CBA Flex Set contains one vial each of Capture Bead and PE Detection Reagent and two vials of Standard. The Capture Bead and PE Detection Reagent components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. The Standard component is lyophilized and should be transferred to a 15 mL polypropylene tube for reconstitution. When reconstituted in 4.0 mL Assay Diluent, the standard has a protein concentration of 2,500 pg/mL. Discard unused reconstituted standard, do not store or reuse. Store lyophilized standard and other components at 4°C. Protect Capture Beads and the PE Detection Reagent from prolonged exposure to light.

Application Notes

Recommended Assay Procedure: The BD CBA Human IL-10 Flex Set must be used in conjunction with a BD CBA Human Soluble Protein Master Buffer Kit (Cat. No. 558264, 100 tests, or 558265, 500 tests), a flow cytometer, and FCAP ArrayTM Software. Detailed instructions on the use of this product can be found in the manual for the BD CBA Human Soluble Protein Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Human IL-10 Flex Set will be 2,500 pg/mL. An example standard curve is shown in Figure 1.

The BD CBA Human IL-10 Flex Set should not be used in the same assay well with any non-BD CBA Human Soluble Protein Flex Set reagents

BD Biosciences

bdbiosciences.com								
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbeau			
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995			
For country cor	ntact informatio	on, visit bdbioscie	ences.com/conta	ct				
of any patents. BL	Delissciences will n	ot be held responsik	ole for patent infrin	gement or other vio	e the above product in violation			
use of our produc	ts. Purchase does n	ot include or carry a	any right to resell or	r transfer this produ	vlations that may occur with the			
product or as a co	mponent of anoth	er product. Any use	of this product oth	er than the permitte	ct either as a stand-alone			
written authoriza	tion of Becton, Dic	kinson and Compan	y is stictly prohibite	d.	ed use without the express			
For Research Use	Only. Not for use ir	diagnostic or thera	peutic procedures.	Not for resale.	an and Company. © 2014 BD			

(such as BD CBA Mouse Soluble Protein or Cell Signaling Flex Sets). For an updated assay compatibility chart for the BD CBA Human Soluble Protein Flex Sets, please refer to the BD CBA Flex Set System homepage at http://www.bdbiosciences.com/cbasetup.

Performance

Limit of Detection: The theoretical limit of detection is 0.13 pg/mL and was determined by evaluating the estimated result of the average MFI of the negative control (0 pg/mL, n=30) + 2 standard deviations.

Specificity		Inter-Assay Reproducibility			Intra-Assay Reproducibility			
		Mean (pg/ml)	Standard Deviation	%CV	Mean (pg/ml)	Standard Deviation	%CV	
Sample 1		39.4	2.2	6%	36.4	2.1	6%	
Human IL-10	Sample 2	146.6	14.5	10%	125.3	2.6	2%	
	Sample 3	585.6	65.8	11%	491.5	17.1	3%	

Reproducibility: The inter-assay and intra-assay reproducibility were determined for the BD CBA Human IL-10 Flex Set by evaluating ten replicates of three different sample levels (intra-assay) and two replicates of three different sample levels from four separate experiments (intra-assay) respectively.

	Cell Culture Supernatant		Serum		Plasma	
Specificity	Average % Recovery	Range	Average % Recovery	Range	Average % Recovery	Range
Human IL-10	83%	77 - 90%	80%	76 - 84%	74%	72 - 77%

Recovery: Cell culture supernatant, serum, or EDTA-treated plasma were spiked with three different levels of protein. The spiked samples were assayed and the results were compared with expected values. Serum and plasma samples were diluted 1:4 before the protein was spiked into each. Serum is a pool of 800 - 1000 donors and the plasma was pooled from at least 20 donors.

Comple Dilution	Cell Culture Supernatant		Ser	um	Plasma		
Sample Dilution	Detected (pg/ml)	% of Expected	Detected (pg/ml)	% of Expected	Detected (pg/ml)	% of Expected	
Spiked sample	525.1	100%	517.6	100%	524.1	100%	
1:2	251.3	96%	250.1	97%	259.6	99%	
1:4	124.9	95%	125.4	97%	133.9	102%	

Linearity: Cell culture supernatant, 1:4 diluted serum, or 1:4 diluted EDTA-treated plasma were spiked with protein and serially diluted. The diluted samples were assayed and the results were compared with the original spiked sample.

Product Notices

- 1. ProClin is a trademark of Rohm and Haas Company.
- 2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- Warning: CBA lyophilized standard contains 0.02% (w/w) of a CMIT/MIT mixture (3:1), which is a mixture of: 5-chloro-2-methyl-4-isothiazolin-3-one [EC No 247-500-7] and 2-methyl-4-isothiazolin-3-one [EC No 220-239-6] (3:1). Hazard statement: May cause an allergic skin reaction.

Precautionary statements: Wear protective gloves/eye protection. Wear protective clothing. Avoid breathing mist/vapours/spray. If skin irritation or rash occurs: Get medical advice/attention. IF ON SKIN: Wash with plenty of water. Dispose of contents/container in accordance with local/regional/national/international regulations.

bdbiosciences.com								
United States 877.232.8995	Canada 800.268.5430	Europe 32.2.400.98.95	Japan 0120.8555.90	Asia Pacific 65.6861.0633	Latin America/Caribbear 55.11.5185.9995			
For country co	ntact informatio	on, visit bdbiosci e	ences.com/conta	ct				
of any patents. BL use of our produc product or as a co	D Biosciences will n ts. Purchase does n imponent of anoth	ot be held responsil ot include or carry a	ble for patent infrin any right to resell or of this product oth	gement or other vio [•] transfer this produ er than the permitte	e the above product in violatior lations that may occur with the ct either as a stand-alone ed use without the express			

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale. Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

