BD™ Cytometric Bead Array (CBA)



Technical Data Sheet Mouse IL-13 Flex Set

Product Information

- Material Number: Size: Bead Position: Assay Range: Reactivity:
- **558349** 100 Tests B8 10-2,500 pg/mL QC Testing: Mouse


Component Description: Component Mat. No: Component Storage Buffer:

Component Description: Component Mat. No: Component Storage Buffer: 51-9004266 Lyophilized in an aqueous buffered solution containing BSA and ProClin[™] 150. Mouse IL-13 PE Detection Reagent 51-9004346

Mouse IL-13 Standard

Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.

Component Description: Component Mat. No: Component Storage Buffer: Mouse IL-13 Capture Bead B8 51-9005260 Aqueous buffered solution containing fetal bovine serum and ≤0.09% sodium azide.

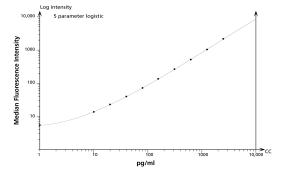

Example BD CBA Mouse IL-13 Flex Set standard curve. Data acquired on a BD FACSArray bioanalyzer and analyzed using FCAP Array Software.

Description

The BDTM CBA Mouse IL-13 Flex Set is a bead-based immunoassay capable of measuring mouse interleukin-13 (IL-13) in serum and cell culture supernatant samples. Mouse reactivity was determined by testing samples with the BD CBA Mouse IL-13 Flex Set. The biology and function of IL13 has been extensively reviewed in the literature. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266 or 558267).

Preparation and Storage

This BDTM CBA Flex Set contains one vial each of Capture Bead and PE Detection Reagent and two vials of Standard. The Capture Bead and PE Detection Reagent components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. The Standard component is lyophilized and should be transferred to a 15 mL polypropylene tube for reconstitution. When reconstituted in 4.0 mL Assay Diluent, the standard has a protein concentration of 2,500 pg/mL. Discard unused reconstituted standard, do not store or reuse. Store lyophilized standard and other components at 4°C. Protect Capture Beads and the PE Detection Reagent from prolonged exposure to light.

Application Notes

Recommended Assay Procedure:

The BD CBA Mouse IL-13 Flex Set must be used in conjunction with a BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266, 100 tests, or 558267, 500 tests), a flow cytometer, and FCAP Array™ Software. Detailed instructions on the use of this product can be found in the manual for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Mouse IL-13 Flex Set will be 2,500 pg/mL. An example standard curve is shown in figure 1.

The BD CBA Mouse IL-13 Flex Set should not be used in the same assay well with any non-BD CBA Mouse Soluble Protein Flex Set reagents

BD Biosciences

bdbiosciences.com								
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbear			
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995			
For country cor	ntact informatio	on, visit bdbiosci e	ences.com/conta	ct				
of any patents. BL	Delissciences will not	ot be held responsik	ble for patent infrin	gement or other vio	e the above product in violation			
use of our produc	ts. Purchase does n	ot include or carry a	any right to resell or	r transfer this produ	vlations that may occur with the			
product or as a co	mponent of anoth-	er product. Any use	of this product oth	er than the permitte	ct either as a stand-alone			
written authorizat	tion of Becton, Dick	kinson and Compan	by is stictly prohibite	d.	ed use without the express			
For Research Use (Only. Not for use in	diagnostic or thera	apeutic procedures.	Not for resale.	an and Company. © 2014 BD			



(such as BD CBA Human or Rat Soluble Protein or Cell Signaling Flex Sets). For an updated assay compatibility chart for the BD CBA Mouse Soluble Protein Flex Sets, please refer to the BD CBA Flex Set System homepage at http://www.bdbiosciences.com/cbasetup.

Performance

Limit of Detection: The theoretical limit of detection is 2.4 pg/mL and was determined by evaluating the estimated result of the average MFI of the negative control (0 pg/mL, n=30) + 2 standard deviations.

Specificity		Inter-Assay Reproducibility			Intra-Assay Reproducibility		
		Mean (pg/ml)	Standard Deviation	%CV	Mean (pg/ml)	Standard Deviation	%CV
	Sample 1	38.2	1.8	5%	36.6	2.3	6%
Mouse IL-13	Sample 2	149.9	4.8	3%	145.5	5.1	3%
	Sample 3	629.9	12.8	2%	621.9	22.8	4%

Reproducibility: The inter-assay and intra-assay reproducibility were determined for the BD CBA Mouse IL-13 Flex Set by evaluating ten replicates of three different sample levels (inter-assay) and two replicates of three different sample levels from four separate experiments (intra-assay) respectively.

		ulture natant	Serum		
Sample Dilution	Average % Recovery	Range	Average % Recovery	Range	
Mouse IL-13	87%	86 - 90%	29%	26 - 31%	

Recovery: Cell culture supernatant and pooled mouse serum were spiked with three different levels of protein. The spiked samples were assayed and the results were compared with expected values. Serum samples were diluted 1:4 before the protein was spiked into each. Serum used was commercially available pooled mouse serum.

		culture natant	Serum		
Sample Dilution	Detected (pg/ml)	% of Expected	Detected (pg/ml)	% of Expected	
Spiked sample	549.2	100%	306.6	100%	
1:2	267.7	97%	191.4	125%	
1:4	131.6	96%	106.5	139%	

Linearity: Cell culture supernatant and 1:4 diluted pooled mouse serum were spiked with protein and serially diluted. The diluted samples were assayed and the results were compared with the original spiked sample.

Product Notices

- 1. ProClin is a trademark of Rohm and Haas Company.
- 2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- Warning: CBA lyophilized standard contains 0.02% (w/w) of a CMIT/MIT mixture (3:1), which is a mixture of: 5-chloro-2-methyl-4-isothiazolin-3-one [EC No 247-500-7] and 2-methyl-4-isothiazolin-3-one [EC No 220-239-6] (3:1). Hazard statement: May cause an allergic skin reaction.

Precautionary statements: Wear protective gloves/eye protection. Wear protective clothing. Avoid breathing mist/vapours/spray. If skin irritation or rash occurs: Get medical advice/attention. IF ON SKIN: Wash with plenty of water. Dispose of contents/container in accordance with local/regional/national/international regulations.

bdbiosciences.com							
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbear		
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995		
For country co	ntact informatio	on, visit bdbiosci e	ences.com/conta	ct			
of any patents. BL use of our produc	D Biosciences will n ts. Purchase does n	ot be held responsil ot include or carry a	ble for patent infrin any right to resell of	gement or other vio • transfer this produ	e the above product in violatior lations that may occur with the ct either as a stand-alone of use without the express		

written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

Page 2 of 2

