Technical Data Sheet

Alexa Fluor® 647 Mouse anti-SLP-76 (pY128)

Product Information

 $\begin{tabular}{llll} \textbf{Material Number:} & 558438 \\ \textbf{Size:} & 50 tests \\ \textbf{Vol. per Test:} & 20 \ \mu l \\ \end{tabular}$

Clone: J141-668.36.58

 Immunogen:
 Phosphorylated Human SLP-76

 Isotype:
 Mouse (BALB/c) IgG1, κ

 Reactivity:
 Confirmed: Human


 Predicted: Mouse, Rat


Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

SLP-76 (SH2 domain-containing Leukocyte Protein of 76 kDa) is a tyrosine phosphoprotein that is involved in the T cell receptor (TCR) -mediated intracellular signaling pathway. It may be involved in the signaling pathways of other peripheral blood leukocytes; thymic/splenic cells; and in human T, B, and monocytic cell lines. SLP-76 consists of several motifs that signify its importance in protein-protein interactions involved in intracellular signaling pathways, such as the SH2 domain in the C-terminus, the three amino-terminus 17-amino acid repeats with conserved tyrosine and acidic residues (DYE(S/P)P), and a proline rich region. SLP-76 has been shown to associate with Gads, Grb2, PLCγ1, SLAP-130, and Vav, all of which are part of the signaling cascade in T lymphocytes. An early event in the T cell activation pathway is the phosphorylation, by the Syk-family kinase ZAP-70, of SLP-76 at the three conserved tyrosine motifs, which then mediate interactions with downstream effectors. The phosphorylated tyrosine 128 (Y128) brings the Rho family guanine nucleotide exchange factor Vav1 and the Nck adapter protein, which binds to p21-activated kinase (PAK1) and Wiskott-Aldrich syndrome protein (WASP), into the activation complex. Vav1, PAK1, and WASP may mediate TCR-stimulated actin cytoskeletal rearrangement.

The J141-688.36.58 monoclonal antibody recognizes the phosphorylated Y128 of activated SLP-76.

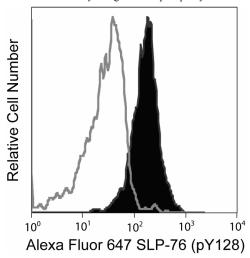

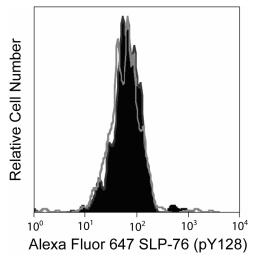

Analysis of SLP-76 (pY128) in human peripheral blood lymphocytes. Human peripheral blood mononuclear cells were either unstimulated (open histograms) or stimulated by cross-linking with CD3 (Purified NA/LE Mouse anti-Human CD3 mAb, clone UCHT/1, Cat. No. 555329) and CD28 (Purified NA/LE Mouse anti-Human CD28, clone CD28.2, Cat. No. 555725) on ice for 15 minutes followed by incubation with the Purified Goat anti-Mouse Ig (Cat. No. 553998) on ice for 15 minutes, and then allowed to undergo phosphorylation at 37°C for 5 minutes (shaded histograms). The cells were fixed (BD Cytofix ™ Fixation buffer Cat. No. 554655) for 10 minutes at 37°C, permeabilized (BD Phosflow™ Perm Buffer III, Cat. No. 558050) on ice for 30 minutes, blocked with normal mouse immunoglobulin, and then stained with Alexa Fluor® 647 Mouse anti-SLP-76 (pY128) (Cat. No. 558438) and PerCP-Cy5.5 Mouse anti-human CD3 mAb SK7 (Cat. No. 340949). The left panel displays the upregulated phosphorylation of SLP-76 in activated CD3-positive T lymphocytes, while the right panel displays the lack of upregulated SLP-76 phosphorylation in CD3-negative non-T cells. Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

BD Biosciences

bdbiosciences.com

United StatesCanadaEuropeJapanAsia PacificLatin America/Caribbean877.232.8995888.268.543032.53.720.5500120.8555.9065.6861.06330800.771.7157For country-specific contact information, visitbdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested	
---	------------------	--

Recommended Assay Procedure:

This antibody conjugate is suitable for intracellular staining of human whole blood (using BD PhosflowTM Lyse/Fix Buffer, Cat. No. 558049) and peripheral blood mononuclear cells (using BD CytofixTM Fixation Buffer, Cat. No. 554655). Any of the 3 BD PhosflowTM Perm Buffers may be used.

Suggested Companion Products

Catalog Number	Name	Size	Clone	
558049	Lyse/Fix Buffer 5X	250 ml	(none)	
558050	Perm Buffer III	125 ml	(none)	
554655	Fixation Buffer	100 ml	(none)	
555329	Purified NA/LE Mouse Anti-Human CD3	0.5 mg	UCHT1	
554656	Stain Buffer (FBS)	500 ml	(none)	
555725	Purified NA/LE Mouse Anti-Human CD28	0.5 mg	CD28.2	
558052	Perm Buffer II	125 ml	(none)	
557885	Perm/Wash Buffer I	125 ml	(none)	

Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 3. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 5. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
- 6. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
- 7. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- 8. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Fang N, Motto DG, Ross SE, Koretzky GA. Tyrosines 113, 128, and 145 of SLP-76 are required for optimal augmentation of NFAT promoter activity. *J Immunol.* 1996; 157:3769-3773. (Biology)

Janssen E, Zhang W. Adaptor proteins in lymphocyte activation. *Curr Opin Immunol*. 2003; 15:269-276. (Biology) Wu JN, Koretzky GA. The SLP-76 family of adapter proteins. *Semin Immunol*. 2004; 16:379-393. (Biology)

558438 Rev. 2 Page 2 of 2