

Technical Data Sheet

Alexa Fluor® 647 Mouse anti-Zap70 (pY292)


Product Information

Material Number:	558515
Alternate Name:	ZAP-70 (pY292)
Size:	50 tests
Vol. per Test:	20 µl
Clone:	J34-602
Immunogen:	Phosphorylated Human ZAP70
Isotype:	Mouse (BALB/c) IgG1, κ
Reactivity:	QC Tested: Human
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The 70-kDa ζ chain-associated protein (ZAP70) is a Syk-family protein tyrosine kinase (PTK) that associates with the ζ subunit of the T cell antigen receptor (TCR) and undergoes tyrosine phosphorylation following TCR stimulation. ZAP70 contains two SH2-like domains with the PTK domain located at the C-terminus. TCR-mediated Lck activity leads to phosphorylation of ZAP70 on Tyrosine 493 in the regulatory loop of the PTK domain leading to upregulation of ZAP70 kinase activity. Tyrosine 292 (Y292), in the linker region between the SH2 and PTK domains, is autophosphorylated by the activated PTK domain. By binding with c-Cbl, the phosphorylated Y292 can down-regulate TCR signaling.

The J34-602 antibody recognizes ZAP70 phosphorylated at Y292.

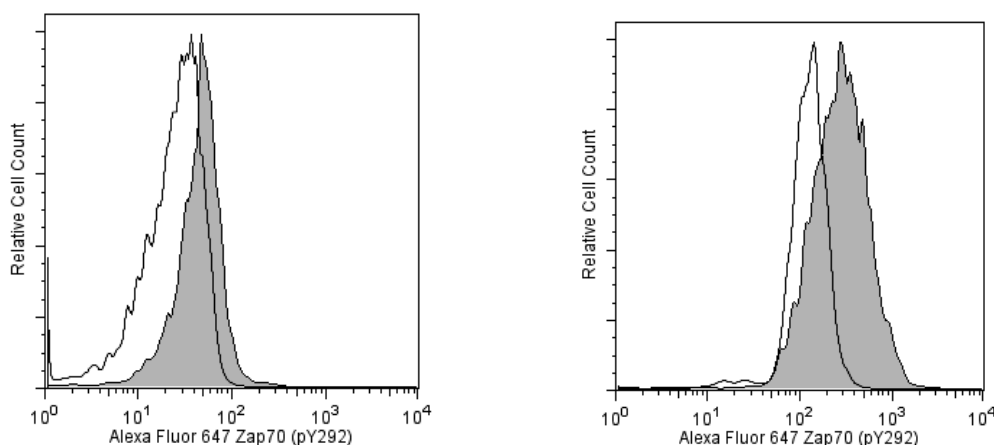

LEFT: Analysis of ZAP70 (pY292) in activated human T lymphocytes. Human whole blood was either stimulated by cross-linking of CD3 with NA/LE Mouse anti-Human CD3 mAb UCHT1 (Cat. No. 555329) on ice for 15 minutes followed by Purified Goat anti-Mouse Ig (Cat. No. 553998) on ice for 15 minutes, and then allowed to undergo phosphorylation at 37°C for 5 minutes (shaded histogram), or unstimulated (open histogram). The cells were lysed and fixed with 1X BD Phosflow™ Lyse/Fix Buffer (Cat. No. 558049) for 10-15 minutes at 37°C, permeabilized (BD Phosflow™ Perm Buffer III, Cat. No. 558050) on ice for 30 minutes, blocked with normal mouse immunoglobulin, and then stained with Alexa Fluor® 647 Mouse anti-ZAP70 (pY292) and PerCP-Cy5.5 Mouse anti-human CD3 mAb SK7 (Cat. No. 340949). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system. For data analysis, CD3-positive and -negative lymphocytes were selected by their scatter and staining profiles. The figure displays the CD3-positive T lymphocytes. Up-regulated phosphorylation of ZAP70 (pY292) was not observed in the CD3-negative lymphocytes (data not shown).

RIGHT: Analysis of ZAP70 (pY292) in activated human T leukemia cells. Jurkat cells (ATCC TIB152) were either stimulated by cross-linking of CD3 and CD28 with NA/LE Mouse anti-Human CD3 mAb UCHT1 (Cat. No. 555329) and NA/LE Mouse anti-Human CD28 mAb CD28.2 (Cat. No. 555725) on ice for 15 minutes followed by Purified Goat anti-Mouse Ig (Cat. No. 553998) on ice for 15 minutes, and then allowed to undergo phosphorylation at 37°C for 5 minutes (shaded histogram), or unstimulated (open histogram). The cells were fixed (BD Cytofix™ buffer, Cat. No. 554655) for 10 minutes at 37°C, permeabilized (BD Phosflow™ Perm Buffer III, Cat. No. 558050) on ice for at least 30 minutes, blocked with normal mouse immunoglobulin, and then stained with Alexa Fluor® 647 Mouse anti-ZAP70 (pY292). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Recommended Assay Procedure:

This antibody conjugate is suitable for intracellular staining of human whole blood (using BD Phosflow™ Lyse/Fix Buffer) and peripheral blood mononuclear cells (using BD Cytotfix™ Fixation Buffer). Any of the three BD Phosflow™ permeabilization buffers may be used.

Suggested Companion Products

Catalog Number	Name	Size	Clone
558049	Lyse/Fix Buffer 5X	250 ml	(none)
558050	Perm Buffer III	125 ml	(none)
554655	Fixation Buffer	100 ml	(none)
558052	Perm Buffer II	125 ml	(none)
557885	Perm/Wash Buffer I	125 ml	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-μl experimental sample (a test).
2. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
3. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
4. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
8. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

Kong G, Dalton M, Wardenburg JB, Straus D, Kurosaki T, Chan AC. Distinct tyrosine phosphorylation sites in ZAP-70 mediate activation and negative regulation of antigen receptor function. *Mol Cell Biol*. 1996; 16:5026-5035. (Biology)

Magnan A, Di Bartolo V, Mura AM, et al. T cell development and T cell responses in mice with mutations affecting tyrosines 292 or 315 of the ZAP-70 protein tyrosine kinase. *J Exp Med*. 2001; 194:491-505. (Biology)

Rao N, Lupher ML Jr, Ota S, Reedquist KA, Druker BJ, Band H. The linker phosphorylation site Tyr292 mediates the negative regulatory effect of Cbl on ZAP-70 in T cells. *J Immunol*. 2000; 164:4616-4626. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

