

Technical Data Sheet

Alexa Fluor® 647 Mouse anti-FAK (pS910)


Product Information

Material Number:	558545
Size:	50 tests
Vol. per Test:	20 µl
Clone:	K73-480
Immunogen:	Phosphorylated Human FAK
Isotype:	Mouse (BALB/c) IgG2b, κ
Reactivity:	QC Testing: Human Predicted: Mouse, Rat
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Focal Adhesion Kinase (FAK) is a cytoplasmic tyrosine kinase that associates with integrins in focal adhesions. Its cellular localization is directed by a 125-amino acid sequence at the C terminus called the "Focal Adhesion Targeting" (FAT) domain, and the phosphorylation state of serine 910 (S910) in the FAT domain may regulate the assembly of focal adhesions. Furthermore, the binding of extracellular matrix ligands to integrins triggers tyrosine phosphorylations near FAK's kinase domain that increase its kinase activity, and additional tyrosine phosphorylations near proline-rich motifs create binding sites for the SH2 domains of various adaptor proteins. FAK's ability to bind numerous structural and signaling proteins via a variety of interactions regulates FAK's targeting to focal adhesions, modulates its kinase activity, and initiates intracellular signaling cascades. Thus, studies suggest that FAK may integrate cellular events controlling cell motility, growth, and invasiveness.

The K73-480 monoclonal antibody recognizes the phosphorylated S910 of human FAK. The orthologous phosphorylation sites in mouse and rat FAK are S948 and S913, respectively.

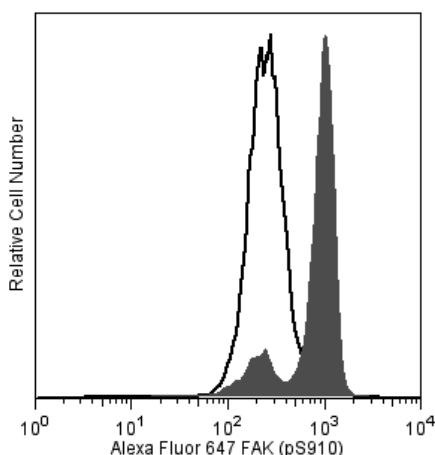

Analysis of FAK (pS910) in human epithelioid carcinoma. *Hela S3 cells (ATCC CCL 2.2) were either stimulated with Nocodazole at 37°C for 16 hours (shaded histogram) or unstimulated (open histogram). The cells were fixed (BD Cytotfix™ buffer, Cat. No. 554655) for 10 minutes at 37°C, then permeabilized (BD™ Phosflow Perm Buffer III, Cat. No. 558050) on ice for at least 30 minutes, and then stained with Alexa Fluor® 647 Mouse anti-FAK (pS910). Flow cytometry was performed on a BD™ FACSCalibur flow cytometry system.*

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Suggested Companion Products

Catalog Number	Name	Size	Clone
558050	Perm Buffer III	125 ml	(none)
554655	Fixation Buffer	100 ml	(none)
554657	Stain Buffer (BSA)	500 ml	(none)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
3. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
4. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
8. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

- Hunger-Glaser I, Perez Salazar E, Sinnott-Smith J, Rozengurt E. Bombesin, lysophosphatidic acid, and epidermal growth factor rapidly stimulate focal adhesion kinase phosphorylation at Ser-910. Requirement for ERK activation. *J Biol Chem.* 2003; 278(25):22631-22643. (Biology)
- Ma A, Richardson A, Schaefer EM, Parsons JT. Serine phosphorylation of focal adhesion kinase in interphase and mitosis: A possible role in modulating binding to p130Cas. *Mol Biol Cell.* 2001; 12:1-12. (Biology)
- Schlaepfer DD, Mitra SK, Ilic D. Control of motile and invasive cell phenotypes by focal adhesion kinase. *Biochim Biophys Acta.* 2004; 1692:77-102. (Biology)
- Yamakita Y, Totsukawa G, Yamashiro S, et al. Dissociation of FAK/c-Src complex during mitosis: Role of mitosis-specific serine phosphorylation of FAK. *J Cell Biol.* 1999; 144(2):315-324. (Biology)