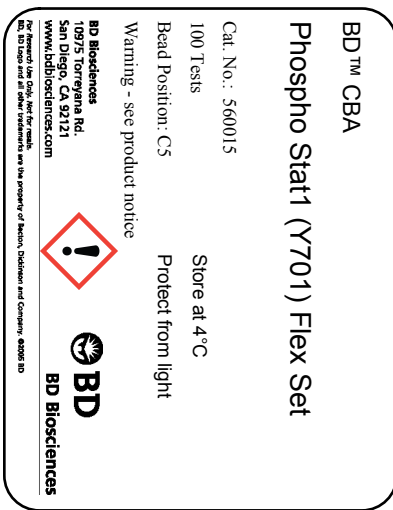




Technical Data Sheet

Phospho Stat1 (Y701) Flex Set

Product Information

Material Number:	560015
Size:	100 Tests
Bead Position:	C5
Assay Range:	3.9-1,000 Units/mL
Reactivity:	QC Testing: Human Tested in Development: Mouse

Component Description:	Phospho Stat1 (Y701) PE Detection Reagent
Component Mat. No:	51-9003792
Component Storage Buffer:	Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.
Component Description:	Phospho Stat1 (Y701) Standard
Component Mat. No:	51-9005201
Component Storage Buffer:	Lyophilized in an aqueous buffered solution containing BSA and ProClin™ 150.
Component Description:	Phospho Stat1 (Y701) Capture Bead C5
Component Mat. No:	51-9005203
Component Storage Buffer:	Aqueous buffered solution containing fetal bovine serum and $\leq 0.09\%$ sodium azide.

Description

The BD™ CBA Phospho Stat1 (Y701) Flex Set is a bead-based immunoassay capable of measuring human or mouse signal transducer and activator of transcription 1 (Stat1) that has been tyrosine-phosphorylated at Tyr-701 in denatured cell lysate samples. Human and mouse reactivity was determined by testing cell lysates with the BD CBA Phospho Stat1 (Y701) Flex Set. The biology and function of Stat1 has been previously reviewed. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Cell Signaling Master Buffer Kit (Cat. No. 560005 or 560006).

Preparation and Storage

This BD™ CBA Flex Set contains one vial of each component listed above. All components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. Store at 4°C. Protect Capture Beads and the PE Detection Reagent from prolonged exposure to light.

The Phospho Stat1 (Y701) Standard provided in this Flex Set is lyophilized and the standard sphere should be transferred to a 1.5 mL microfuge tube for reconstitution. Reconstitute the standard with 0.1 mL Assay Diluent from the BD CBA Cell Signaling Master Buffer Kit (Cat. No. 560005/560006), warm to 37 °C and vortex prior to use. After reconstitution, the standard concentration is 50,000 Units/mL and is stable for 3 months when stored at 4 °C. When using reconstituted standard, warm to 37 °C and vortex to mix thoroughly before use.

Application Notes

Recommended Assay Procedure: The BD CBA Phospho Stat1 (Y701) Flex Set must be used in conjunction with a BD CBA Cell Signaling Master Buffer Kit (Cat. No. 560005, 100 tests, or 560006, 500 tests) and a flow cytometer. Detailed instructions on the use of this product can be found in the manual for the BD CBA Cell Signaling Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Phospho Stat1 (Y701) Flex Set will be 1,000 Units/mL. An example standard curve is shown in Figure 1.

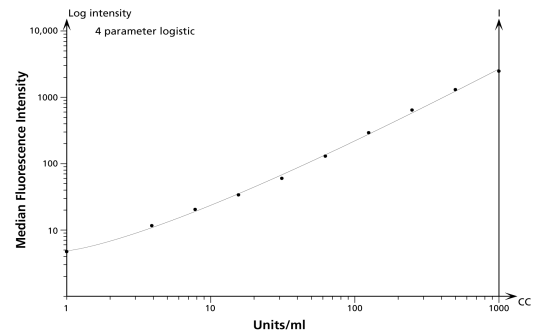


Figure 1. Example BD CBA Phospho Stat1 (Y701) Flex Set standard curve. Data acquired on a BD FACSAria bioanalyzer and analyzed using FCAP Array Software.

BD Biosciences

bdbiosciences.com					
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



The BD CBA Phospho Stat1 (Y701) Flex Set should not be used in the same assay well with any non-BD CBA Cell Signaling Flex Set reagents (such as BD CBA Human or Mouse Soluble Protein Flex Sets) **nor with any BD CBA Total Protein Cell Signaling Flex Set reagents**. For an updated assay compatibility chart for the BD CBA Cell Signaling Flex Sets, please refer to the BD CBA Flex Set System homepage at <http://www.bdbiosciences.com/cbasetup>.

Performance

Limit of Detection: The theoretical limit of detection is 1.36 Units/mL and was determined by evaluating the estimated result of the average MFI of the negative control (0 Units/mL, $n=30$) + 2 standard deviations.

Specificity		Inter-Assay Reproducibility			Intra-Assay Reproducibility		
		Mean (Units/ml)	Standard Deviation	%CV	Mean (Units/ml)	Standard Deviation	%CV
Phospho Stat1 (Y701)	Sample 1	30.2	2.4	8%	26.1	1.5	6%
	Sample 2	115.1	8.1	7%	104.7	3.2	3%
	Sample 3	492.6	21.1	4%	470.7	17.2	4%

Reproducibility: The inter-assay and intra-assay reproducibility were determined for the BD CBA Phospho Stat1 (Y701) Flex Set by evaluating ten replicates of three different sample levels (intra-assay) and two replicates of three different sample levels from four separate experiments (inter-assay).

Lysate Dilution	Phospho Stat1 (Y701)	
	Detected (Units/ml)	% of Expected
Neat	663.4	100%
1 : 2	334.3	101%
1 : 4	174.3	105%
1 : 8	88.4	107%

Linearity: An activated cell lysate was serially diluted to determine the linearity of the assay.

Product Notices

- ProClin is a trademark of Rohm and Haas Company.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- Warning: CBA lyophilized standard contains 32.07% sodium dodecyl sulfate (w/w) and 0.01% (w/w) of a CMIT/MIT mixture (3:1), which is a mixture of: 5-chloro-2-methyl-4-isothiazolin-3-one [EC No 247-500-7] and 2-methyl-4-isothiazolin-3-one [EC No 220-239-6] (3:1). Hazard statement: May be harmful if swallowed. Causes skin irritation. Causes serious eye irritation. May cause an allergic skin reaction. Precautionary statements: Wear protective gloves/eye protection. Wear protective clothing. Avoid breathing mist/vapours/spray. If skin irritation or rash occurs: Get medical advice/attention. IF ON SKIN: Wash with plenty of water. Dispose of contents/container in accordance with local/regional/national/international regulations.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 800.268.5430 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

