


## Technical Data Sheet

# Human Soluble TNFRII Flex Set

### Product Information

<b>Material Number:</b>	<b>560155</b>
<b>Size:</b>	100 Tests
<b>Bead Position:</b>	C5
<b>Assay Range:</b>	10-2,500 pg/mL
<b>Reactivity:</b>	QC Testing: Human

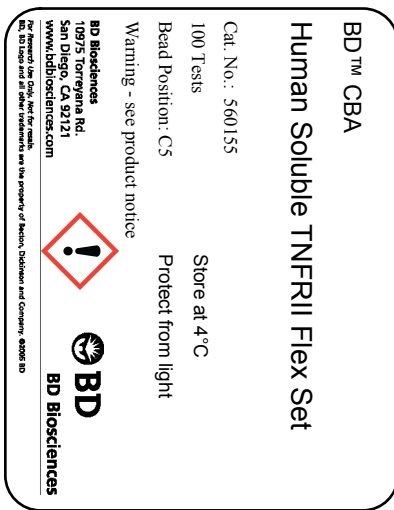

Component Description: Human Soluble TNFRII Standard  
 Component Mat. No: 51-9005507  
 Component Storage Buffer: Lyophilized in an aqueous buffered solution containing BSA and ProClin™ 150.

Component Description: Human Soluble TNFRII PE Detection Reagent  
 Component Mat. No: 51-9005605  
 Component Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Component Description: Human Soluble TNFRII Capture Bead C5  
 Component Mat. No: 51-9005607  
 Component Storage Buffer: Aqueous buffered solution containing fetal bovine serum and ≤0.09% sodium azide.

### Description


The BD™ CBA Human TNFRII Flex Set is a bead-based immunoassay capable of measuring soluble human type 2 tumor necrosis factor receptor (TNFRII), also known as CD120b, in serum, plasma, and cell culture supernatant samples. Human reactivity was determined by testing samples with the BD CBA Human TNFRII Flex Set. The biology and function of TNFRII has been extensively reviewed in the literature. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Human Soluble Protein Master Buffer Kit (Cat. No. 558264 or 558265).

This BD™ CBA Flex Set contains one vial each of Capture Bead and PE Detection Reagent and two vials of Standard. The Capture Bead and PE Detection Reagent components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. The Standard component is lyophilized and should be transferred to a 15 ml polypropylene tube for reconstitution. When reconstituted in 4.0 ml Assay Diluent, the standard has a protein concentration of 2,500 pg/ml. Discard unused reconstituted standard, do not store or reuse. Store lyophilized standard and other components at 4°C. Protect PE Detection Reagent from prolonged exposure to light.

### Application Notes

**Recommended Assay Procedure:** The BD CBA Human TNFRII Flex Set must be used in conjunction with a BD CBA Human Soluble Protein Master Buffer Kit (Cat. No. 558264, 100 tests, or 558265, 500 tests), a flow cytometer, and the FCAP Array™ Software (Cat. No. 641488). Detailed instructions on the use of this product can be found in the manual for the BD CBA Human Soluble Protein Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Human TNFRII Flex Set will be 2,500 pg/ml. An example standard curve is shown in Figure 1.

The BD CBA Human TNFRII Flex Set should not be used in the same assay well with any non-BD CBA Human Soluble Receptor or CD Marker Flex Set reagents (such as BD CBA Mouse Soluble Protein or Cell Signaling Flex Sets) **nor with any BD CBA Human Soluble Protein Flex**

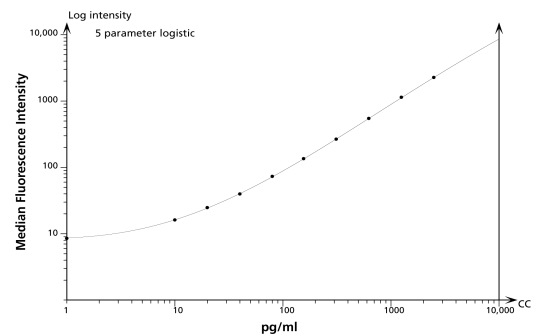

**Figure 1. Example BD CBA Human Soluble TNFRII Flex Set standard curve.** Data acquired on a BD FACSAry bioanalyzer and analyzed using the FCAP Array Software (Cat. No. 641488).

### BD Biosciences

bdbiosciences.com

<b>United States</b>	<b>Canada</b>	<b>Europe</b>	<b>Japan</b>	<b>Asia Pacific</b>	<b>Latin America/Caribbean</b>
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

*Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.*

For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



**Set reagents.** For an updated assay compatibility chart for the BD CBA Human Soluble Protein Flex Sets, please refer to the BD CBA Flex Set System homepage at <http://www.bdbiosciences.com/flexset>.

**Performance**

**Limit of Detection:** The theoretical limit of detection is 1.4 pg/ml and was determined by evaluating the estimated result of the average MFI of the negative control (0 pg/ml,  $n=30$ ) + 2 standard deviations.

Specificity		Inter-Assay Reproducibility			Intra-Assay Reproducibility		
		Mean (pg/mL)	Standard Deviation	%CV	Mean (pg/mL)	Standard Deviation	%CV
Human Soluble TNFRII	Sample 1	38.0	3.0	8%	34.2	0.5	2%
	Sample 2	150.4	9.5	6%	137.7	2.3	2%
	Sample 3	605.8	33.3	5%	590.7	16.7	3%

**Reproducibility:** The intra-assay and inter-assay reproducibility were determined for the BD CBA Human Soluble TNFRII Flex Set by evaluating ten replicates of three different sample levels (intra-assay) and three replicates of three different sample levels from four separate experiments (inter-assay).

Specificity	Cell Culture Supernatant		Serum		Plasma	
	Average % Recovery	Range	Average % Recovery	Range	Average % Recovery	Range
Human Soluble TNFRII	90%	88 - 91%	85%	65 - 99%	91%	83 - 96%

**Recovery:** Cell culture supernatant, serum, or EDTA-treated plasma were spiked with three different levels of protein. The spiked samples were assayed and the results were compared with expected values. Serum and plasma samples were diluted 1:4 before the protein was spiked into each. Serum is a pool of 800 - 1000 donors and the plasma was pooled from at least 20 donors.

Sample Dilution	Cell Culture Supernatant		Serum		Plasma	
	Detected (pg/mL)	% of Expected	Detected (pg/mL)	% of Expected	Detected (pg/mL)	% of Expected
Spiked sample	549.3	100%	631.1	100%	584.3	100%
1 : 2	282.2	103%	280.3	89%	265.8	91%
1 : 4	143.9	105%	126.3	80%	124.3	85%

**Linearity:** Cell culture supernatant, 1:4 diluted serum, or 1:4 diluted EDTA-treated plasma were spiked with protein and serially diluted. The diluted samples were assayed and the results were compared with the original spiked sample.

**Product Notices**

1. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
2. ProClin is a trademark of Rohm and Haas Company.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Warning: CBA lyophilized standard contains 0.02% (w/w) of a CMIT/MIT mixture (3:1), which is a mixture of: 5-chloro-2-methyl-4-isothiazolin-3-one [EC No 247-500-7] and 2-methyl-4-isothiazolin-3-one [EC No 220-239-6] (3:1).  
Hazard statement: May cause an allergic skin reaction.  
Precautionary statements: Wear protective gloves/eye protection. Wear protective clothing. Avoid breathing mist/vapours/spray. If skin irritation or rash occurs: Get medical advice/attention. IF ON SKIN: Wash with plenty of water. Dispose of contents/container in accordance with local/regional/national/international regulations.

**BD Biosciences**

[bdbiosciences.com](http://bdbiosciences.com)

United States 877.232.8995    Canada 800.268.5430    Europe 32.2.400.98.95    Japan 0120.8555.90    Asia Pacific 65.6861.0633    Latin America/Caribbean 55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

*Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.*

For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

