# **Technical Data Sheet**

# APC Rat anti-Mouse CD229.1

#### **Product Information**


**Material Number:** 560188 Alternate Name: Ly-9.1 0.1 mg Size 0.2 mg/ml Concentration: 30C7 Clone:

Immunogen: SJL/J mouse splenocytes Isotype: Rat (LOU) IgG2a, ĸ Reactivity: QC Testing: Mouse

Storage Buffer: Aqueous buffered solution containing ≤0.09% sodium azide.

#### Description

The 30C7 antibody reacts with CD229.1, the Ly-9.1 (Lgp-100a) surface glycoprotein alloantigen found on most thymocytes, peripheral T and B cells, bone marrow lymphoid cells, hematopoietic progenitors, and lymphomas of most strains (e.g., AKR, Balb/c, CBA, SJL). 30C7 mAb does not react with the Ly-9.2 alloantigen expressed on mice derived from the C57/J line (e.g., C57BL, C57BR, C57L, C58).

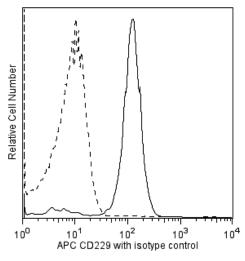

#### Flow cytometric analysis of APC-conjugated anti-mouse CD229.1 on mouse splenocytes.

Splenocytes from BALB/c mice were stained with APC anti-mouse CD229.1 (clone 30C7, Cat. No. 560188) and compared to splenocytes stained with an APC rat IgG2a isotype control (clone R35-95, Cat. No. 553932). The isotype control is represented by a dashed line and the APC anti-mouse CD229.1 by the solid line. Flow cytometry was performed on a BD FACSCalibur™ System.

#### **Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

#### **Application Notes**

Application

Flow cytometry Routinely Tested

## **Suggested Companion Products**

Catalog Number	Name	Size	Clone
553932	APC Rat IgG2a κ Isotype Control	0.1 mg	R35-95

### **Product Notices**

- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser. 3.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

#### References

## **BD Biosciences**

bdbiosciences.com United States Asia Pacific Europe 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how\_to\_order/ Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation

cof any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Lanier LL, Warner NL, Ledbetter JA, Herzenberg LA. Quantitative immunofluorescent analysis of surface phenotypes of murine B cell lymphomas and plasmacytomas with monoclonal antibodies. *J Immunol.* 1981; 127(4):1691-1697. (Clone-specific)

Ledbetter JA, Goding JW, Tsu TT, Herzenberg LA. A new mouse lymphoid alloantigen (Lgp100) recognized by a monoclonal rat antibody. *Immunogenetics.* 1979; 8(1):347-360. (Immunogen)

Ledbetter JA, Herzenberg LA. Xenogeneic monoclonal antibodies to mouse lymphoid differentiation antigens. Immunol Rev. 1979; 47:63-90. (Biology)

560188 Rev. 1 Page 2 of 2