BD™ Cytometric Bead Array (CBA)



Technical Data Sheet Mouse Soluble CD62L (L-Selectin) Flex Set

Product Information

- Material Number: Size: Bead Position: Assay Range: Reactivity:
- 560385 100 Tests E7 40-10,000 pg/mL QC Testing: Mouse

Component Description: Component Mat. No: Component Storage Buffer:

Component Description: Component Mat. No: Component Storage Buffer:

Component Description: Component Mat. No: Component Storage Buffer: Mouse Soluble CD62L (L-Selectin) Standard 51-9006074 Lyophilized in an aqueous buffered solution containing BSA and ProClin[™] 150. Mouse Soluble CD62L (L-Selectin) PE Detection Reagent 51-9006089

Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Mouse Soluble CD62L (L-Selectin) Capture Bead E7 51-9006091 Aqueous buffered solution containing fetal bovine serum and ≤0.09% sodium azide.

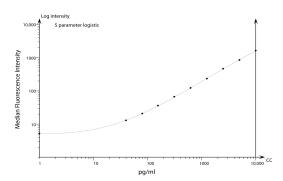


Figure 1. Example BD CBA Mouse Soluble CD62L (L-Selectin) Flex Set standard curve. Data acquired on a BD FACSArray bioanalyzer and analyzed using FCAP Array Software.

Description

The BDTM CBA Mouse Soluble CD62L (L-Selectin) Flex Set is a bead-based immunoassay capable of measuring mouse soluble CD62L (L-Selectin) in cell culture supernatant samples. Mouse reactivity was determined by testing samples with the BD CBA Mouse Soluble CD62P (L-Selectin) Flex Set. The biology and function of CD62L (L-Selectin) has been extensively reviewed in the literature. For more information on bead-based immunoassays, refer to the product insert for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266/558267).

Preparation and Storage

This BD[™] CBA Flex Set contains one vial each of Capture Bead and PE Detection Reagent and two vials of Standard. The Capture Bead and PE Detection Reagent components of this flex set have been formulated to a 50x concentration to ensure product performance when multiplexed. The Standard component is lyophilized and should be transferred to a 15 mL polypropylene tube for reconstituted standard, do not reconstituted in 4.0 mL Assay Diluent, the standard has a protein concentration of 10,000 pg/mL. Discard unused reconstituted standard, do not store or reuse. Store lyophilized standard and other components at 4°C. Protect Capture Beads and the PE Detection Reagent from prolonged exposure to light.

Application Notes

Recommended Assay Procedure: The BD CBA Mouse Soluble CD62L (L-Selectin) Flex Set must be used in conjunction with a BD CBA Mouse/Rat Soluble Protein Master Buffer Kit (Cat. No. 558266, 100 tests, or 558267, 500 tests), a flow cytometer, and FCAP Array[™] Software. Detailed instructions on the use of this product can be found in the manual for the BD CBA Mouse/Rat Soluble Protein Master Buffer Kit. When following the directions in the Master Buffer Kit, the top standard point for the BD CBA Mouse Soluble CD62L (L-Selectin) Flex Set will be 10,000 pg/mL. An example standard curve is shown in figure 1.

The BD CBA Mouse Soluble CD62L (L-Selectin) Flex Set should not be used in the same assay well with any non-BD CBA Mouse Soluble

BD Biosciences

bdbiosciences.	om				
United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbea
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995
For country cor	ntact informatio	on, visit bdbiosci e	ences.com/conta	ct	
of any patents. BL	D Biosciences will n	ot be held responsik	ble for patent infrin	gement or other vio	e the above product in violation
use of our product	ts. Purchase does n	ot include or carry a	any right to resell or	r transfer this produ	vlations that may occur with the
product or as a co	mponent of anoth	er product. Any use	of this product oth	er than the permitte	ct either as a stand-alone
written authoriza	tion of Becton, Dic	kinson and Compan	by is stictly prohibite	d.	ed use without the express
For Research Use	Only. Not for use ir	diagnostic or thera	apeutic procedures.	Not for resale.	an and Company. © 2014 BD

Protein Flex Set reagents (such as BD CBA Human or Rat Soluble Protein or Cell Signaling Flex Sets). For an updated assay compatibility chart for the BD CBA Mouse Soluble Protein Flex Sets, please refer to the BD CBA Flex Set System homepage at http://www.bdbiosciences.com/cbasetup.

Performance

Limit of Detection: The theoretical limit of detection is 8.3 pg/mL and was determined by evaluating the estimated result of the average MFI of the negative control (0 pg/mL, n=30) + 2 standard deviations.

	Inter-Assay Reproducibility			Intra-Assay Reproducibility			
Specificity	Specificity		Standard Deviation	%CV	Mean (pg/ml)	Standard Deviation	%CV
	Sample 1	158.3	13.1	8%	150.4	5.6	4%
Mouse sCD62L	Sample 2	609.7	17.7	3%	589.5	11.6	2%
(L-Selectin)	Sample 3	2449.3	85.9	4%	2338.4	47.0	2%

Reproducibility: The intra-assay and inter-assay reproducibility were determined for the BD CBA Mouse Soluble CD62L (L-Selectin) Flex Set by evaluating ten replicates of three different sample levels (intra-assay) and three replicates of three different sample levels from four separate experiments (inter-assay).

		ulture natant	Ser	um	Plasma		
Specificity	Average % Recovery	Range	Average % Recovery	Range	Average % Recovery	Range	
Mouse sCD62L (L-Selectin)	96%	94 - 98%	ND*	ND*	ND*	ND*	

Recovery: Cell culture supernatant and pooled mouse serum were spiked with three different levels of protein. The spiked samples were assayed and the results were compared with expected values. Serum samples were diluted 1:4 before the protein was spiked into each. Serum used was commercially available pooled mouse serum

*ND = Not Determined. Due to high levels of CD62L (L-Selectin) in serum and EDTA-plasma samples, recovery information could not be generated.

	Cell Culture Supernatant		Ser	um	Plasma		
Sample Dilution	Detected (pg/mL)	% of Expected	Detected (pg/mL)	% of Expected	Detected (pg/mL)	% of Expected	
Spiked sample	2329.6	100%	ND*	ND*	ND*	ND*	
1:2	1254.8	108%	ND*	ND*	ND*	ND*	
1:4	665.2	114%	ND*	ND*	ND*	ND*	

Linearity: Cell culture supernatant and 1:4 diluted pooled mouse serum were spiked with protein and serially diluted. The diluted samples were assayed and the results were compared with the original spiked sample

*ND = Not Determined. Due to high levels of CD62L (L-Selectin) in serum and EDTA-plasma samples, linearity information could not be generated.


Product Notices

- ProClin is a trademark of Rohm and Haas Company. 1.
- 2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before 3. discarding to avoid accumulation of potentially explosive deposits in plumbing.

4. Warning: CBA lyophilized standard contains 0.02% (w/w) of a CMIT/MIT mixture (3:1), which is a mixture of: 5-chloro-2-methyl-4-isothiazolin-3-one [EC No 247-500-7] and 2-methyl-4-isothiazolin-3-one [EC No 220-239-6] (3:1). Hazard statement: May cause an allergic skin reaction.

Precautionary statements: Wear protective gloves/eye protection. Wear protective clothing. Avoid breathing mist/vapours/spray. If skin irritation or rash occurs: Get medical advice/attention. IF ON SKIN: Wash with plenty of water. Dispose of contents/container in accordance with local/regional/national/international regulations.

BD Biosciences bdbiosciences.com United States Canada Asia Pacific Latin America/Caribbean Europe Japan 800.268.5430 32.2.400.98.95 0120.8555.90 877.232.8995 65.6861.0633 55.11.5185.9995 For country contact information, visit bdbiosciences.com/contact Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.

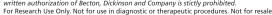

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

