

Technical Data Sheet


APC Mouse Anti-Human IL-4

Product Information

Material Number:	560671
Size:	50 tests
Vol. per Test:	5 µl
Clone:	8D4-8
Immunogen:	Recombinant Human IL-4
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The 8D4-8 monoclonal antibody reacts with human interleukin-4 (IL-4). The immunogen used to raise the 8D4-8 hybridoma was recombinant human IL-4. The 8D4-8 antibody binds to an epitope that is different than the epitope recognized by the MP4-25D2 antibody (Cat. No. 554485). Clone 8D4-8 displays an increased amount of non-specific binding to dead cells when compared to the clone MP4-25D2. It is recommended to use a fixable viability dye in conjunction with this clone.

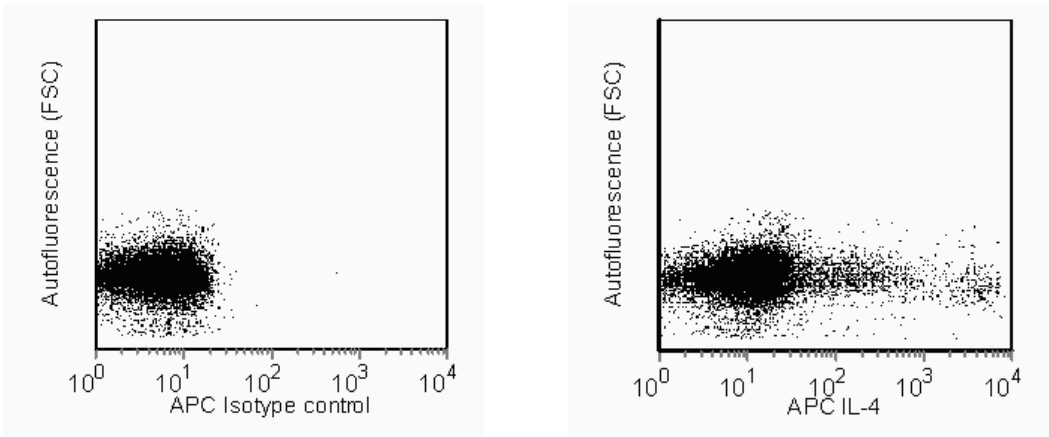

Flow cytometric analysis for IL-4 in stimulated human peripheral blood mononuclear cells (PBMC).
Human PBMC were stimulated for at least 6 hours with 50 ng/mL PMA (Sigma-Aldrich Cat. No. P-8139) and 250 ng/mL calcium ionophore A23187 (Sigma-Aldrich Cat. No. C-9275) in the presence of 2 µM BD GolgiStop™ (Cat. No. 554724). Cells were then fixed and permeabilized using BD Cytotfix/Cytoperm™ (Cat. No. 554714) followed by staining with either a APC Mouse IgG1, κ isotype control (left panel) or with the APC Mouse Anti-Human IL-4 antibody (right panel). Dot plots were derived from gated events based on light scattering characteristics for lymphocytes. Flow cytometry was performed on a BD™ LSR II flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.
The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.
The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Recommended Assay Procedure:

Flow cytometry: The 8D4-8 antibody is useful for immunofluorescent staining and flow cytometric analysis to identify and enumerate IL-4 producing cells within mixed cell populations. A useful control investigators may consider using for demonstrating specificity of staining, is to pre-block with one of the following reagents: (1) recombinant human IL-4 (Cat. No. 554605) or (2) unlabeled 8D4-8 antibody (Cat. No. 554515), prior to staining.

BD Biosciences
bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



Suggested Companion Products

Catalog Number	Name	Size	Clone
554681	APC Mouse IgG1 κ Isotype Control	0.1 mg	MOPC-21
554714	BD Cytofix/Cytoperm™ Fixation/Permeabilization Kit	250 tests	(none)
554724	Protein Transport Inhibitor (Containing Monensin)	0.7 ml	(none)
555062	HiCK-2 Human Cytokine Positive Control Cells	1.0 ml	(none)
554605	Recombinant Human IL-4	5 μ g	(none)
554515	Purified Mouse Anti-Human IL-4	0.5 mg	8D4-8

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
5. This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

Bird C, Wadhwa M, Thorpe R. Development of immunoassays for human interleukin 3 and interleukin 4, some of which discriminate between different recombinant DNA-derived molecules. *Cytokine*. 1991; 3(6):562-567. (Biology)

Prussin C, Metcalfe DD. Detection of intracytoplasmic cytokine using flow cytometry and directly conjugated anti-cytokine antibodies. *J Immunol Methods*. 1995; 188(1):117-128. (Methodology: Flow cytometry)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

