

Technical Data Sheet

V500 Mouse IgG1, κ Isotype Control**Product Information**


Material Number:	560787
Size:	0.1 mg
Concentration:	0.2 mg/ml
Clone:	X40
Immunogen:	Keyhole limpet hemocyanin (KLH)
Isotype:	Mouse (BALB/c) IgG1, κ
Reactivity:	QC Testing: Human and rat
Storage Buffer:	Aqueous buffered solution containing protein stabilizer, glycerol and ≤0.09% sodium azide.

Description

The X40 hybridoma was derived from the fusion of mouse Sp2/0-Ag 14 myeloma cells with spleen cells from BALB/c mice immunized with keyhole limpet hemocyanin (KLH). The hybridoma secretes monoclonal antibody that binds specifically to KLH, an antigen not expressed on human cells or human cell lines.

The antibody is conjugated to BD Horizon™ V500, which has been developed for use in multicolor flow cytometry experiments and is available exclusively from BD Biosciences. It is excited by the Violet laser with an Ex max of 415 nm and Em Max at 500 nm. BD Horizon V500 conjugates emit at a similar wavelength to Amcyan yet exhibit reduced spillover into the FITC channel. For more information on BD Horizon V500, visit bdbiosciences.com/colors.

When compensating dyes in this spectral range (such as Horizon™ V500 and AmCyan), the most accurate compensation can be obtained using single stained cellular controls. Due to spectral differences between cells and beads in this channel, using BD CompBeads can result in spillover errors for V500 and AmCyan reagents. Therefore, the use of BD CompBeads or BD CompBeads Plus to determine spillover values for these reagents is not recommended. Different V500 reagents (e.g. CD4 vs. CD45) can have slightly different fluorescence spillover therefore, it may also be necessary to use clone specific compensation controls when using these reagents.

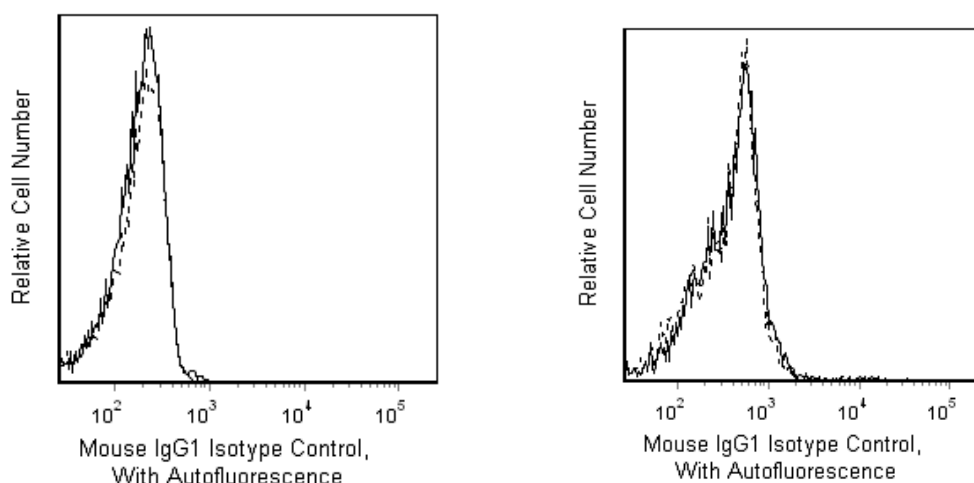

Profile of BD Horizon™ V500 Mouse IgG1, κ Isotype Control (X40) staining on human lysed whole blood and rat bone marrow. Human whole blood (left panel) or bone marrow cells from a Lew rat (right panel) were stained with BD Horizon™ V500 Mouse IgG1, κ Isotype Control (solid line) or left unstained (dotted line). Flow Cytometry was performed on a BD FACSCanto™ II flow cytometry system.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ V500 under optimum conditions, and unreacted BD Horizon™ V500 was removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Flow cytometry	Routinely Tested
Isotype control	Routinely Tested

Product Notices

1. Please refer to www.bdbiosciences.com/pharming/en/protocols for technical protocols.
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
4. BD Horizon™ V500 has a maximum absorption of 415 nm and maximum emission of 500 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

