


## Technical Data Sheet

**FITC Mouse Anti-Human CD144****Product Information**

<b>Material Number:</b>	<b>560874</b>
<b>Alternate Name:</b>	VE-cadherin; Cadherin-5; CDH5; Vascular endothelial cadherin
<b>Size:</b>	25 tests
<b>Vol. per Test:</b>	20 µl
<b>Clone:</b>	55-7H1
<b>Isotype:</b>	Mouse IgG1, κ
<b>Reactivity:</b>	QC Testing: Human
<b>Workshop:</b>	VI E092
<b>Storage Buffer:</b>	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

**Description**

The 55-7H1 antibody reacts with a calcium-independent epitope on cadherin-5, a member of the cadherin family of calcium-dependent adhesion molecules. Cadherin-5 is expressed on endothelial cells *in vivo* and *in vitro*. It may play a role in the organization of lateral endothelial junctions and in the control of permeability properties of vascular endothelium.

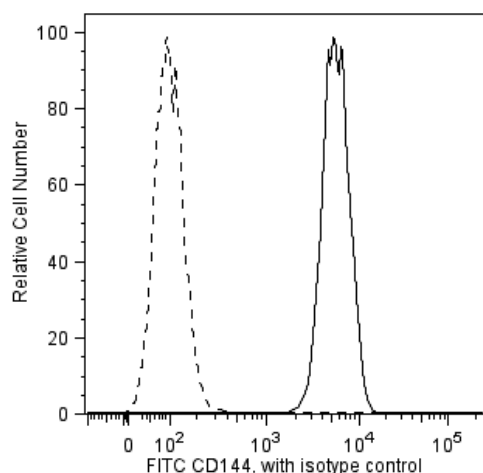

**Analysis of CD144 on Human Umbilical Vein Endothelial Cells (HUVECs).** HUVECs were stained with either FITC Mouse anti-Human CD144 (solid line) or FITC Mouse IgG1, κ Isotype Control (clone MOPC-21, Cat. No. 555748, dotted line). Flow cytometry was performed on a BD LSR™ II flow cytometry system.

**Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed.

**Application Notes****Application**

Flow cytometry

Routinely Tested

**Suggested Companion Products**

Catalog Number	Name	Size	Clone
555748	FITC Mouse IgG1, κ Isotype Control	100 tests	MOPC-21

**Product Notices**

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
5. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).

**BD Biosciences**

[bdbiosciences.com](http://bdbiosciences.com)

United States 877.232.8995	Canada 800.979.9408	Europe 32.53.720.550	Japan 0120.8555.90	Asia Pacific 65.6861.0633	Latin America/Caribbean 55.11.5185.9995
-------------------------------	------------------------	-------------------------	-----------------------	------------------------------	--

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



6. Please refer to [www.bdbiosciences.com/pharmingen/protocols](http://www.bdbiosciences.com/pharmingen/protocols) for technical protocols.

## References

- Breier G, Breviario F, Caveda L, et al. Molecular cloning and expression of murine vascular endothelial-cadherin in early stage development of cardiovascular system. *Blood*. 1996; 87(2):630-641. (Biology)
- Lampugnani MG, Resnati M, Raiteri M, et al. A novel endothelial-specific membrane protein is a marker of cell-cell contacts. *J Cell Biol*. 1996; 118(6):1511-1522. (Biology)
- Vincent PA, Xiao K, Buckley KM, Kowalczyk AP. VE-cadherin: adhesion at arm's length. *Am J Physiol Cell Physiol*. 2004; 286:C987-C997. (Biology)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

*Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.*

*For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.*

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

