Technical Data Sheet

PE Mouse Anti-Human CD152

Product Information

560939 **Material Number:** CTLA-4 Alternate Name: 25 tests Size: 20 µl Vol. per Test: BNI3 Clone:


Isotype: Mouse IgG2a, κ Reactivity: QC Testing: Human

Aqueous buffered solution containing BSA and ≤0.09% sodium azide. Storage Buffer:

Description

Reacts with the "cytolytic T lymphocyte-associated antigen", CTLA-4. CTLA-4 is transiently expressed on activated CD28+ T cells and binds to CD80 and CD86 present on antigen presenting cells (APC) with high avidity. This interaction appears to deliver a negative regulatory signal to the T cell. There are recent reports that indicate that CTLA-4 is also expressed on B cells when cultured with activated T cells, suggesting a possible role of CTLA-4 in the regulation of B-cell response. Immobilized BNI3.1 enhances T-cell proliferation induced by CD3 and CD28.

Recent studies have showed that CD152 can be expressed by regulatory T (Treg) cells. It has been found this antibody can stain the intracellular CD152 on the Treg cells after fixation and permeabilization of cells.

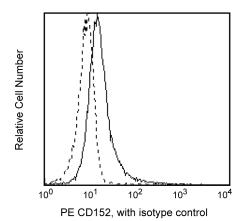

Profile of Concanavalin A-stimulated peripheral blood mononuclear cells analyzed on a FACScan (BDIS, San Jose, CA).

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

Application Notes

Αı	nr	olia	at	ion

Flow cytometry	Routinely Tested
Intracellular staining (flow cytometry)	Tested During Development

Suggested Companion Products

Catalog Number	Name	Size	Clone
555574	PE Mouse IgG2a, κ Isotype Control	100 tests	G155-178
554714	BD Cytofix/Cytoperm™ Fixation/Permeablization Kit	250 tests	(none)

BD Biosciences

bdbiosciences.com

United States 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.
For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Page 1 of 2 560939 Rev. 1

Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-μl experimental sample (a test).
- 2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 5. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997. (Biology)

Kuiper HM, Brouwer M, Linsley PS, van Lier RA. Activated T cells can induce high levels of CTLA-4 expression on B cells. *J Immunol*. 1995; 155(4):1776-1783.

(Riology)

Lindsten T, Lee KP, Harris ES, et al. Characterization of CTLA-4 structure and expression on human T cells. *J Immunol.* 1993; 151(7):3489-3499. (Biology) Morton PA, Fu XT, Stewart JA, et al. Differential effects of CTLA-4 substitutions on the binding of human CD80 (B7-1) and CD86 (B7-2). *J Immunol.* 1996; 156(3):1047-1054. (Biology)

560939 Rev. 1 Page 2 of 2