

Technical Data Sheet

V500 Mouse Anti-Human CD16

Product Information


Material Number:	561393
Alternate Name:	FcRIII; Fc-gamma RIII; FCG3; FCGR3; FCGRIII; FcγRIII; IGFR3
Size:	25 tests
Vol. per Test:	5 µl
Clone:	3G8
Immunogen:	Human polymorphonuclear leukocytes
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human
Workshop:	IV N409
Storage Buffer:	Aqueous buffered solution containing protein stabilizer, glycerol and ≤0.09% sodium azide.

Description

The 3G8 monoclonal antibody specifically binds to the 50-65 kDa transmembrane form of the IgG Fc Receptor (FcγRIII), a human NK-cell-associated antigen. CD16 is expressed on NK cells as well as macrophages and granulocytes. Reports indicate that CD16 plays a role in signal transduction and NK cell activation. The 3G8 antibody blocks the binding of soluble immune complexes to granulocytes.

The antibody is conjugated to BD Horizon™ V500, which has been developed for use in multicolor flow cytometry experiments and is available exclusively from BD Biosciences. It is excited by the Violet laser with an Ex max of 415 nm and Em Max at 500 nm. BD Horizon V500 conjugates emit at a similar wavelength to Amcyan yet exhibit reduced spillover into the FITC channel. For more information on BD Horizon V500, visit bdbiosciences.com/colors.

When compensating dyes in this spectral range (such as Horizon™ V500 and AmCyan), the most accurate compensation can be obtained using single stained cellular controls. Due to spectral differences between cells and beads in this channel, using BD CompBeads can result in spillover errors for V500 and AmCyan reagents. Therefore, the use of BD CompBeads or BD CompBeads Plus to determine spillover values for these reagents is not recommended. Different V500 reagents (e.g. CD4 vs. CD45) can have slightly different fluorescence spillover therefore, it may also be necessary to use clone specific compensation controls when using these reagents.

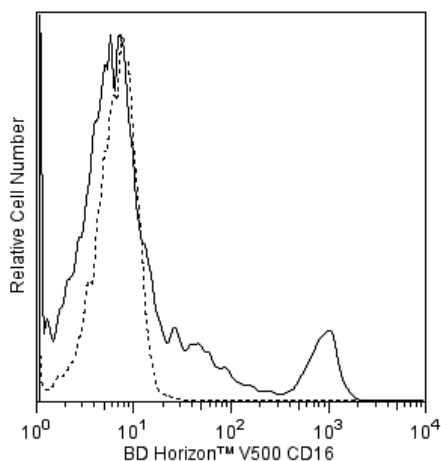

Flow cytometric analysis of CD16 expression on human peripheral blood lymphocytes. Human whole blood was stained with the BD Horizon™ V500 Mouse anti-Human CD16 antibody (Cat. No. 561394/561393; solid line histogram) or with a BD Horizon™ V500 Mouse IgG1, κ Isotype Control (Cat. No. 560787; dashed line histogram). The erythrocytes were lysed with BD PharmLyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ V500 under optimum conditions, and unreacted BD Horizon™ V500 was removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.268.5430	32.53.720.550	0120.8555.90	65.6861.0633	0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



Application Notes

Application

Flow cytometry

Routinely Tested

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
561394	V500 Mouse Anti-Human CD16	100 tests	3G8
560787	V500 Mouse IgG1, κ Isotype Control	0.1 mg	X40
555899	Lysing Buffer	100 ml	(none)
554656	Stain Buffer (FBS)	500 ml	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Please refer to www.bdbiosciences.com/pharming/en/protocols for technical protocols.
4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
5. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
6. BD Horizon™ V500 has a maximum absorption of 415 nm and maximum emission of 500 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.

References

- Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997. (Biology)
- Fleit HB, Wright SD, Unkeless JC. Human neutrophil Fc gamma receptor distribution and structure. *Proc Natl Acad Sci U S A*. 1982; 79(10):3275-3279. (Biology)
- Knapp W, Dorken B, Rieber EP, et al, ed. *Leucocyte Typing IV*. New York: Oxford University Press; 1989:1-1208. (Clone-specific)
- Stroncek DF, Skubitz KM, Plachta LB, et al. Alloimmune neonatal neutropenia due to an antibody to the neutrophil Fc-gamma receptor III with maternal deficiency of CD16 antigen. *Blood*. 1991; 77(7):1572-1580. (Biology)
- Wirthmueller U, Kurosaki T, Murakami MS, Ravetch JV.. Signal transduction by Fc gamma RIII (CD16) is mediated through the gamma chain. *J Exp Med*. 1992; 175(5):1381-1390. (Biology)