Technical Data Sheet


PE Mouse Anti-Human CD1d

Product Information

Material Number:	561757
Size:	25 tests
Vol. per Test:	20 µl
Clone:	CD1d42
Isotype:	Mouse IgG1, ĸ
Reactivity:	QC Testing: Human
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Reacts with a glycosylated type I transmembrane α chain (43-49 kDa) non- covalently associated with β 2-microglobulin. CD1d is a member of the CD1 family of molecules, which belong to the immunoglobulin superfamily. Sequence homology data classifies the CD1 molecules into two groups. Group 1 includes CD1a, CD1b and CD1c molecules; group 2 includes CD1d molecules and their homologues in other species. CD1d is expressed on human intestinal epithelial cells (IEC). CD1d42 monoclonal antibody reacts with peripheral blood B cells coexpressing CD19 and CD45RA. It is not seen on CD3[+], or CD45RO[+] cells. Studies suggest that CD1d participates in lipid antigen presentation.

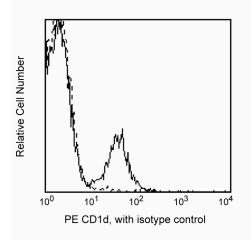

Profile of peripheral blood lymphocytes analyzed by flow cytometry

Preparation and Storage

Store undiluted at 4°C.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

Application Notes

Application							
Flow cytometry	Routinely Tested						
Suggested Companion Products							
<u>Catalog Number</u> 555749	<u>Name</u> PE Mouse IgG1, κ Isotype Control		<u>Size</u> 100 tests	Clone MOPC-21			
555747	i E Mouse igoi, k isotype control		100 10313	MOI C 21			

Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^{6} cells in a 100-µl experimental sample (a test).
- 2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 5. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- 6. An isotype control should be used at the same concentration as the antibody of interest.

BD Biosciences

bdbiosciences.com								
United States 877,232,8995	Canada 888 268 5430	Europe 32,53,720,550	Japan 0120.8555.90	Asia Pacific	Latin America/Caribbean 0800 771 7157			
For country-specific contact information, visit bdbiosciences.com/how_to_order/								
Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to result or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale. BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD								



References Hong S, Scherer DC, Singh N. Lipid antigen presentation in the immune system: lessons learned from CD1d knockout mice. *Immunol Rev.* 1999; 169:31-44. (Biology) Kishimoto T, von dem Borne AEG, Goyert SM,et al., ed. *Leucocyte Typing VI: White Cell Differentiation Antigens*. London: Garland Publishing; 1997. (Biology)