Technical Data Sheet

APC Mouse Anti-Human CD27

Product Information

Material Number: 561786

Alternate Name: TNFRSF7; TNF receptor superfamily, member 7; T14; Tp55; S152


Size 25 tests Vol. per Test: 20 µl M-T271 Clone:

Mouse (BALB/c) IgG1, κ Isotype: Reactivity: QC Tested: Human Workshop: V 5T CD27.03

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The M-T271 monoclonal antibody specifically reacts with CD27. CD27 presents as a type I transmembrane, disulphide-linked 110 kDa homodimer comprised of two polypeptide chains. The CD27 molecule is a lymphocyte-specific member of the TNF/NGF-R family, and is expressed on a subset of human thymocytes and on the majority of mature T lymphocytes, activated B cells and NK cells. CD27 is highly induced on T cells after TCR stimulation. CD27 binds to CD70 (also known as, CD27 ligand or CD27L) and may be involved in cellular interaction of T and B lymphocytes.

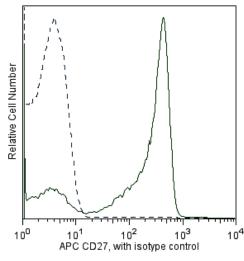

Flow cytometric analysis of anti-human CD27 on human whole blood lymphocytes. Human lymphocytes were stained with APC-anti human CD27 (solid line, Cat. No. 558664) or APC mouse IgG1 isotype control (dashed line, Cat. No. 555751/554681). Erythrocytes were lysed using BD Pharm Lyse™ lysis buffer (Cat. No. 555899). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone
554681	APC Mouse IgG1 κ Isotype Control	0.1 mg	MOPC-21
555751	APC Mouse IgG1, κ Isotype Control	100 tests	MOPC-21
554656	Stain Buffer (FBS)	500 ml	(none)
555899	Lysing Buffer	100 ml	(none)

BD Biosciences

bdbiosciences.com

United States Europe Japan 877.232.8995 888.268.5430 32.53.720.550 0120.8555.90 65.6861.0633 0800.771.7157

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2011 BD



561786 Rev. 1 Page 1 of 2

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-μl experimental sample (a test).
- 2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- 3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 4. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 7. This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.
- 8. An isotype control should be used at the same concentration as the antibody of interest.

References

Bigler RD, Bushkin Y, Chiorazzi N. S152 (CD27). A modulating disulfide-linked T cell activation antigen. *J Immunol*. 1988; 141(1):21-28. (Biology)
Bigler RD, Donat TL, Boselli CM. Definition of three epitopes of the CD27 molecule [P 120->55] present on activated normal lymphocytes. In: Knapp W, Dorken B, Rieber EP, et al, ed. *Leukocyte Typing IV: White Cell Differentiation Antigens*. New York: Oxford University Press; 1989:351-352. (Clone-specific)
Schlossman S, Boumell L, et al, ed. *Leucocyte Typing V*. New York: Oxford University Press; 1995. (Biology)
Watts TH. TNF/TNFR family members in costimulation of T cell responses. *Annu Rev Immunol*. 2005; 23:23-68. (Biology)

561786 Rev. 1 Page 2 of 2