Technical Data Sheet

PE-Cy[™]5 Mouse Anti-Human CD28

Product Information

Material Number: 561791

Alternate Name: CD28 antigen; T44; Tp44; TP44

25 Tests Size Vol. per Test: 20 µl CD28.2 Clone:


Immunogen: Human CD28 Transfected Cell Line Isotype: Mouse (C3H x BALB/c) IgG1, κ

Reactivity: QC Testing: Human Workshop: V 5T CD28.05

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The CD28.2 monoclonal antibody specifically binds to CD28, a 44 kDa homodimeric transmembrane glycoprotein present on most mature T cells, thymocytes and plasma cells. CD28 is a costimulatory receptor that binds CD80 and CD86 as ligands and plays a very important role in T cell-B cell interactions. It has been suggested that CD28 initiates and regulates a separate and distinct signal transduction pathway from those stimulated by the TCR complex. Additionally, it has been reported that CD28 antibody clones vary in their ability to stimulate T cells to produce IL-2 and increase intracellular Ca2+ concentration. This finding suggests the existence of functionally distinct subregions on the CD28 molecule. CD28.2 has been demonstrated to bind to the same molecule as clone L293, another CD28 mAb, and has been reported to induce Ca2+ influx in Jurkat T cells.

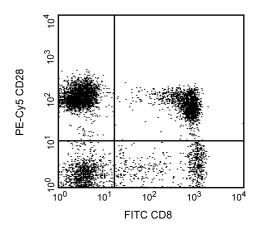

Profile of peripheral blood lymphocytes analyzed on a BD FACScan™ flow cytometer. (BDIS, San Jose, CA)

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with PE-Cy5 (formerly known as BD Cy-ChromeTM) under optimum conditions, and unconjugated antibody and free PE-Cy5 were removed.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name Name	Size	Clone	
555750	PE-Cy TM 5 Mouse IgG1 κ Isotype Control	100 Tests	MOPC-21	
554656	Stain Buffer (FBS)	500 mL	(none)	
555899	Lysing Buffer	100 mL	(none)	

BD Biosciences

bdbiosciences.com

United States Asia Pacific 866.979.9408 32.2.400.98.95 0120.8555.90 65.6861.0633 877.232.8995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD



Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-μl experimental sample (a test).
- 2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- 3. PE-Cy5 is a tandem fluorochrome composed of R-phycoerythrin (PE), which is excited by the 488 nm light of an Argon ion laser and serves as an energy donor, coupled to the cyanine dye Cy5, which acts as an energy acceptor and fluoresces at 670 nm. BD Biosciences Pharmingen has maximized the fluorochrome energy transfer in PE-Cy5, thus maximizing its fluorescence emission intensity, minimizing residual emission from PE, and minimizing lot-to-lot variation.
- 4. PE-Cy5 is optimized for use with a single argon ion laser emitting 488-nm light. Because of the broad absorption spectrum of the PE-Cy5 tandem fluorochrome, extra care must be taken when using dual-laser cytometers which may directly excite both PE and Cy5TM.
- 5. PE-Cy5 tandem fluorochromes have been reported to bind some classes of human macrophages and granulocytes via Fc receptors, and PE has been reported to bind to mouse B lymphocytes via Fc receptors. Preincubation of mouse leukocytes with Mouse BD Fc BlockTM purified anti-mouse CD16/CD32 mAb 2.4G2 can reduce the non-specific binding of PE-Cy5-conjugated reagents to mouse B cells. However, PE-Cy5 conjugated reagents should not be used to stain splenocytes of SJL, NOD, and MRL mice as B lymphocytes and/or other leukocytes have been reported to non-specifically stain regardless of the use of Mouse BD Fc BlockTM (the CD72c complex has been implicated for PE-Cy5 binding in these strains). Reagents conjugated to PE, PerCP, PerCP-Cy5.5, APC, and APC-Cy7 tandem fluorochrome can be used on leukocytes from these mouse strains.
- 6. Cy is a trademark of Amersham Biosciences Limited.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- 8. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- 9. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- 10. Please observe the following precautions: Absorption of visible light can significantly alter the energy transfer occurring in any tandem fluorochrome conjugate; therefore, we recommend that special precautions be taken (such as wrapping vials, tubes, or racks in aluminum foil) to prevent exposure of conjugated reagents, including cells stained with those reagents, to room illumination.
- 11. This product is sold under license. Purchase of this product does not include rights to (i) incorporate this product into the purchaser's own products for resale to end-users, or (ii) use this product to conduct for-profit research for or on behalf of another party. For information on obtaining a license to this product for such prohibited uses, contact INSERM, 7 rue Watt, 75013 Paris. Telephone: +33 1 55 03 01 60. Facsimile: +33 1 55 03 01 18. Email: techtransfert@inserm-transfert.fr
- 12. An isotype control should be used at the same concentration as the antibody of interest.

References

Barclay NA, Brown MH, Birkeland ML, et al, ed. The Leukocyte Antigen FactsBook. San Diego, CA: Academic Press; 1997. (Biology)

June CH, Bluestone JA, Nadler LM, Thompson CB. The B7 and CD28 receptor families. Immunol Today. 1994; 15(7):321-331. (Biology)

Kuiper H, Brouwer M, Vermeire S, van Lier R. Analysis of the Workshop CD28 Panel mAb: distinct signalling pathways coupled to CD28. In: Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leucocyte Typing V: White Cell Differentiation Antigens*. Oxford: Oxford University Press; 1995:373-374. (Clone-specific: Activation, Calcium Flux, (Co)-stimulation)

Nunes J, Klasen S, Franco MD, et al. Signalling through CD28 T-cell activation pathway involves an inositol phospholipid-specific phospholipase C activity. Biochem J. 1993; 293(3):835-842. (Biology)

Nunes J, Klasen S, Ragueneau M, et al. CD28 mAbs with distinct binding properties differ in their ability to induce T cell activation: analysis of early and late activation events. *Int Immunol.* 1993; 5(3):311-315. (Biology)

Olive D, Cerdan C, Costello R, Sielleur I, Ragueneau M, Pages F, Klasen S, Nunes J, Imbert J. CD28 and CTLA-4 cluster report. In: Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leucocyte Typing V: White Cell Differentiation Antigens*. Oxford: Oxford University Press; 1995:360-370. (Clone-specific: (Co)-stimulation, Flow cytometry, Functional assay, Inhibition, Stimulation)

Schlossman SF, Boumsell L, Gilks W, et al, ed. Leukocyte Typing V: White Cell Differentiation Antigens. Oxford: Oxford University Press; 1995. (Clone-specific) Verwilghen J, Vandenberghe P, Wallays G, et al. Simultaneous ligation of CD5 and CD28 on resting T lymphocytes induces T cell activation in the absence of T cell receptor/CD3 occupancy. J Immunol. 1993; 150(3):835-846. (Biology)

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 866.979.9408
 32.2.400.98.95
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit ${\bf bdbiosciences.com/contact}$

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.



Page 2 of 2